



US008262645B2

(12) **United States Patent**
Bagwell et al.

(10) **Patent No.:** **US 8,262,645 B2**

(45) **Date of Patent:** **Sep. 11, 2012**

(54) **DEVICES FOR CLEARING BLOCKAGES IN IN-SITU ARTIFICIAL LUMENS**

(75) Inventors: **Roger B. Bagwell**, Bellefonte, PA (US); **Paul L. Frankhouser**, Miami Beach, FL (US); **Maureen L. Mulvihill**, Bellefonte, PA (US); **Josue R. Crespo**, University Park, PA (US); **Brian M. Park**, Bellefonte, PA (US); **Dana B. Mallen**, State College, PA (US); **Gabriela Hernandez Meza**, State College, PA (US); **Debora L. Demers**, State College, PA (US)

(73) Assignee: **Actuated Medical, Inc.**, Bellefonte, PA (US)

(*) Notice: Subject to any disclaimer, the term of this patent is extended or adjusted under 35 U.S.C. 154(b) by 0 days.

(21) Appl. No.: **12/964,252**

(22) Filed: **Dec. 9, 2010**

(65) **Prior Publication Data**

US 2011/0106019 A1 May 5, 2011

Related U.S. Application Data

(63) Continuation-in-part of application No. 12/274,937, filed on Nov. 20, 2008, now abandoned.

(60) Provisional application No. 61/099,737, filed on Sep. 24, 2008, provisional application No. 60/989,484, filed on Nov. 21, 2007.

(51) **Int. Cl.**
A61M 5/00 (2006.01)
A61M 5/46 (2006.01)

(52) **U.S. Cl.** **604/540**; 604/541; 604/35; 604/43; 15/104.05

(58) **Field of Classification Search** 604/540, 604/541, 35, 43

See application file for complete search history.

(56) **References Cited**

U.S. PATENT DOCUMENTS

3,765,420 A 10/1973 Felczak
(Continued)

FOREIGN PATENT DOCUMENTS

EP 1666168 6/2006
(Continued)

OTHER PUBLICATIONS

Dolmatch, et al., "Synthetic Dialysis Shunts: Thrombolysis with the Cragg Thrombolytic Brush Catheter." Radiology 1999; 213:180-184.

(Continued)


Primary Examiner — Michael Kornakov
Assistant Examiner — Natasha Campbell

(74) *Attorney, Agent, or Firm* — Metz Lewis Brodman Must O'Keefe LLC; Barry L. Friedman

(57) **ABSTRACT**

Devices and methods for the effective clearing of artificial tubes, especially in-situ clearing of artificial tubes in a living being are covered in this disclosure. The devices and methods provide an elongated clearing member having a first end that is coupled to a driving mechanism and having a second working end that is subjected to repetitive motion for clearing blockages within the artificial tube. The elongated clearing member includes either a fixed or an adjustable element that selectively defines the portion of the elongated clearing member that is insertable within the artificial tube. The proximal end of the clearing member is releasably secured to the driving mechanism and the driving mechanisms may comprise a wide variety of repetitive motion drivers such as voice coil motors, piezoelectric actuators, pneumatic actuators, DC motors, etc. These devices/methods may comprise a free-standing console for hands-free operation or may comprise hand-held versions. The distal working end of the clearing member may comprise tips of differing functions, including an irrigation/aspiration feature.

70 Claims, 44 Drawing Sheets

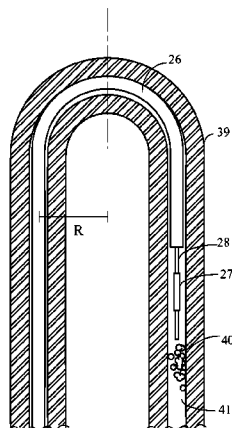

U.S. PATENT DOCUMENTS

3,818,913 A 6/1974 Wallach
 4,363,980 A 12/1982 Petersen
 4,698,058 A 10/1987 Greenfeld et al.
 4,800,104 A 1/1989 Cruickshank
 4,845,688 A 7/1989 Butler
 4,850,071 A 7/1989 Lawrence
 4,894,056 A 1/1990 Bommarito
 4,999,819 A 3/1991 Newnham et al.
 5,197,946 A 3/1993 Tachibana
 5,240,675 A 8/1993 Wilk et al.
 5,251,356 A 10/1993 Oaki et al.
 5,269,291 A 12/1993 Carter
 5,271,735 A 12/1993 Greenfeld et al.
 5,276,657 A 1/1994 Newnham et al.
 5,315,998 A 5/1994 Tachibana et al.
 5,318,014 A 6/1994 Carter
 5,324,282 A 6/1994 Dodick
 5,362,309 A 11/1994 Carter
 5,375,589 A 12/1994 Bhatta
 5,395,398 A 3/1995 Rogozinski
 5,424,299 A 6/1995 Monte
 5,431,640 A 7/1995 Gabriel
 5,458,633 A 10/1995 Bailey
 5,512,048 A 4/1996 Slettenmark
 5,615,439 A 4/1997 Bourrelly
 5,636,403 A 6/1997 Grimsley et al.
 5,664,594 A 9/1997 Kee
 5,725,494 A 3/1998 Briskin
 5,729,077 A 3/1998 Newnham et al.
 5,868,708 A 2/1999 Hart et al.
 5,897,534 A 4/1999 Heim et al.
 5,906,611 A 5/1999 Dodick et al.
 6,047,431 A 4/2000 Canonica
 6,126,647 A 10/2000 Posey et al.
 6,173,199 B1 1/2001 Gabriel
 6,187,266 B1 2/2001 Lin et al.
 6,326,340 B1 12/2001 Labib et al.
 6,439,246 B2 8/2002 Stanley
 6,460,231 B2 10/2002 Bourgerie
 6,465,936 B1 10/2002 Knowles et al.
 6,524,251 B2 2/2003 Rabiner et al.
 6,551,337 B1 4/2003 Rabiner et al.
 6,619,302 B2 9/2003 Labib et al.
 6,652,547 B2 11/2003 Rabiner et al.
 6,660,013 B2 12/2003 Rabiner et al.
 6,681,783 B2 1/2004 Kawazoe
 6,699,331 B1 3/2004 Kritzler
 6,725,492 B2 4/2004 Moore et al.
 6,866,670 B2 3/2005 Rabiner et al.
 7,055,203 B1 6/2006 Franzino
 7,081,133 B2 7/2006 Chinn et al.
 7,393,501 B2 7/2008 Zumeris et al.
 2002/0069893 A1* 6/2002 Kawazoe 134/1
 2002/0099387 A1* 7/2002 Gauderer et al. 606/106
 2003/0109837 A1 6/2003 McBride-Sakal
 2003/0140946 A1 7/2003 Coats
 2003/0181876 A1 9/2003 Ahn et al.
 2003/0215433 A1 11/2003 Kokai-Kun et al.
 2004/0084070 A1 5/2004 Sasaki et al.

2004/0147595 A1 7/2004 Kjelleberg et al.
 2004/0181194 A1 9/2004 Perkins
 2004/0230117 A1 11/2004 Tosaya et al.
 2006/0293559 A1 12/2006 Grice, III et al.
 2007/0049904 A1 3/2007 Deutsch
 2007/0179071 A1 8/2007 Thoele
 2007/0179513 A1 8/2007 Deutsch
 2007/0215190 A1 9/2007 Efinger et al.
 2007/0225615 A1* 9/2007 Chechelski et al. 600/585
 2007/0244423 A1 10/2007 Zumeris et al.
 2008/0141473 A1 6/2008 Arai et al.
 2008/0267688 A1 10/2008 Busted
 2008/0283065 A1 11/2008 Chang et al.
 2009/0036341 A1 2/2009 Carney et al.
 2009/0054881 A1 2/2009 Krespi
 2009/0264833 A1 10/2009 Boyle, Jr.
 2010/0186783 A1 7/2010 Knaff et al.

FOREIGN PATENT DOCUMENTS

EP 1829474 9/2007
 WO WO 2004098654 A2 * 11/2004
 WO 2006/063176 6/2006
 WO 2007/110870 10/2007

OTHER PUBLICATIONS

Cauda, et al., "Heparin Coating on Ureteral Double J Stents Prevents Encrustations: An in Vivo Case Study", Journal of Endourology, Mar. 2008, 22(3): 465-472. DOI: 10.1089/end.2007.0218.
 Kite, et al., "Use of in Vivo-Generated Biofilms from Hemodialysis Catheters to Test the Efficacy of a Novel Antimicrobial Catheter Lock for Biofilm Eradication In Vitro", J Clin Microbiol. Jul. 2004; 42(7): 3073-3076. doi: 10.1128/JCM.42.7.3073-3076.2004. PMID: PMC446317 Copyright © 2004, American Society for Microbiology. Viasys (Corpak) Medsystems, "Clog Zapper™", (Wheeling, IL; 800-323-6305), Dated 2008.
 Fish, J., "Tube feeding in the ICU: Overcoming the Obstacles", www.rosslearningcenter.com/library/Tube%20Feeding%20in%20the%20ICU_RX.pdf (Oct. 14, 2004).
 Nicholson, L., "Declogging Small-Bore Feeding Tubes", Journal of Parenteral and Enteral Nutrition, vol. 11, No. 6, 594-597 (1987) DOI: 10.1177/0148607187011006594.
 Smith, R., et al., (Web Article) "2 Devices That Unclog Feeding Tubes", http://rn.modernmedicine.com/rnweb/article/articleDetail.jsp?id=142656, Published Jan. 1, 2005.
 Barnadas, G., "Navigating Home Care: Enteral Nutrition—Part One", Practical Gastroenterology, Oct. 2003. www.oralcancerfoundation.org/dental.tube-feeding.htm, "Tube Feeding", Dated 2001.
 Schaffner, et al., "Long-Term Enteral Nutrition Facilities Optimization of Body Weight", Journal of Parenteral and Enteral Nutrition, vol. 29, No. 3, 2005.
 Vickery, et al., "Removal of Biofilm from Endoscopes: Evaluation of Detergent Efficiency", American Journal of Infection Control, vol. 32, No. 3, May 2004.
 International Search Report for corresponding PCT Application No. PCT/US2010/061900, dated May 5, 2011.

* cited by examiner

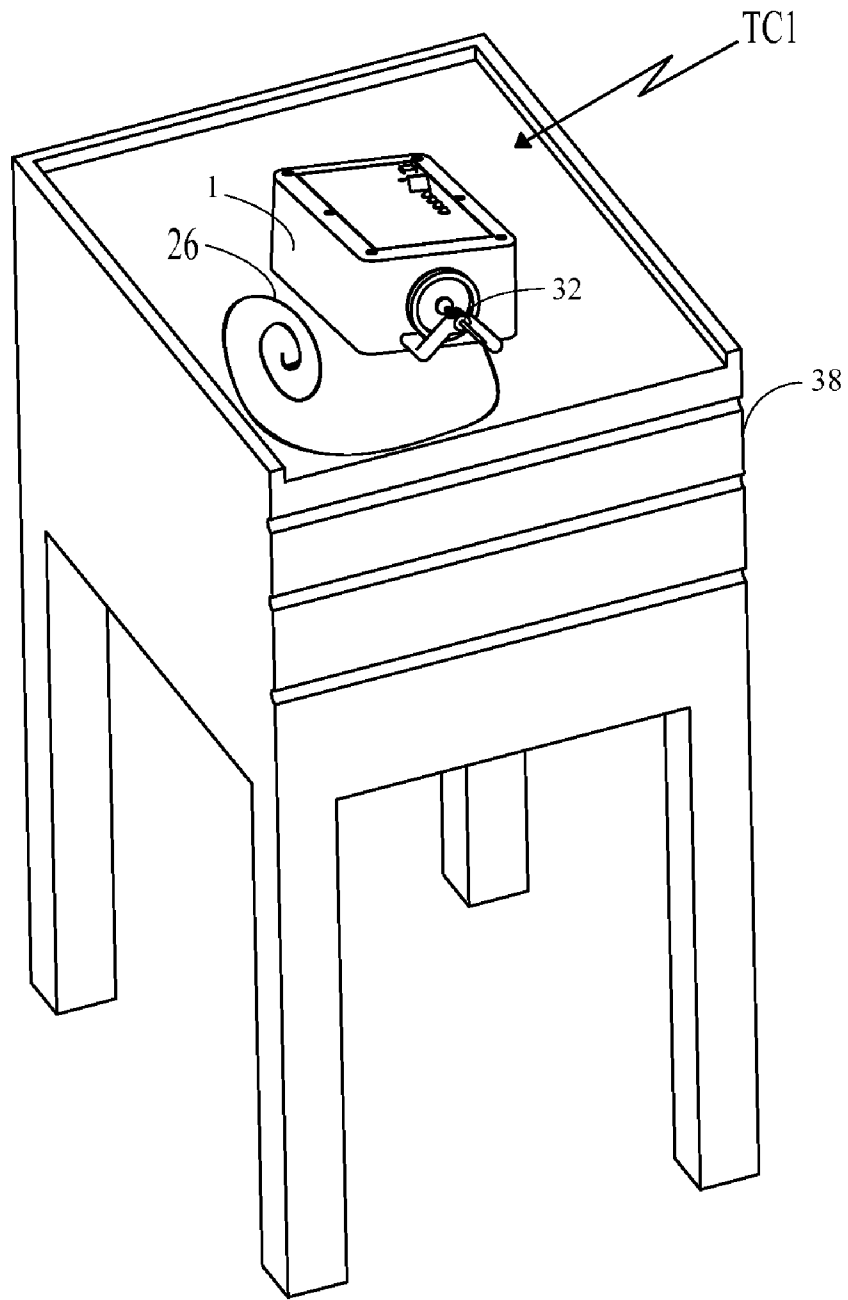


FIG. 1

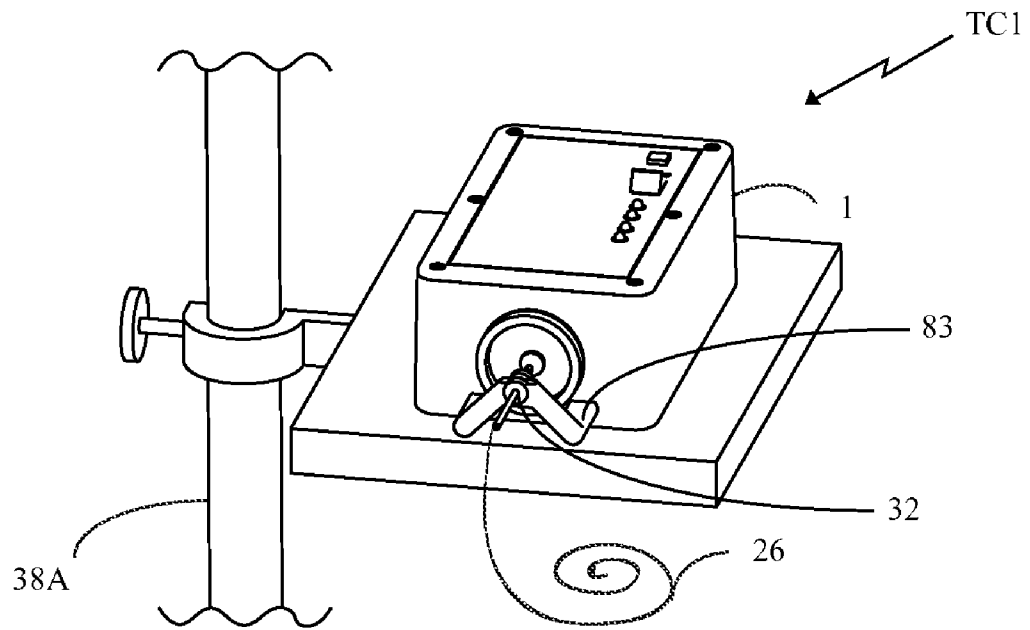


FIG. 1A

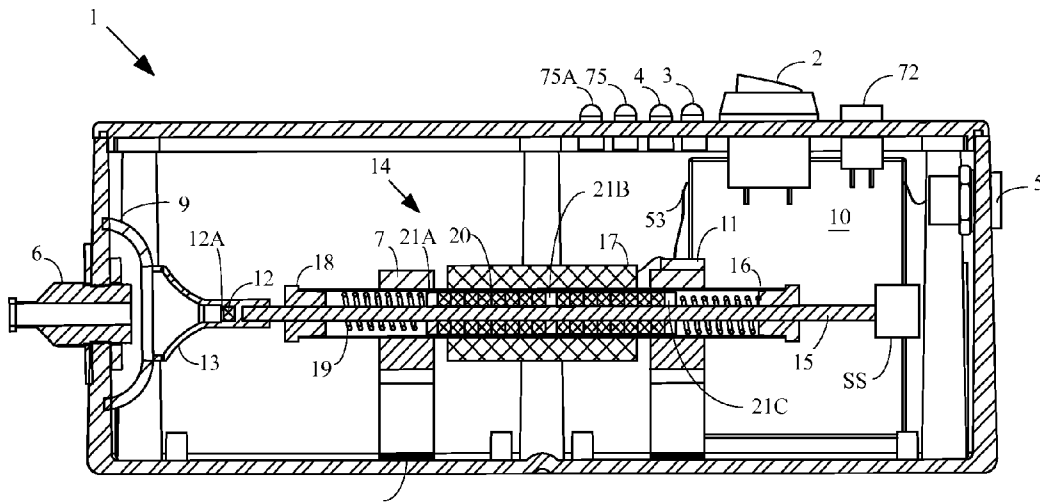


FIG. 2A

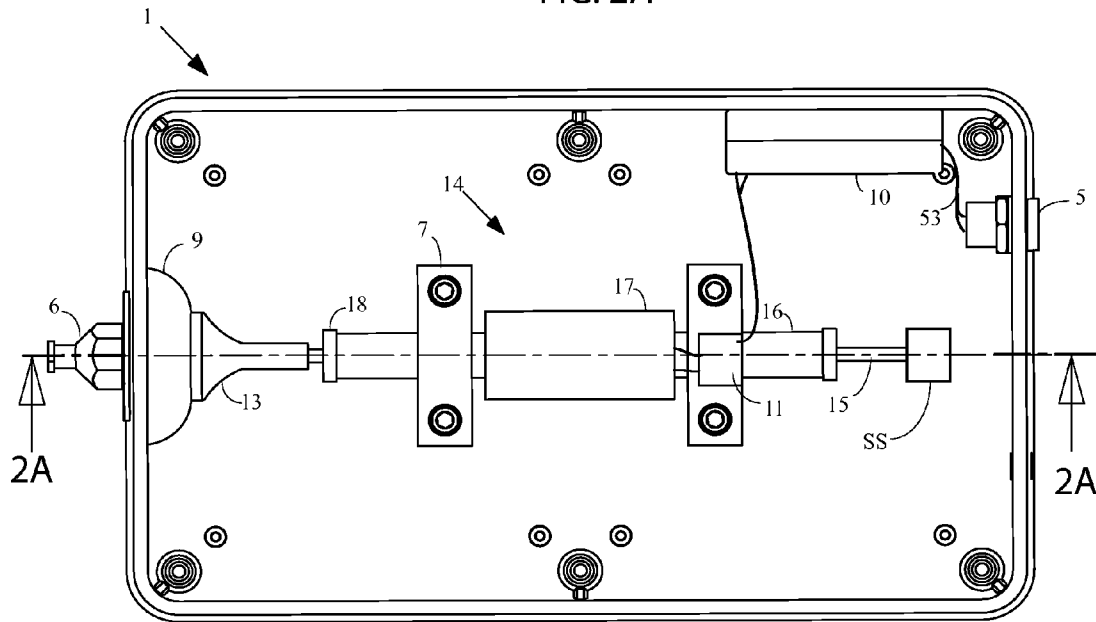


FIG. 2

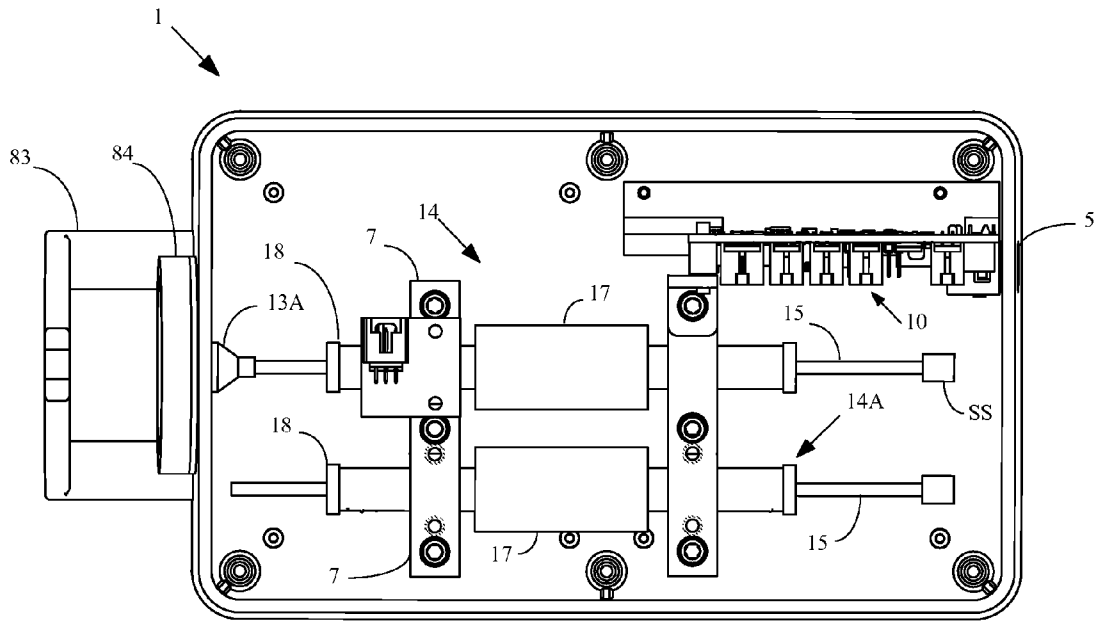


FIG. 2B

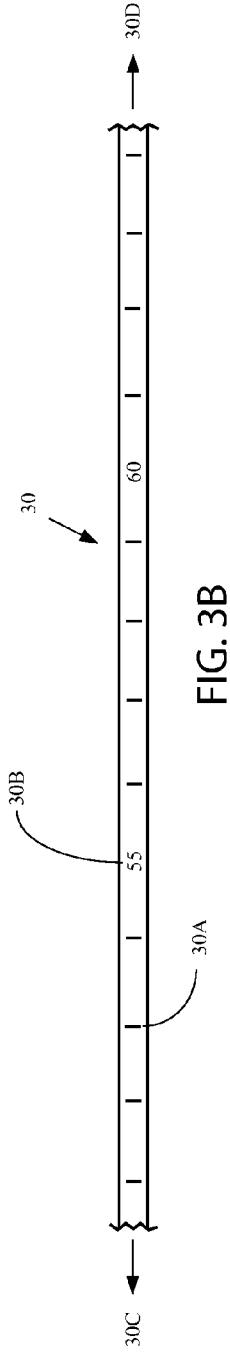


FIG. 3B

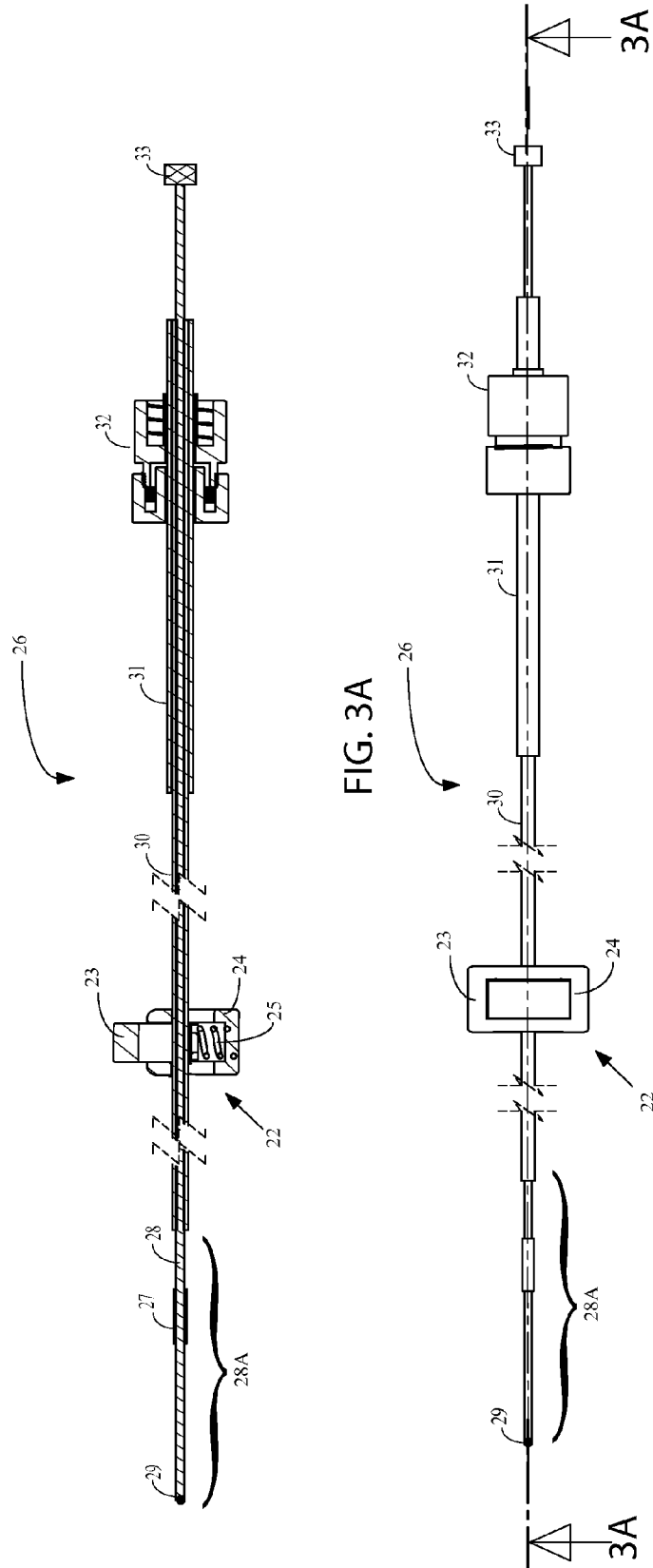


FIG. 3A

FIG. 3

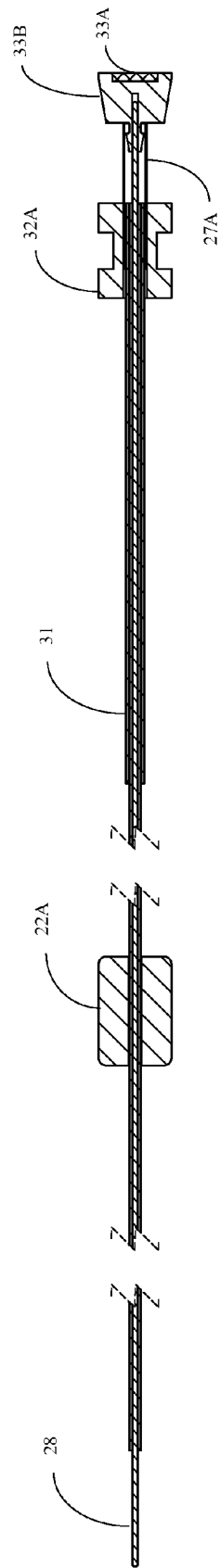


FIG. 3D

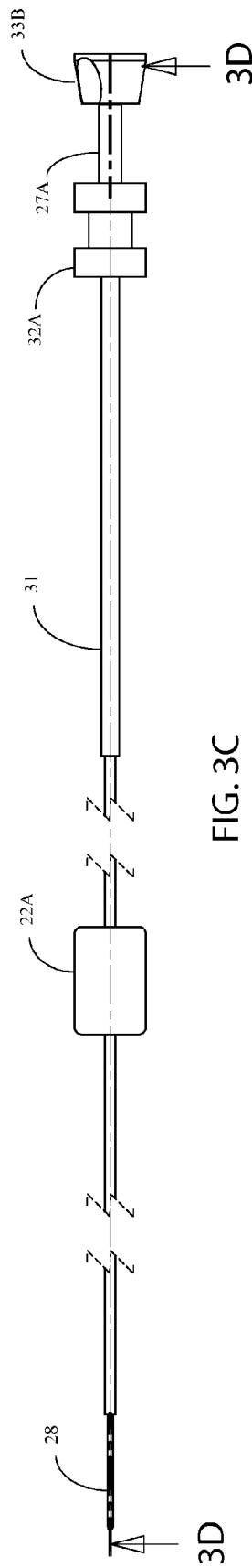


FIG. 3C

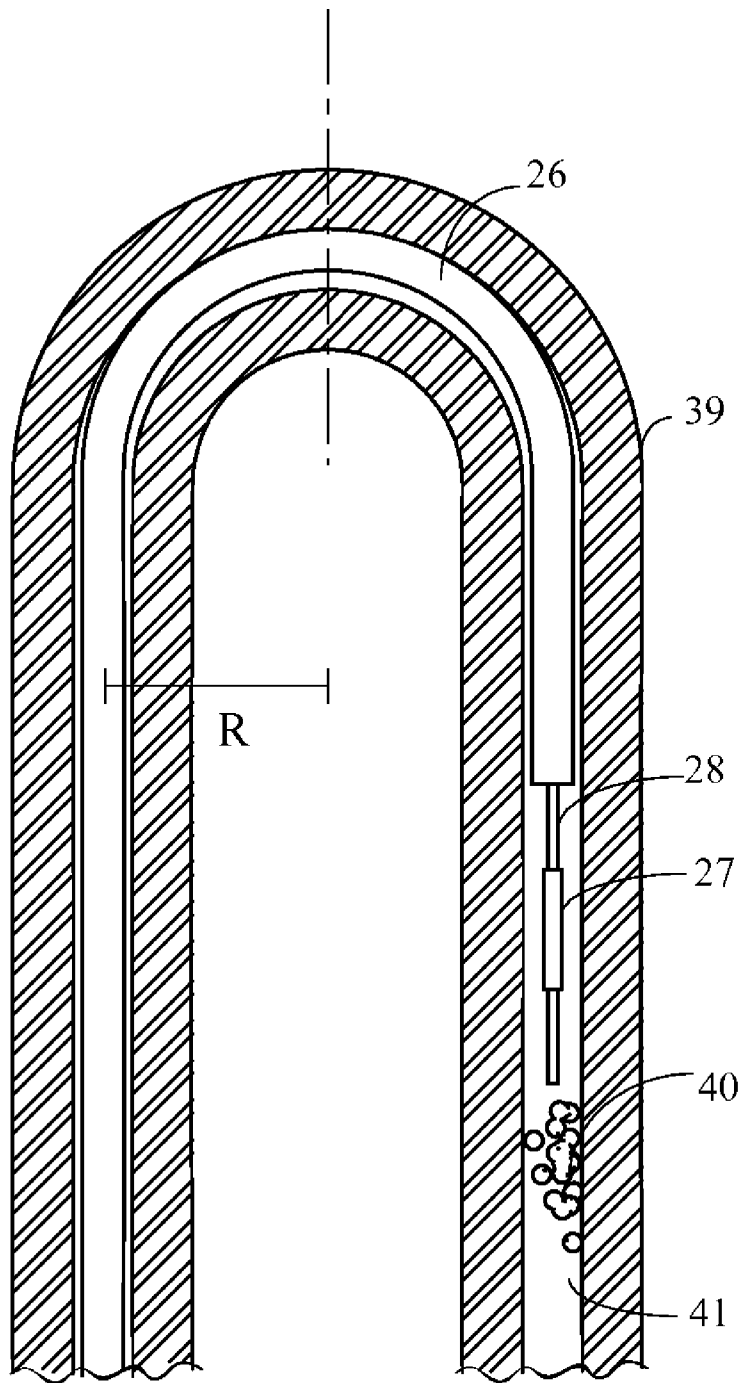


FIG. 4

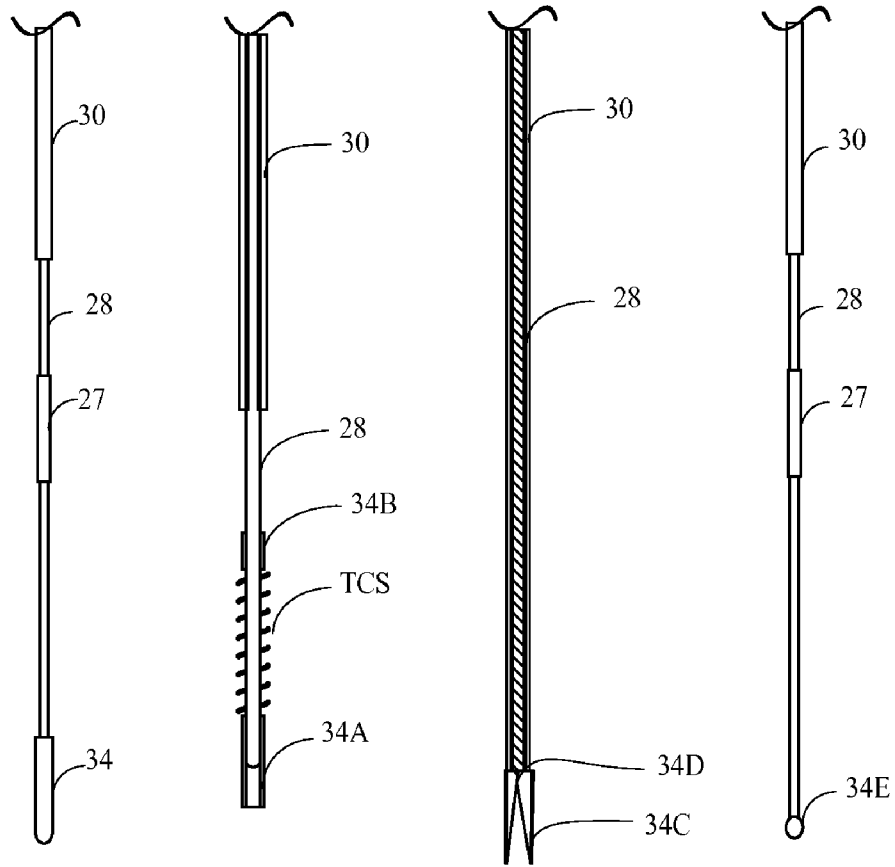


FIG. 5A

FIG. 5B

FIG. 5C

FIG. 5D

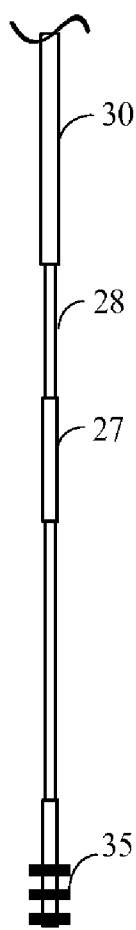


FIG. 6

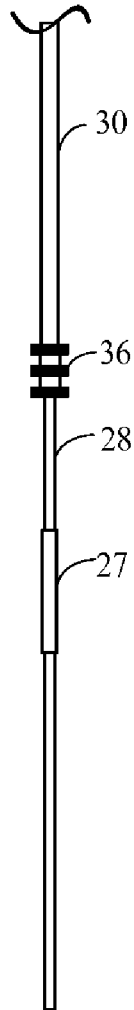


FIG. 7

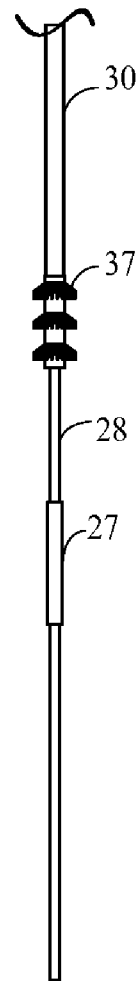


FIG. 8

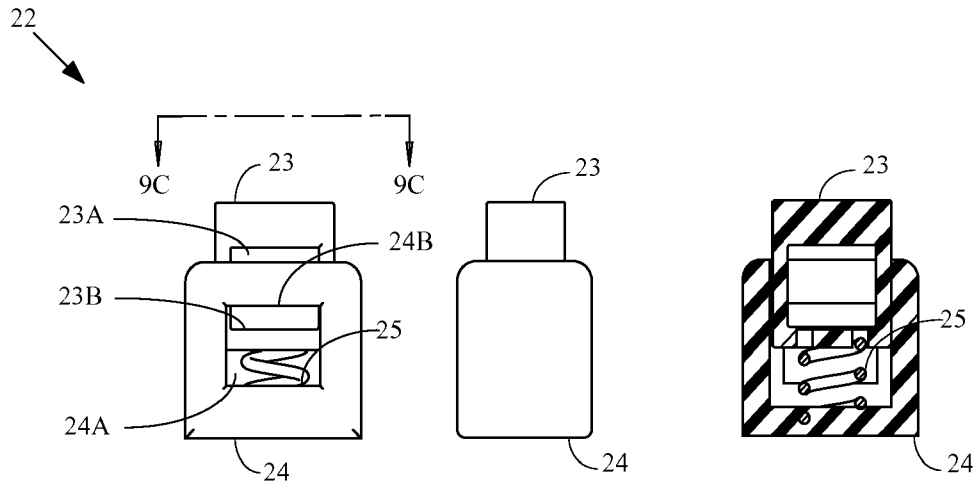


FIG. 9A

FIG. 9B

FIG. 9C

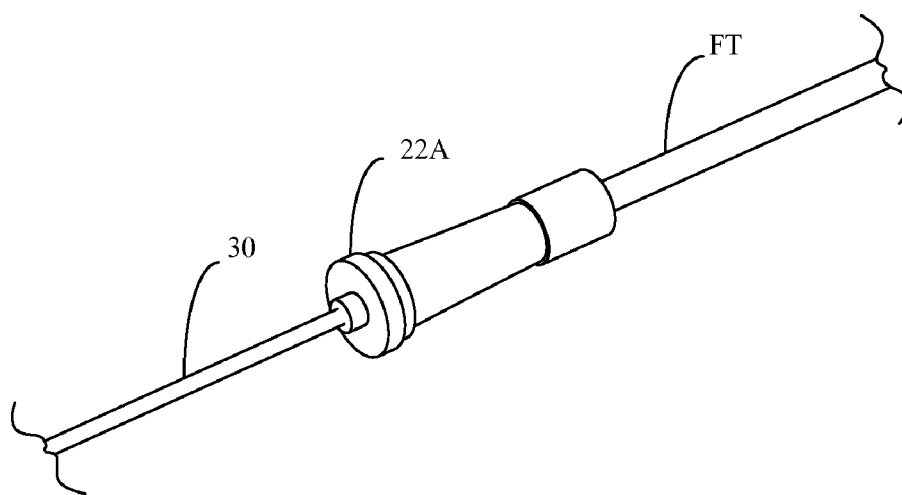


FIG. 9D

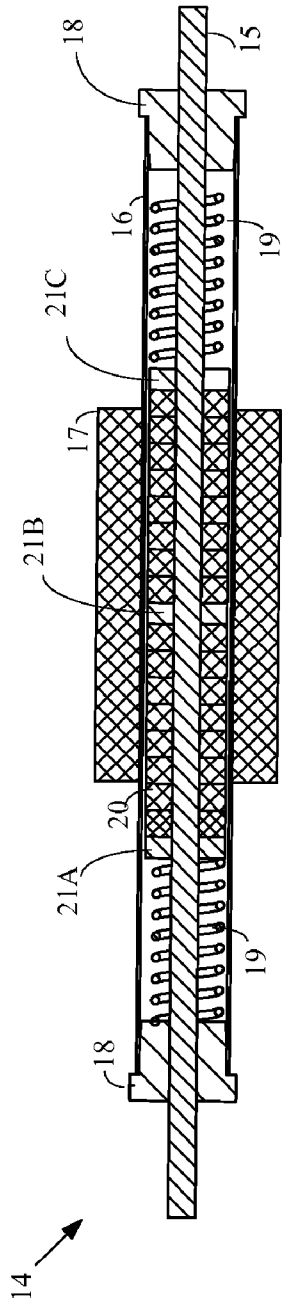


FIG. 10A

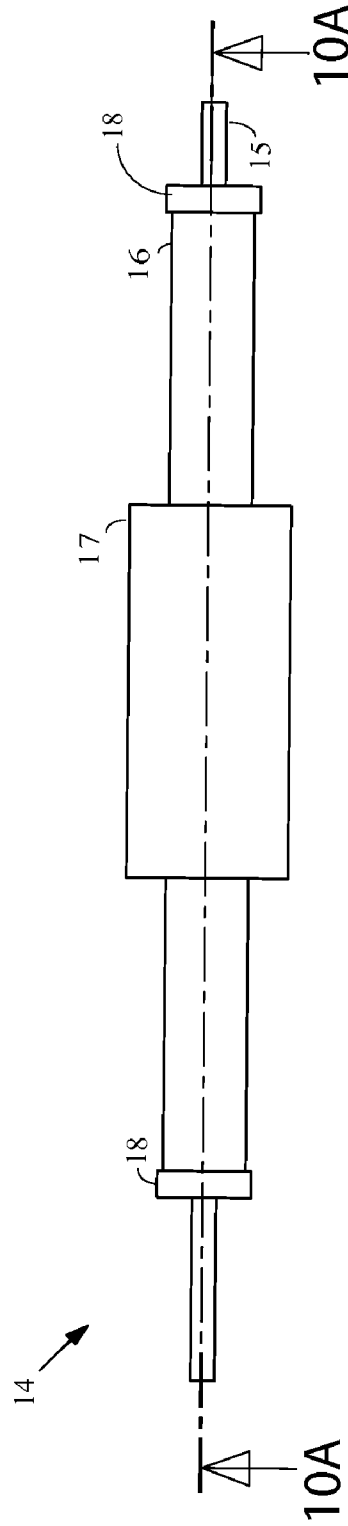


FIG. 10

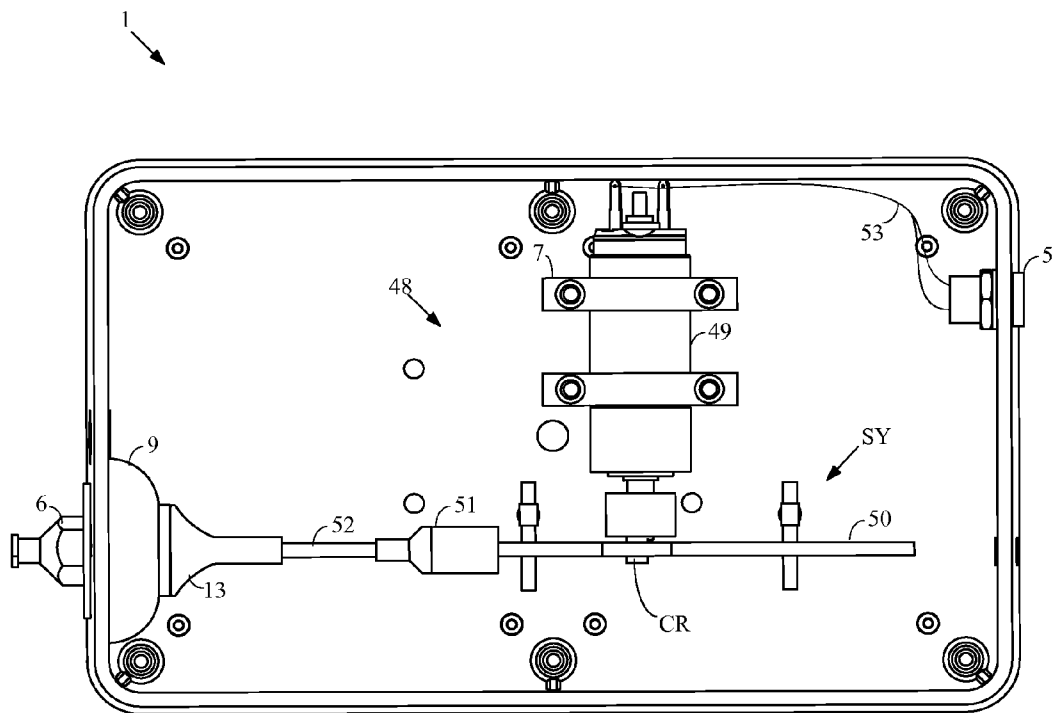


FIG. 11

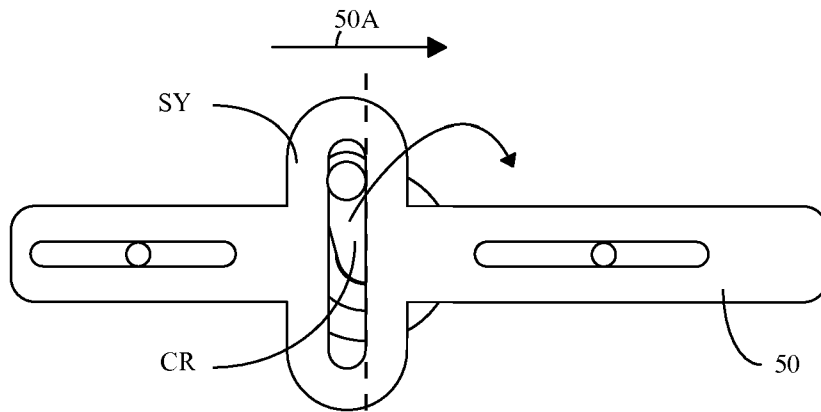


FIG. 11A

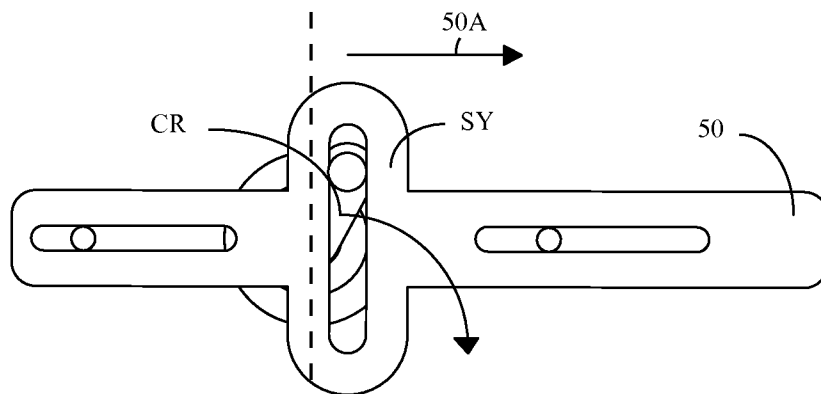


FIG. 11B

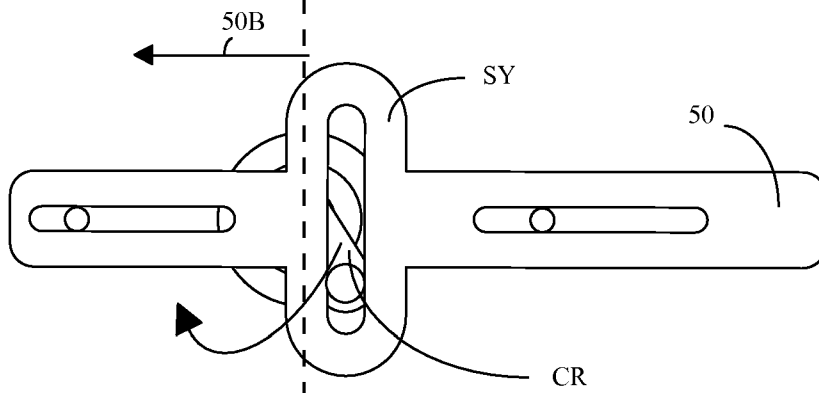


FIG. 11C

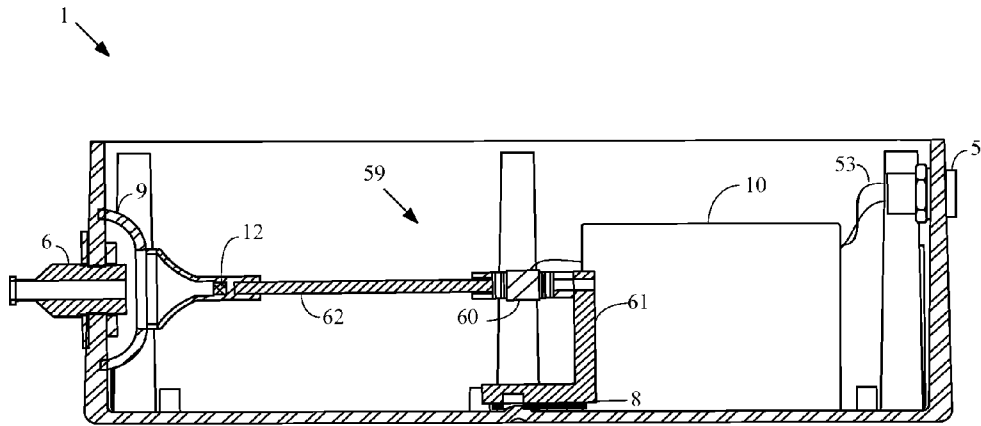


FIG. 12A

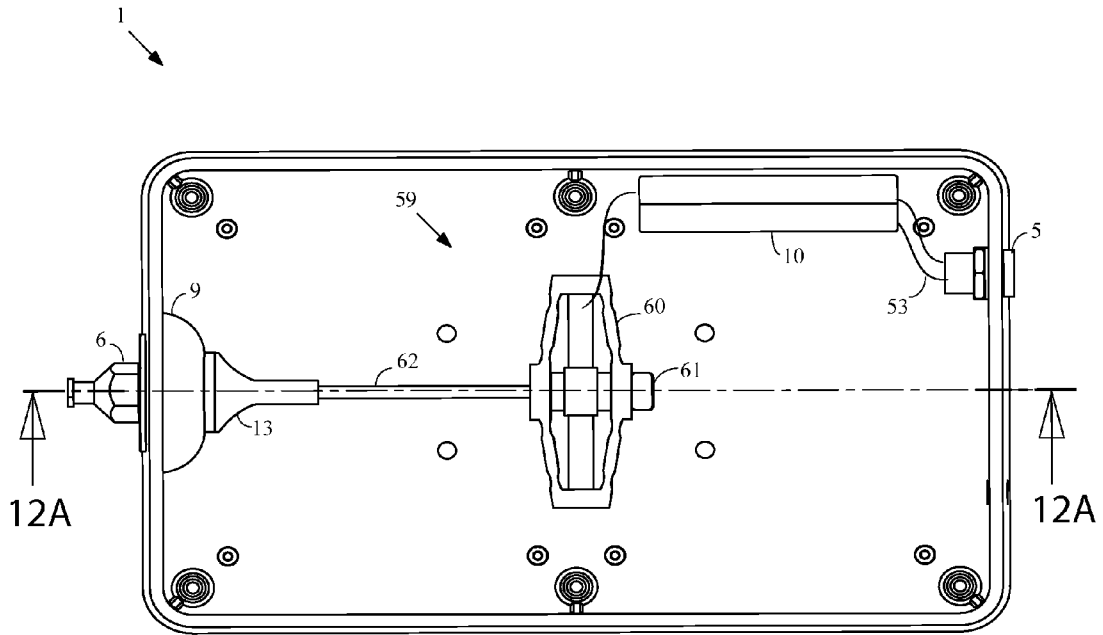


FIG. 12

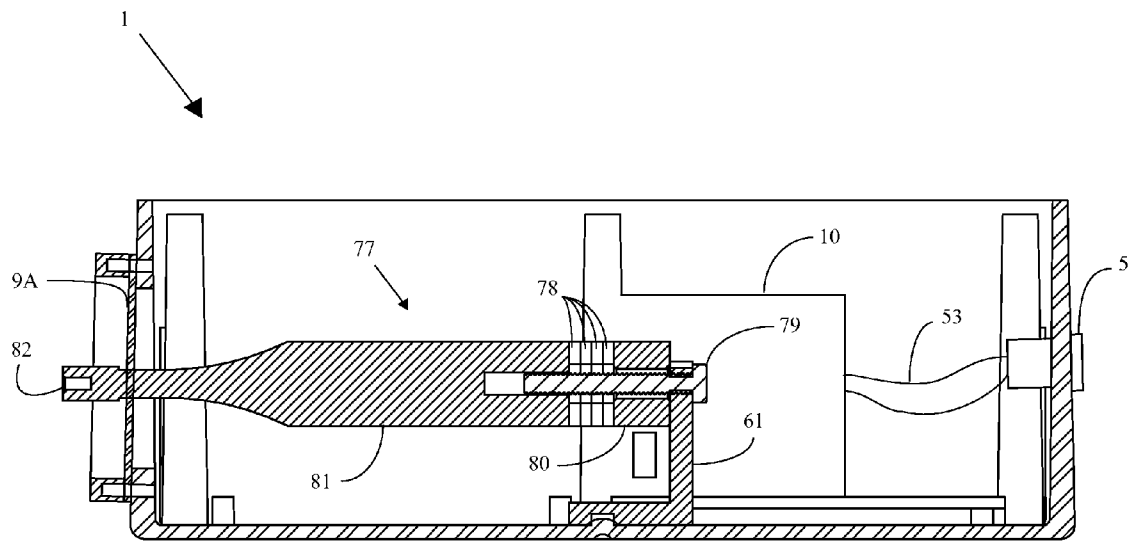


FIG. 12B

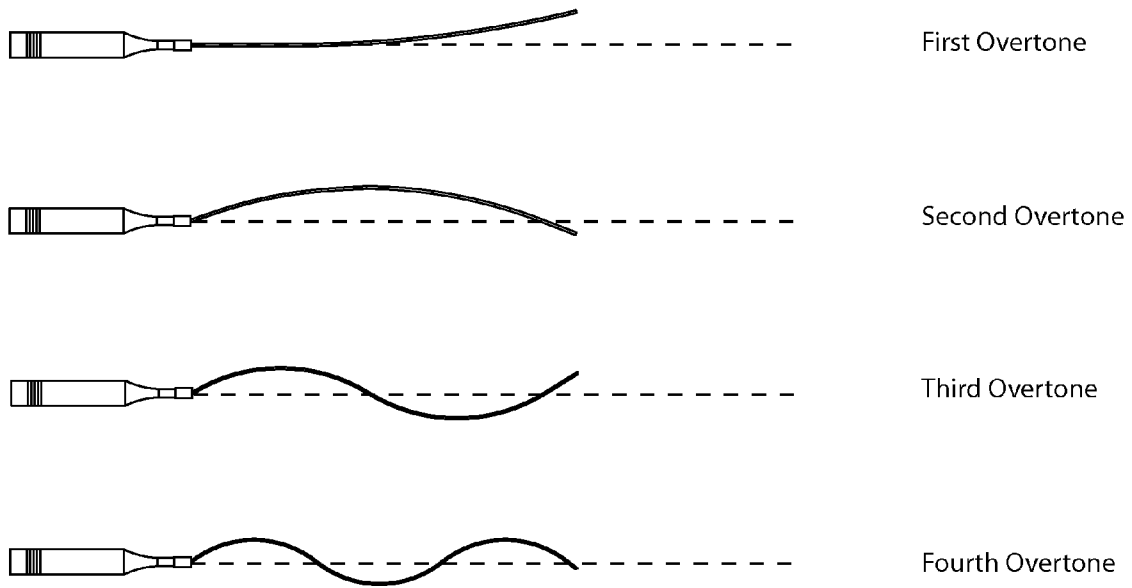


FIG. 12C

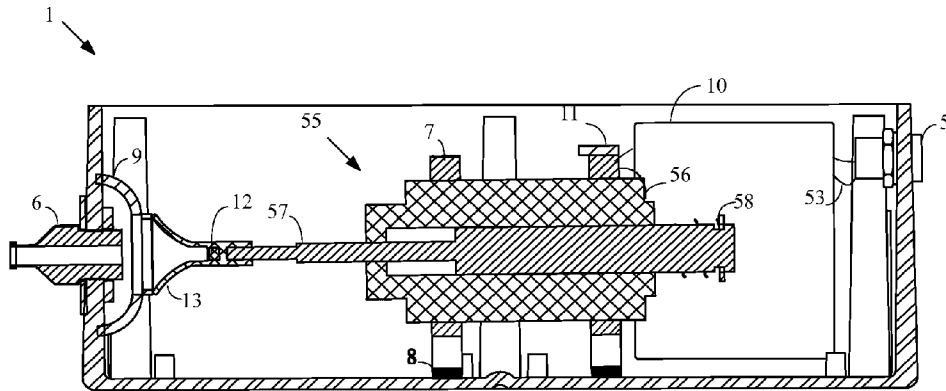


FIG. 13A

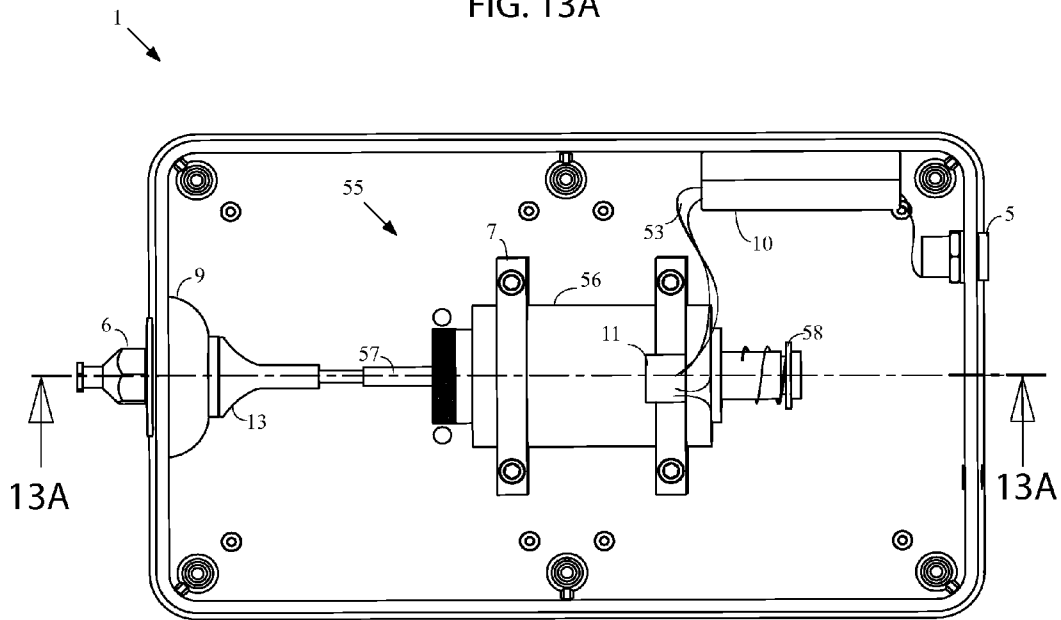


FIG. 13

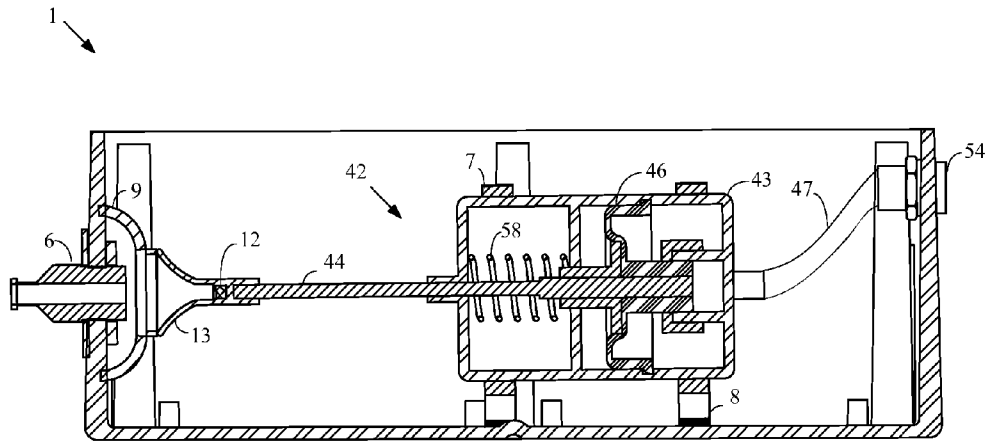


FIG. 14A

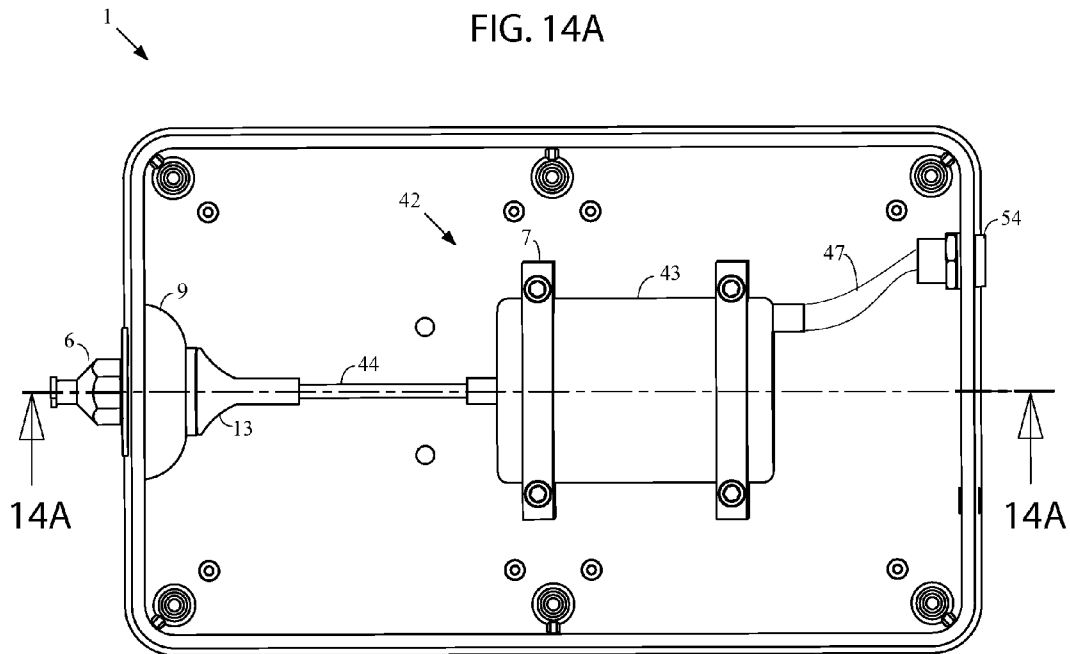


FIG. 14

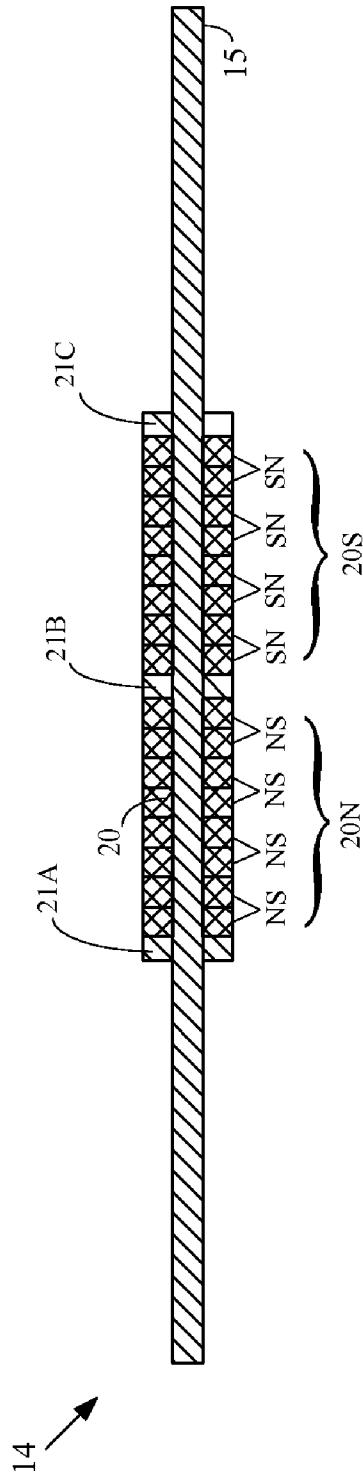


FIG. 15

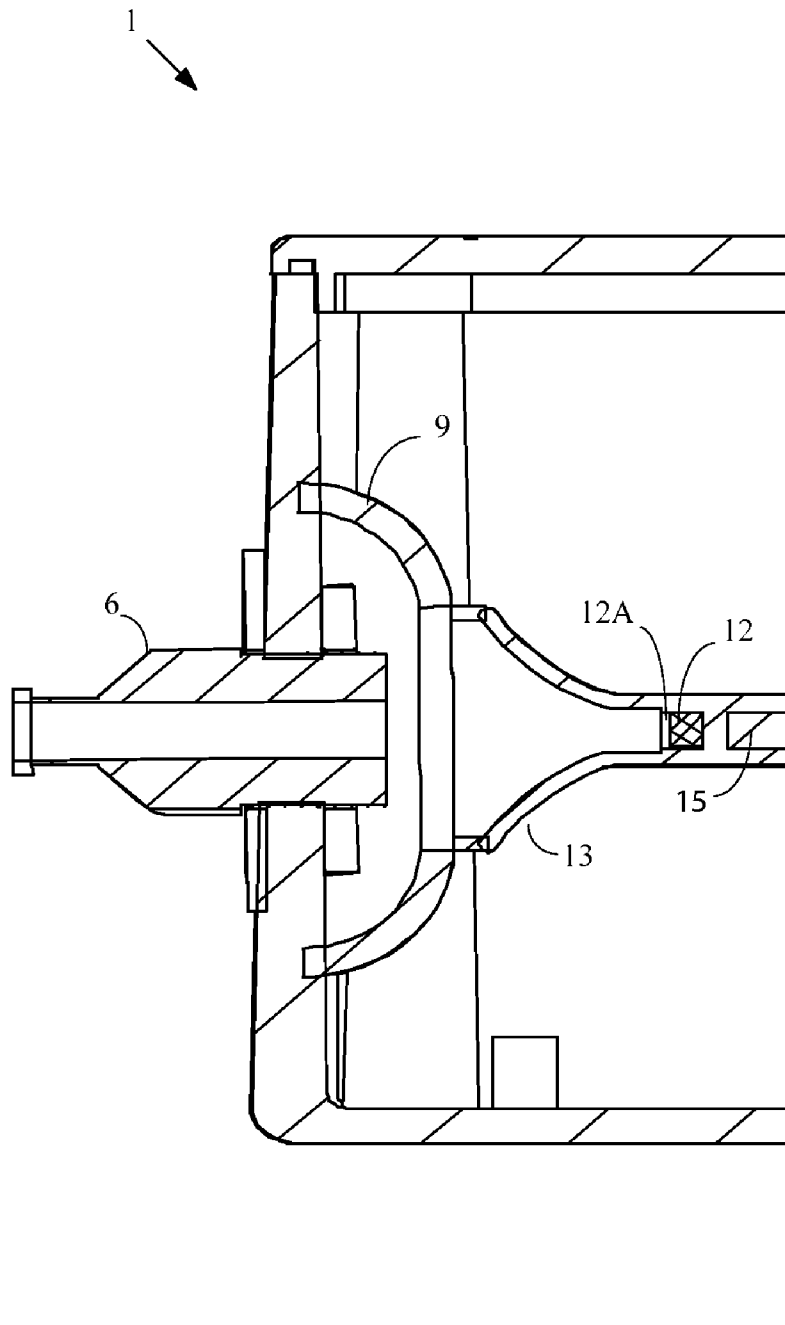


FIG. 16A

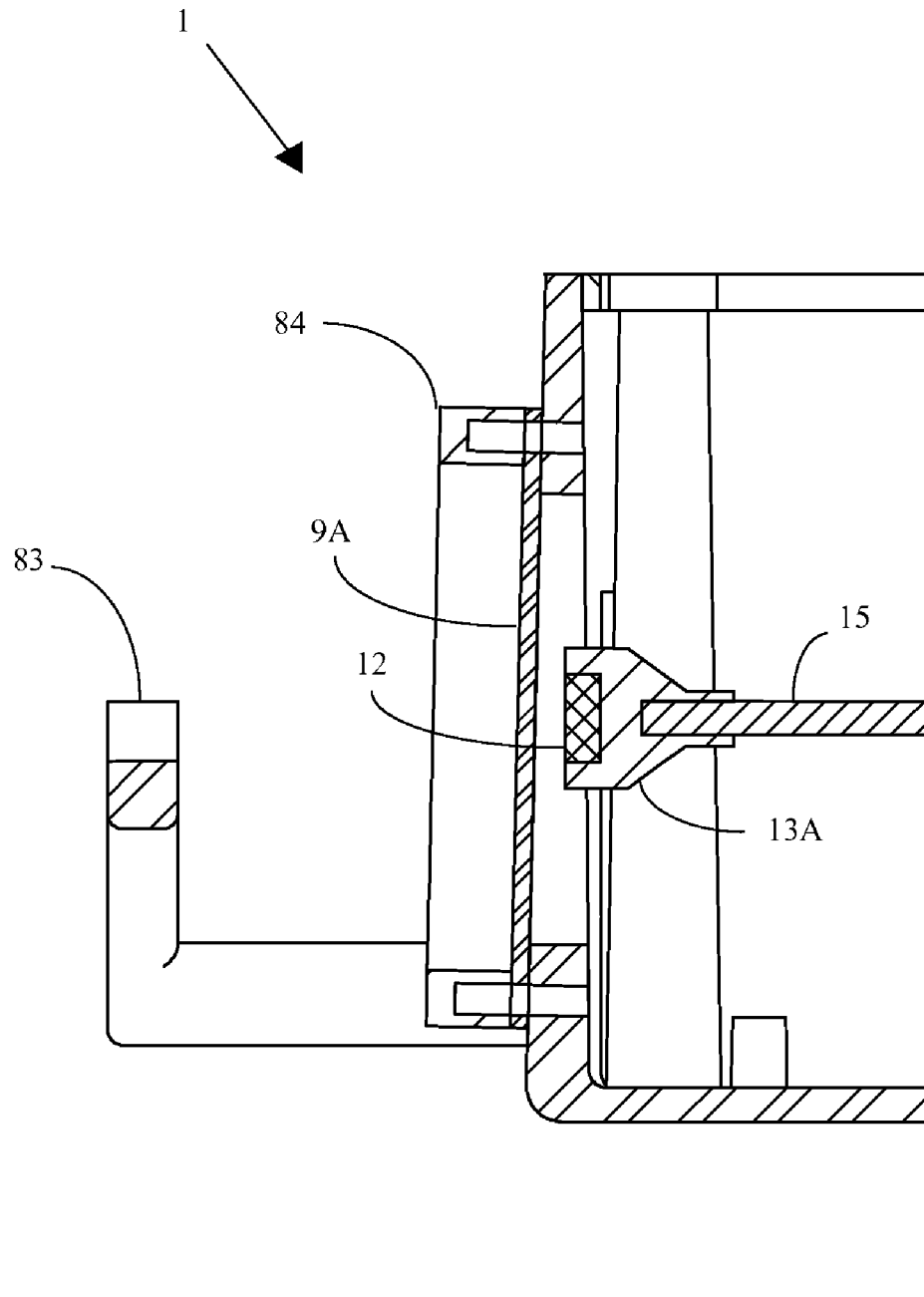


FIG. 16B

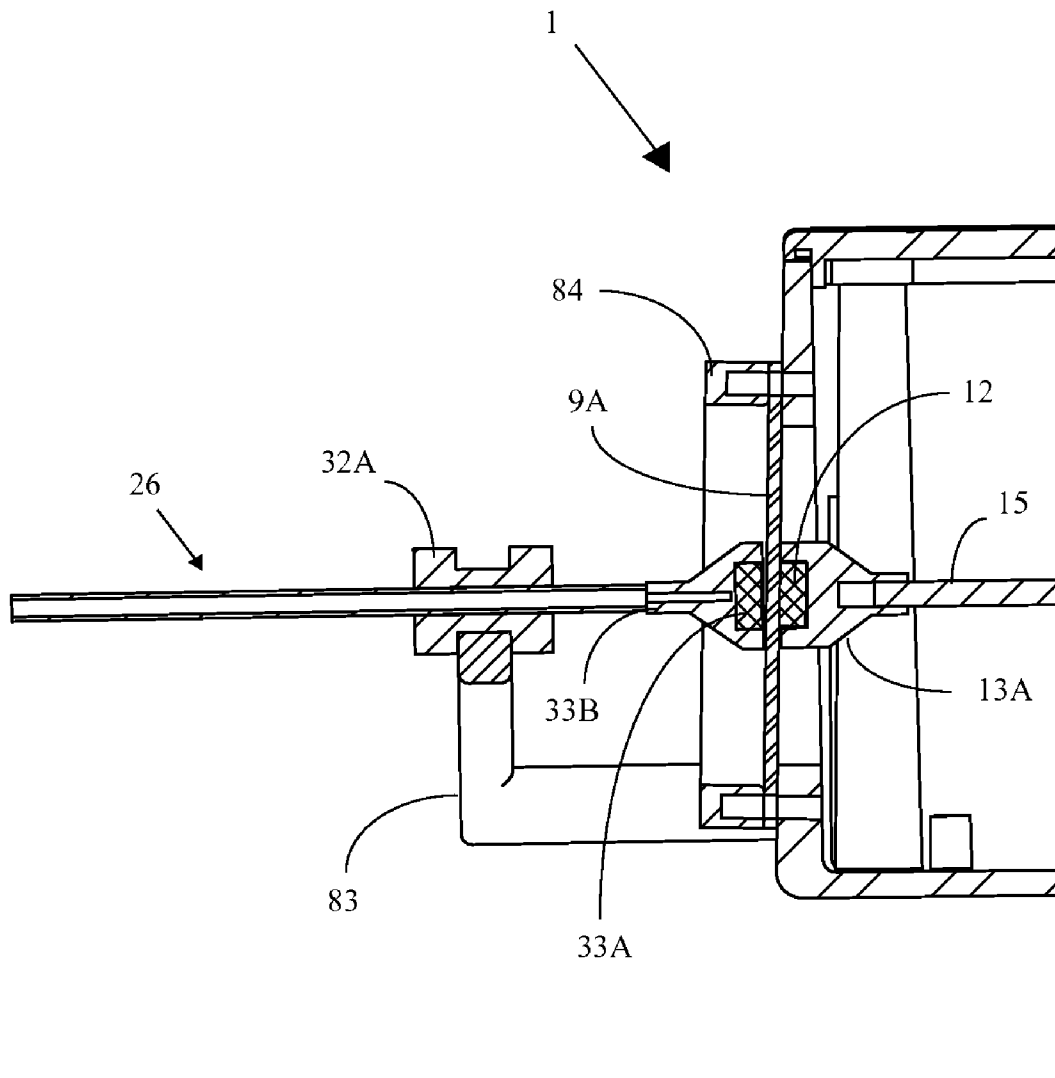


FIG. 16C

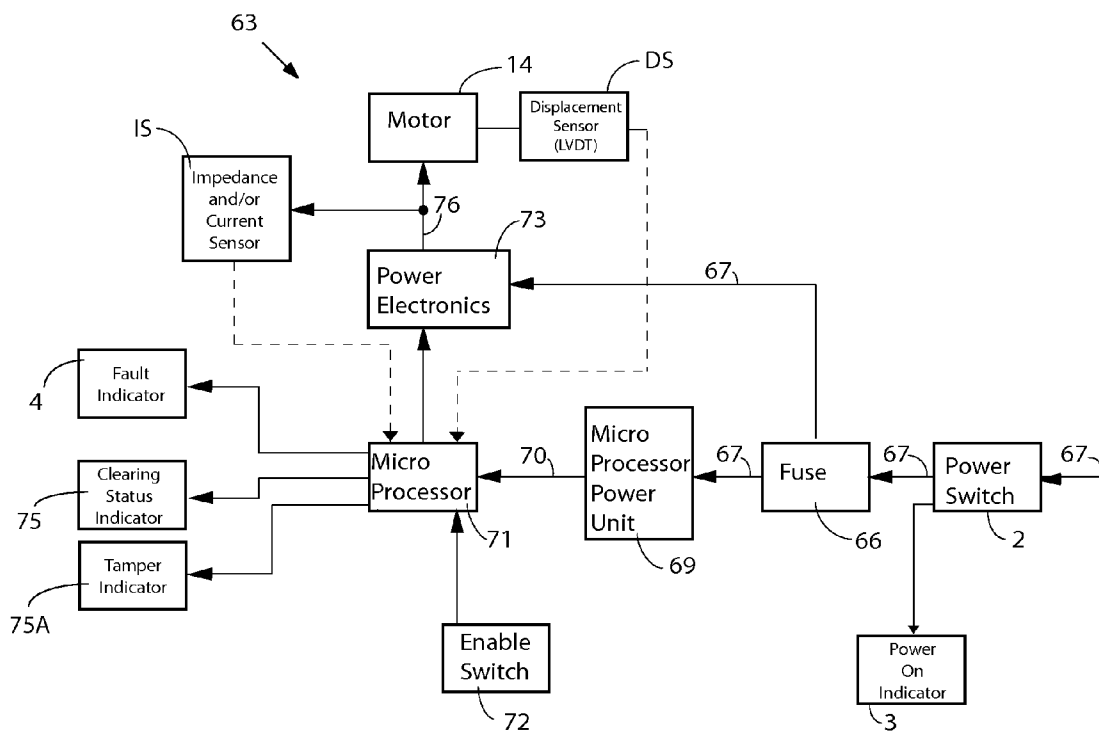


FIG. 17A

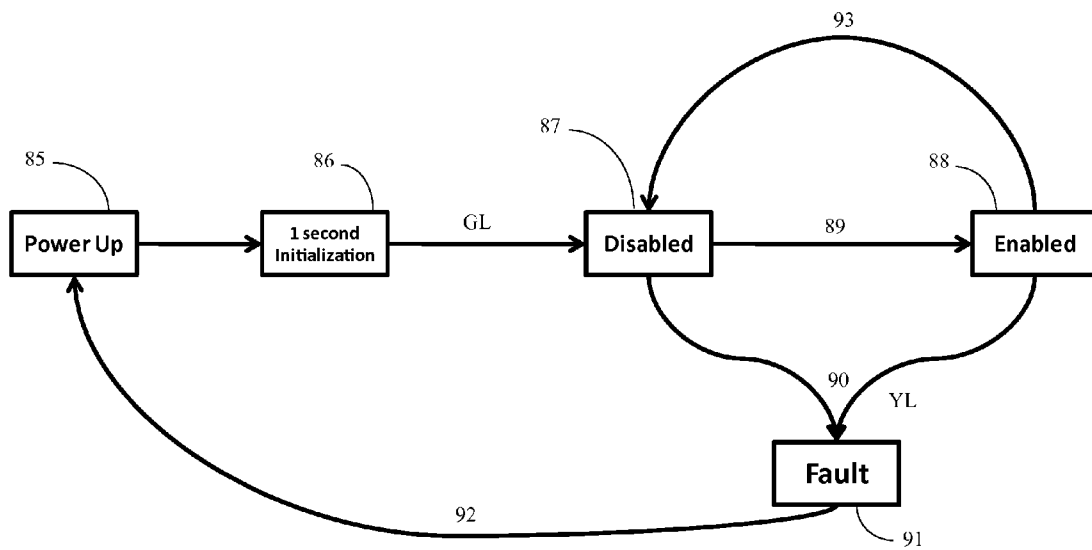


FIG. 17B

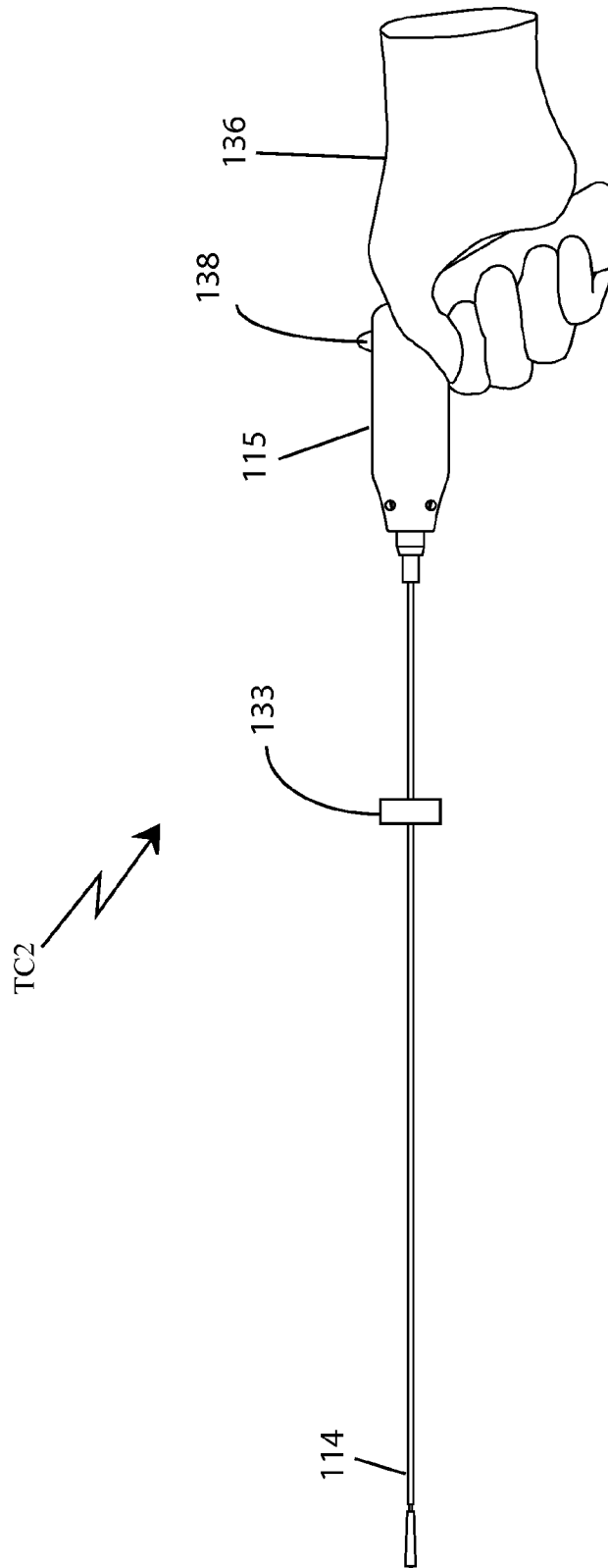


FIG. 18A

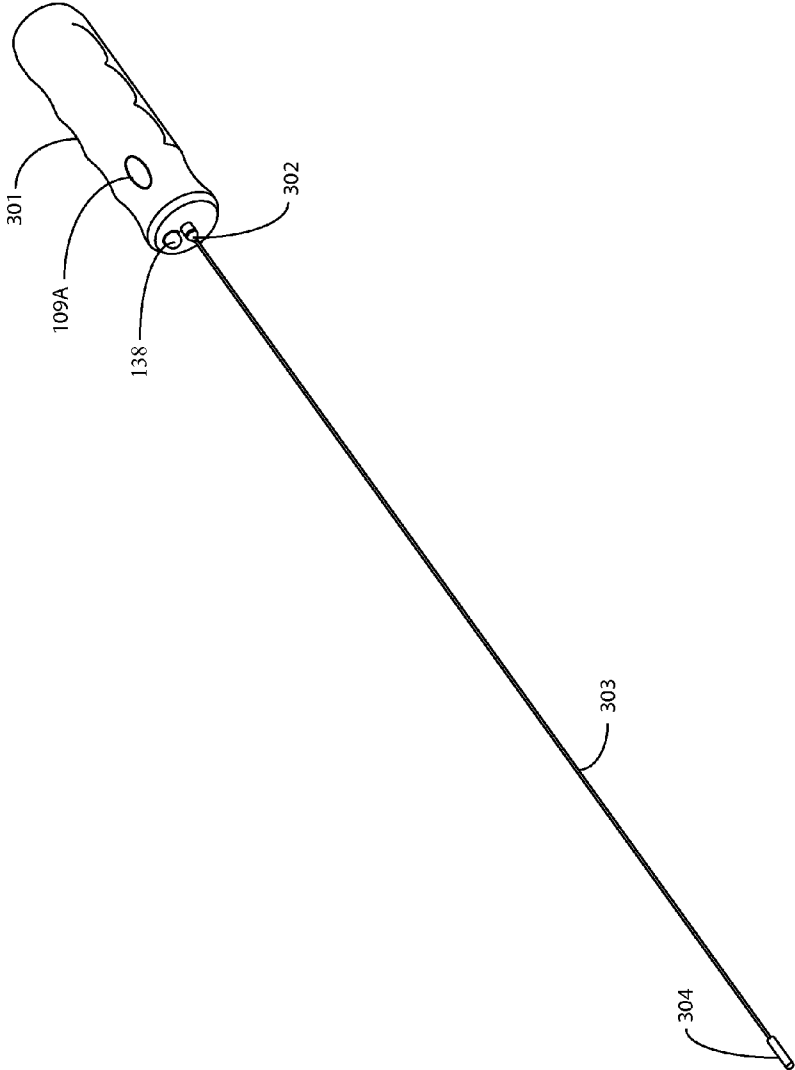


FIG. 18B

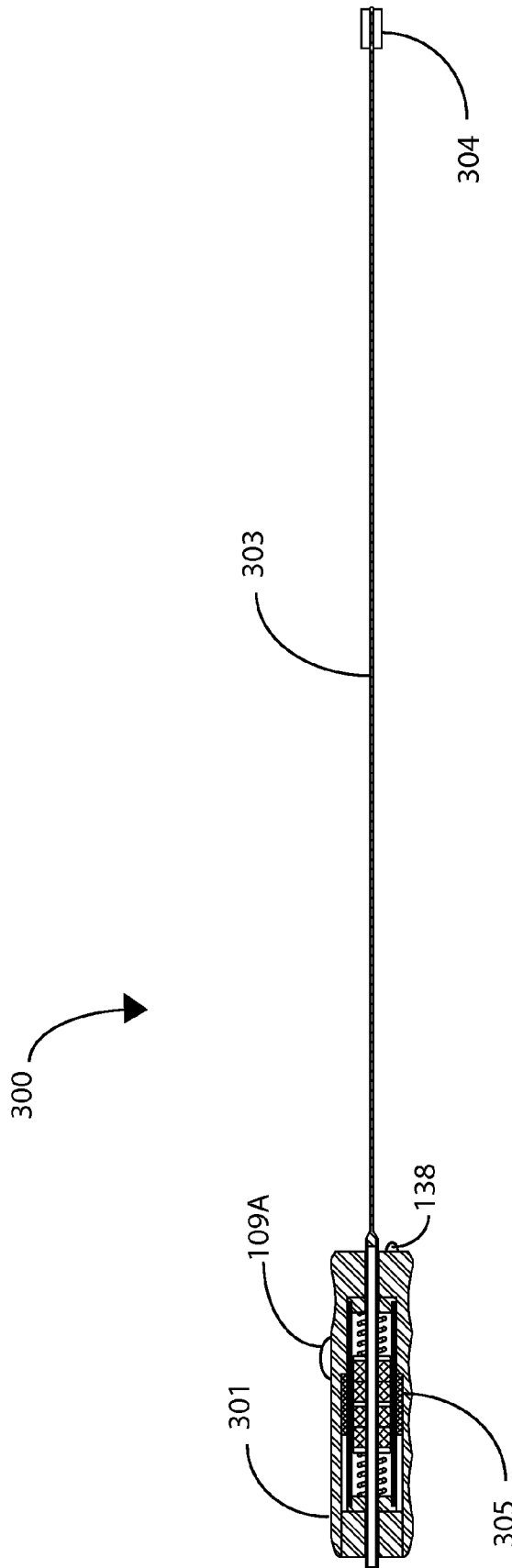


FIG. 18C

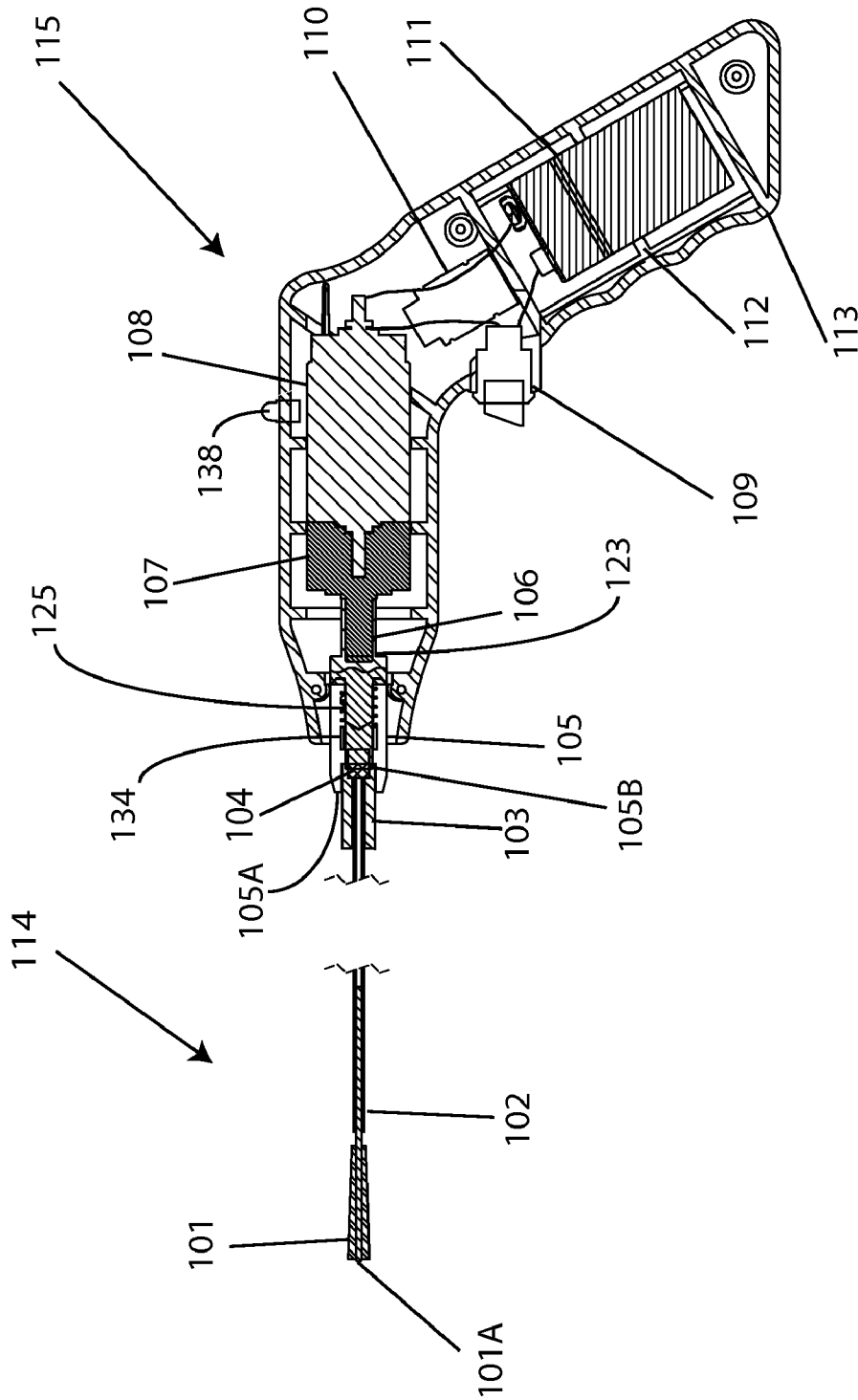


FIG. 19

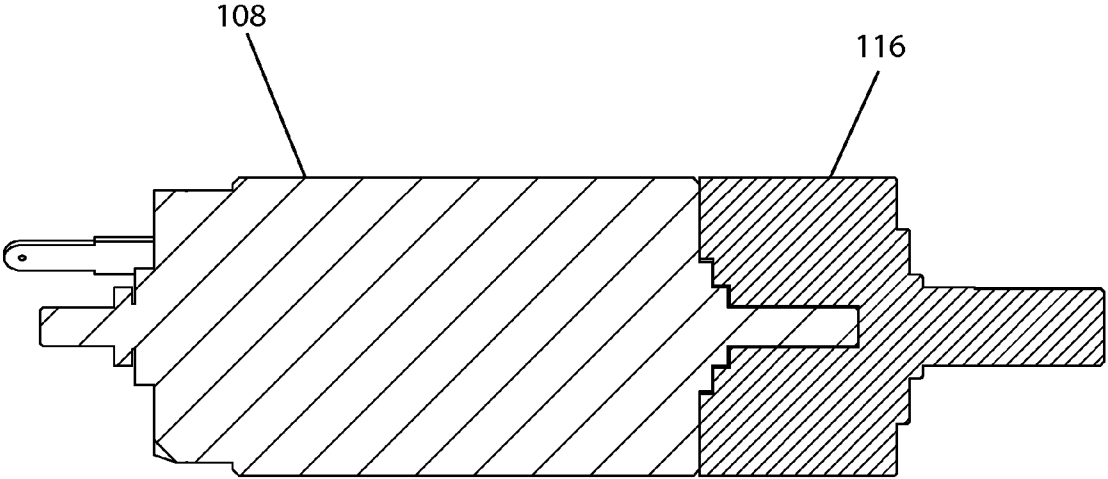


FIG. 20

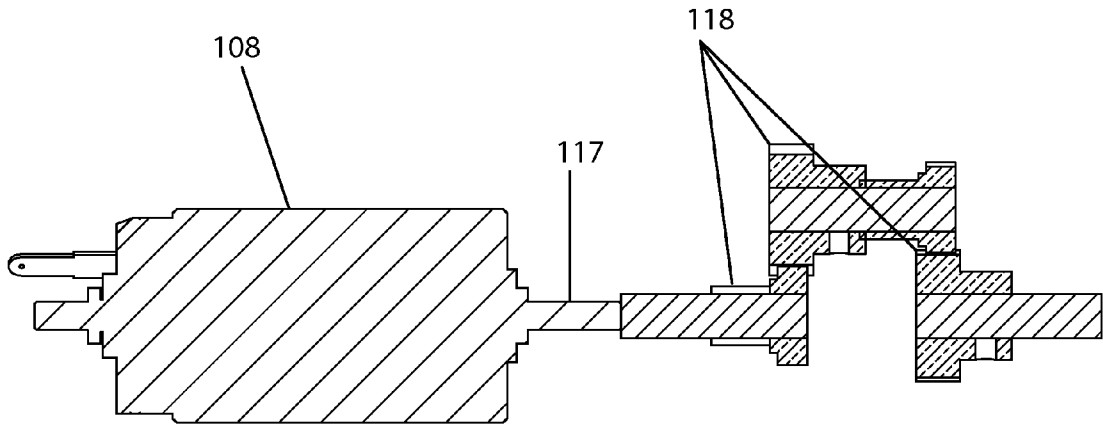


FIG. 21

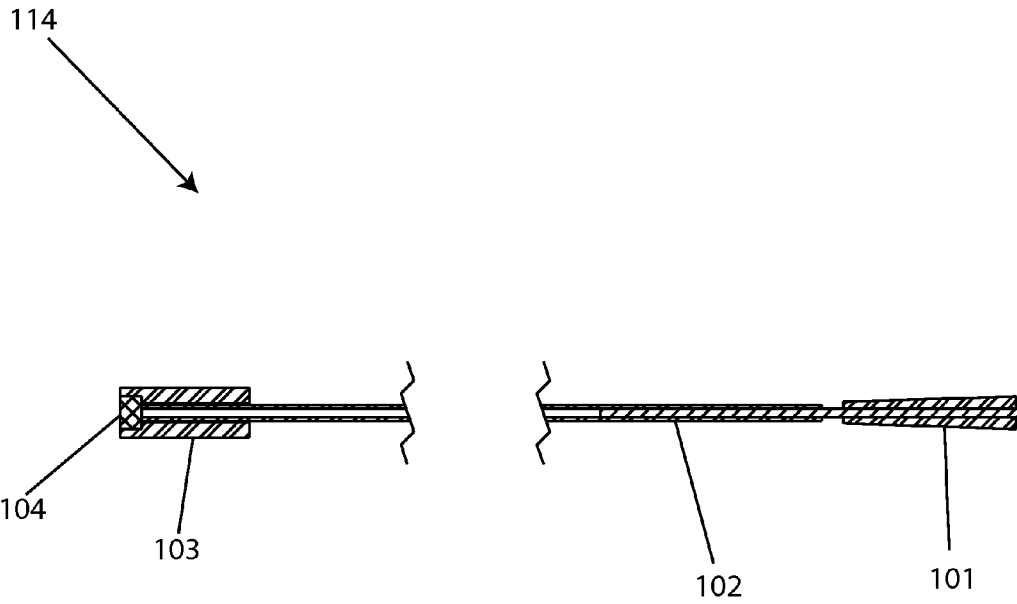


FIG. 22

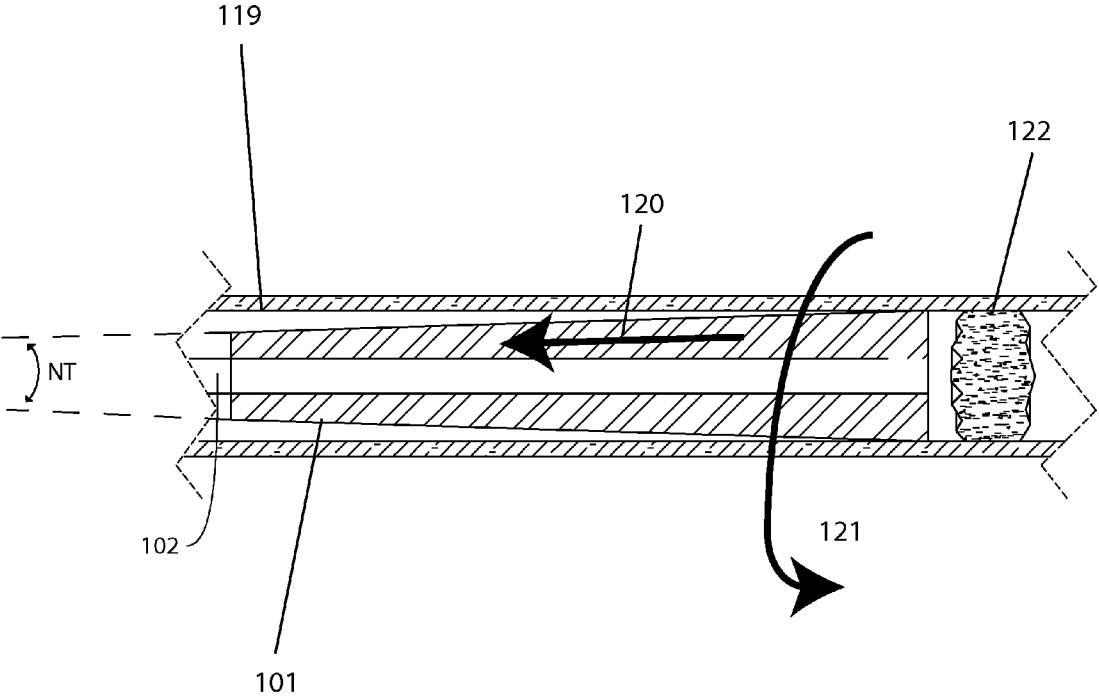


FIG. 23

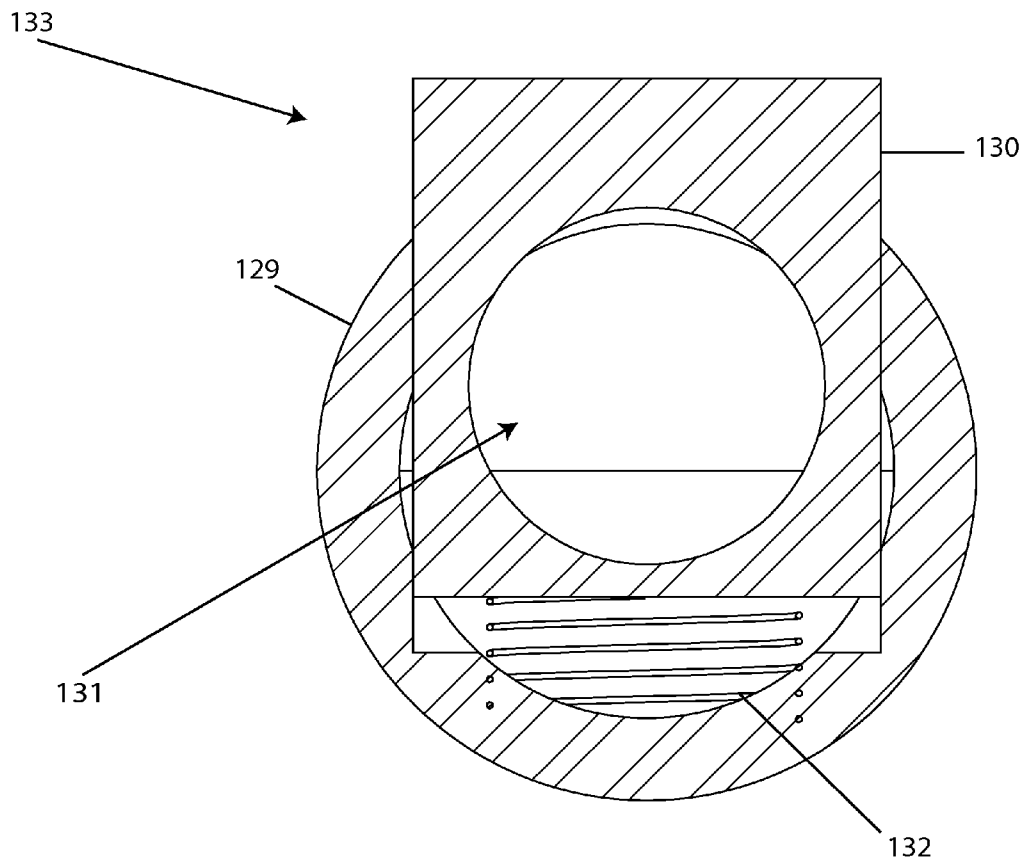


FIG. 24

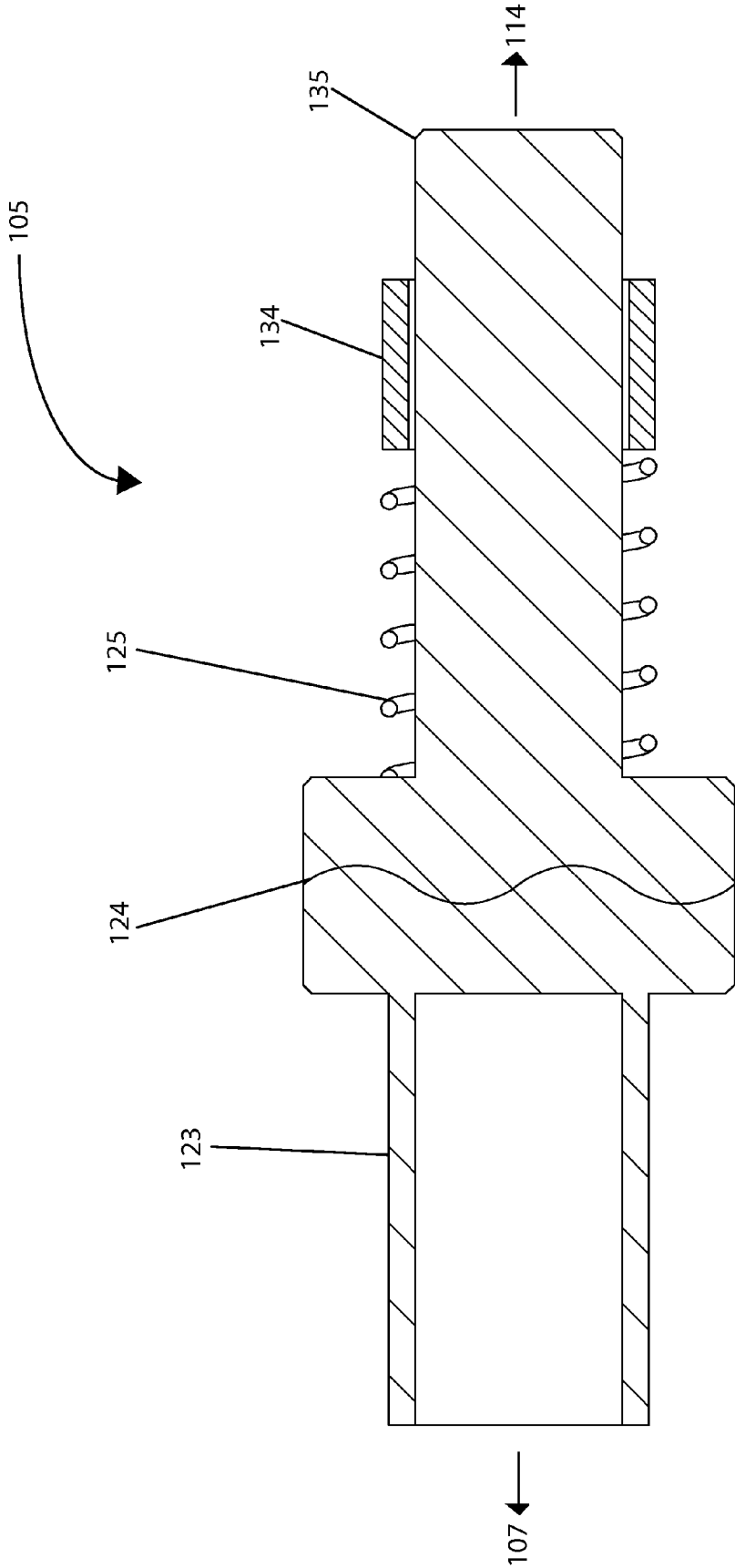


FIG. 25

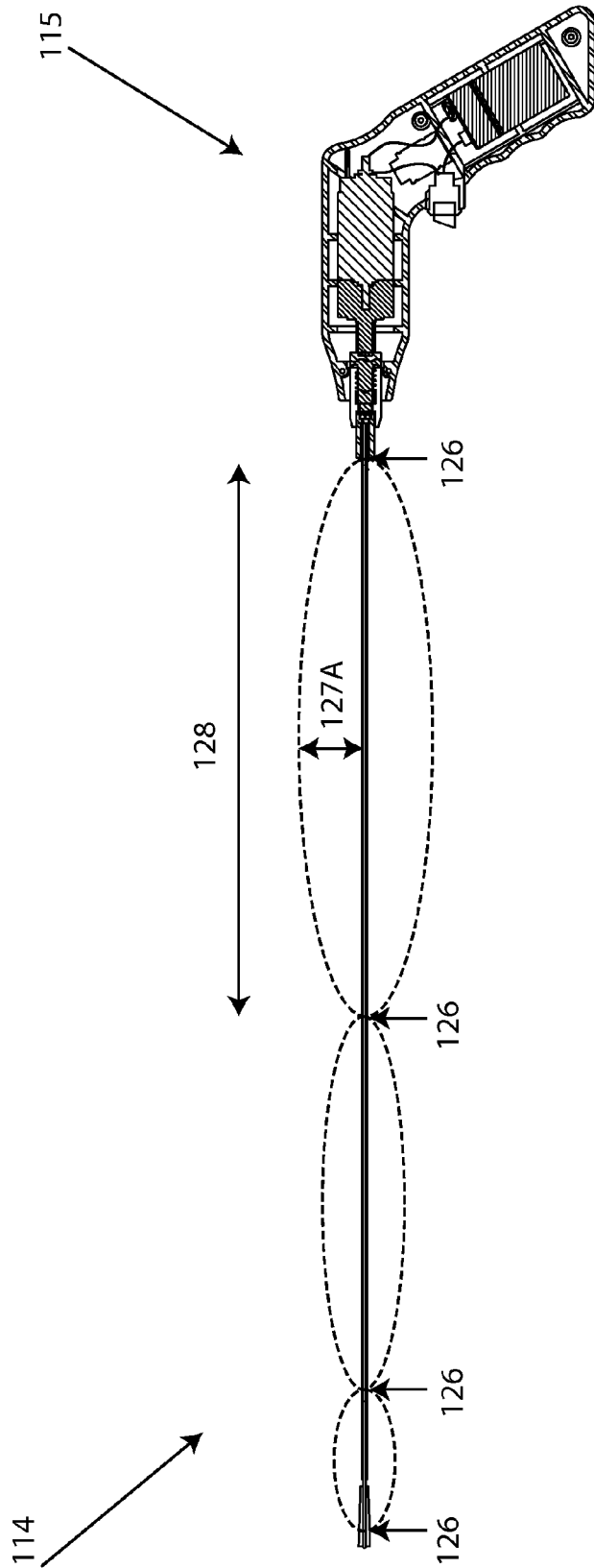


FIG. 26

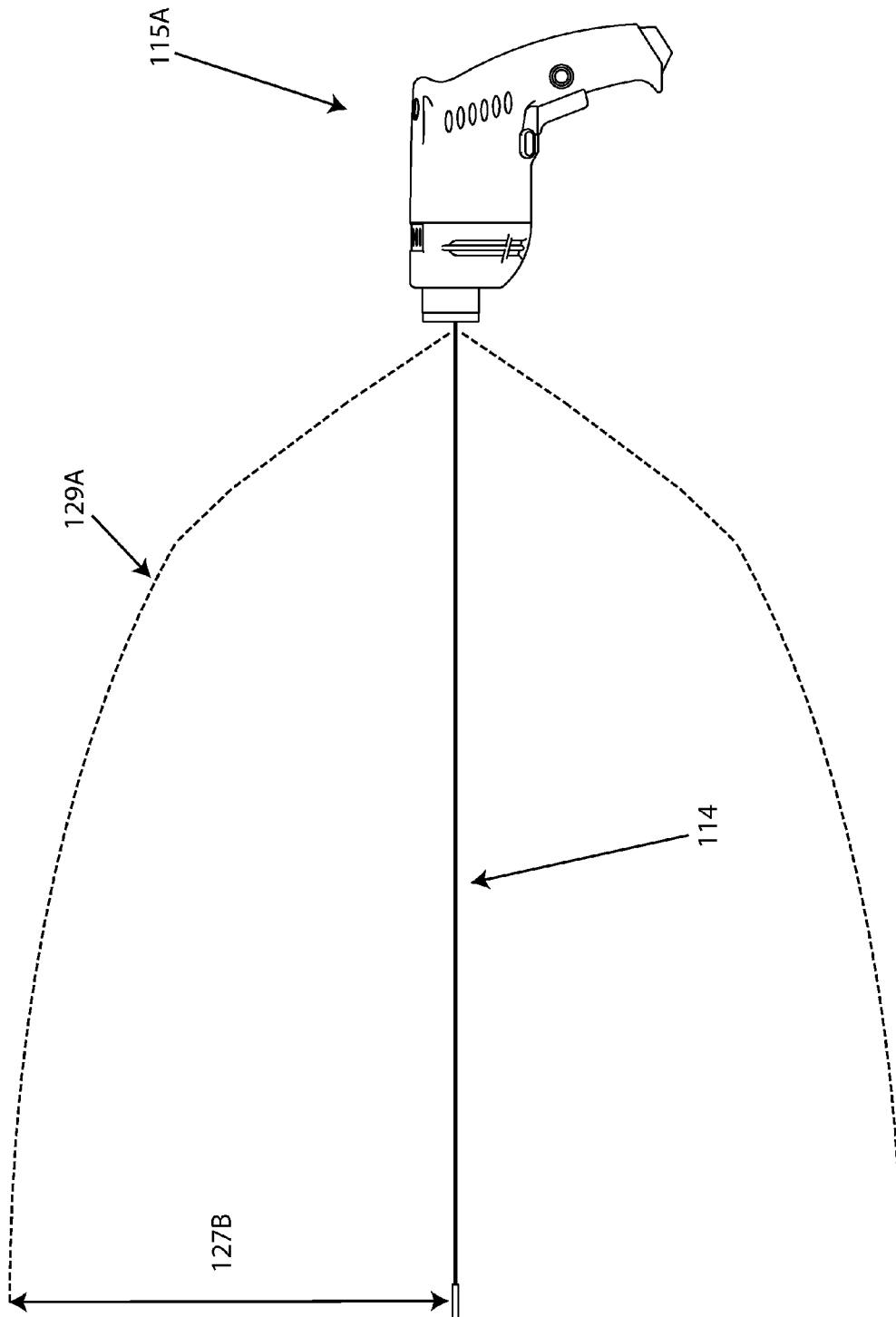


FIG. 27

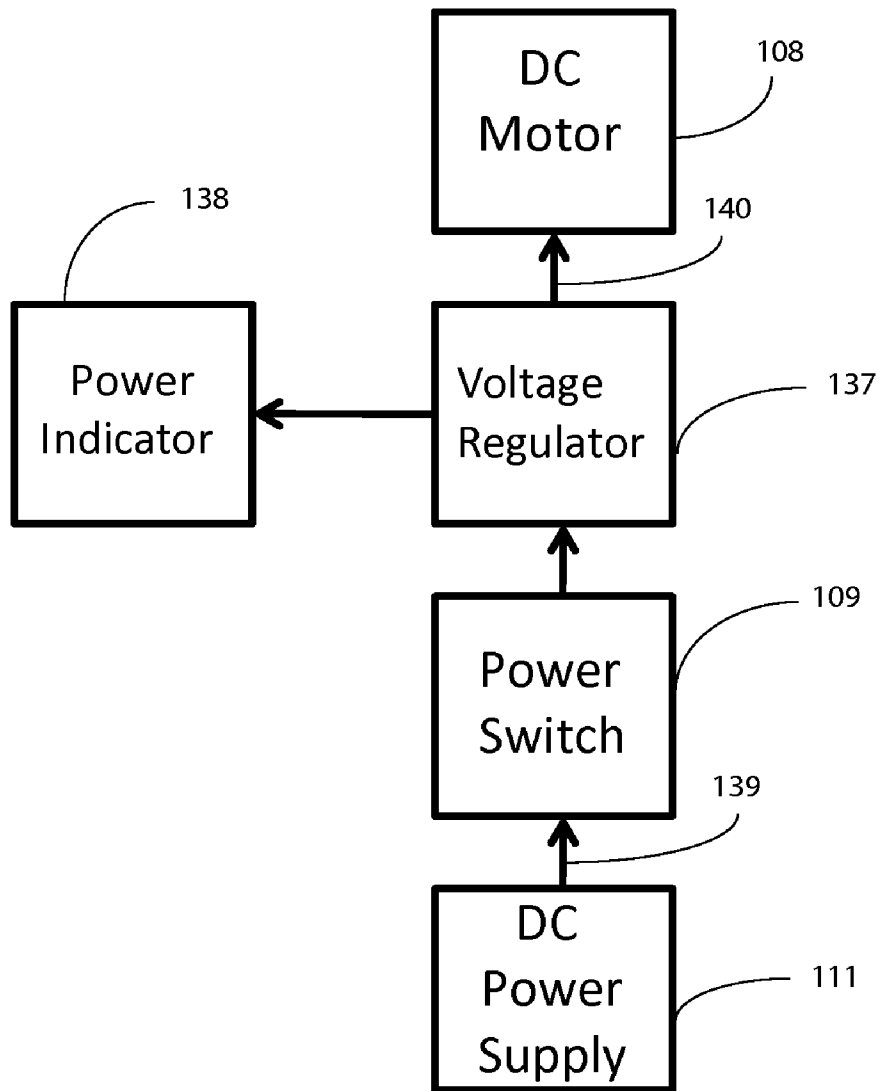


FIG. 28

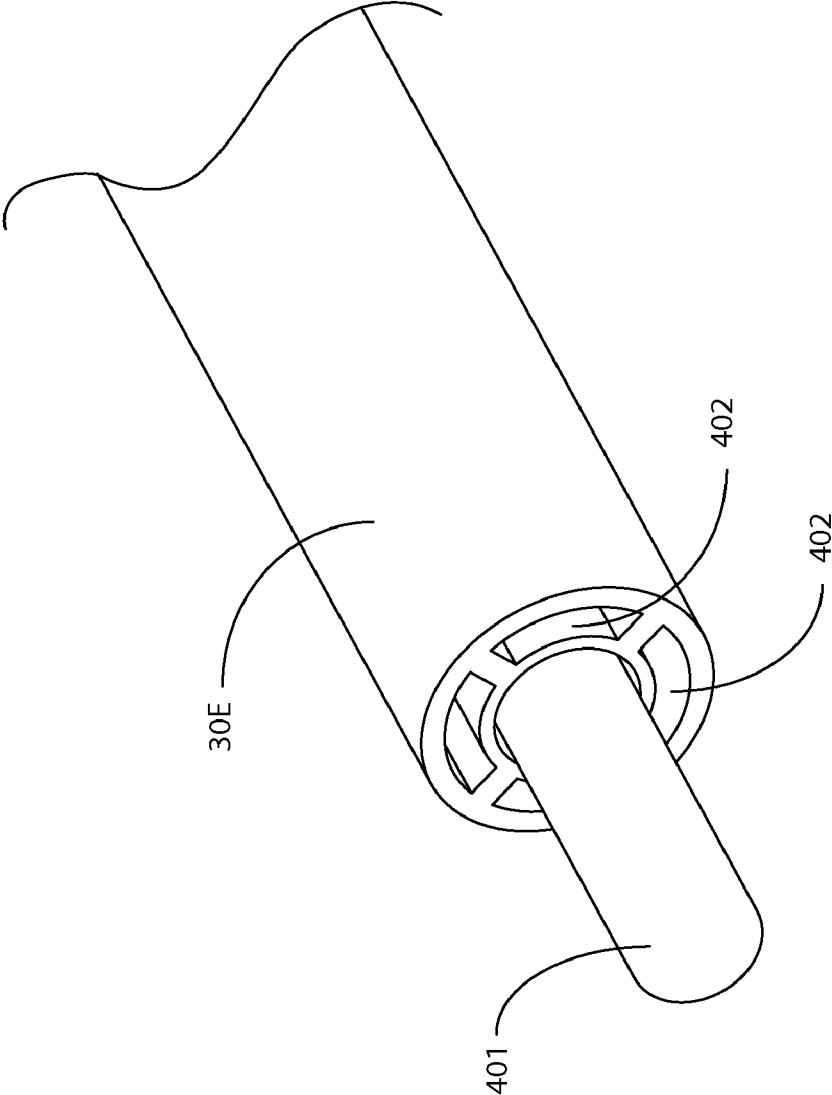


FIG. 29

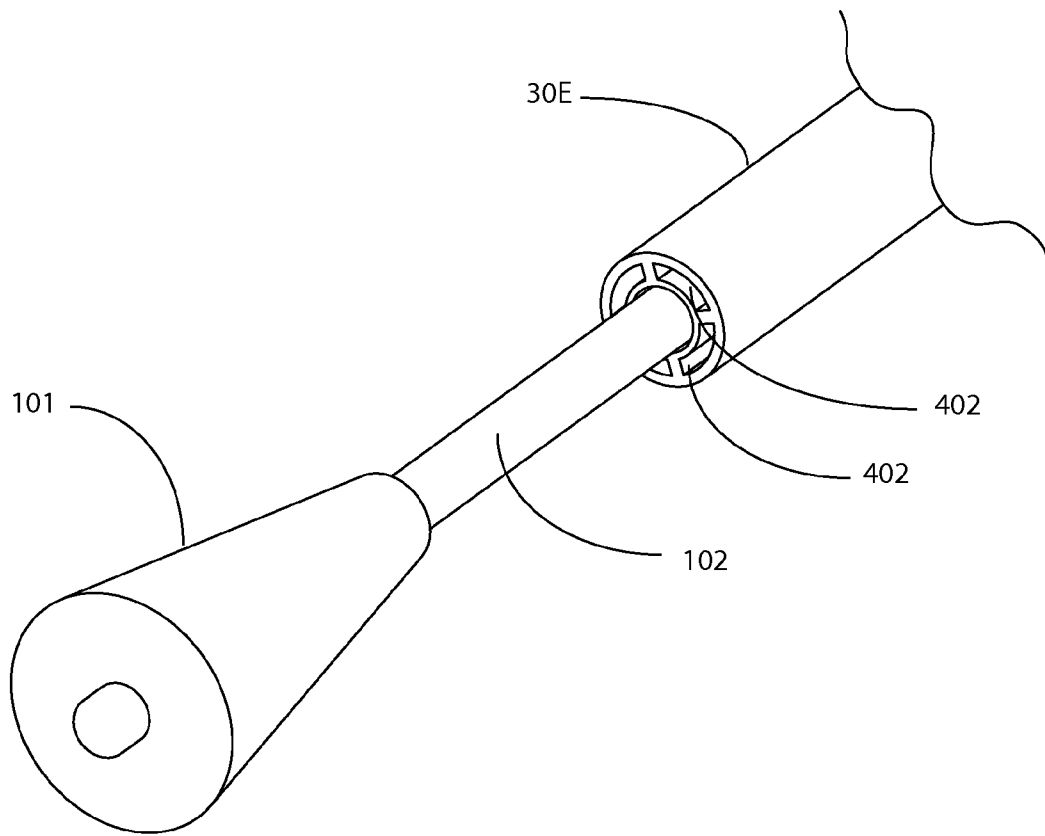


FIG. 29A

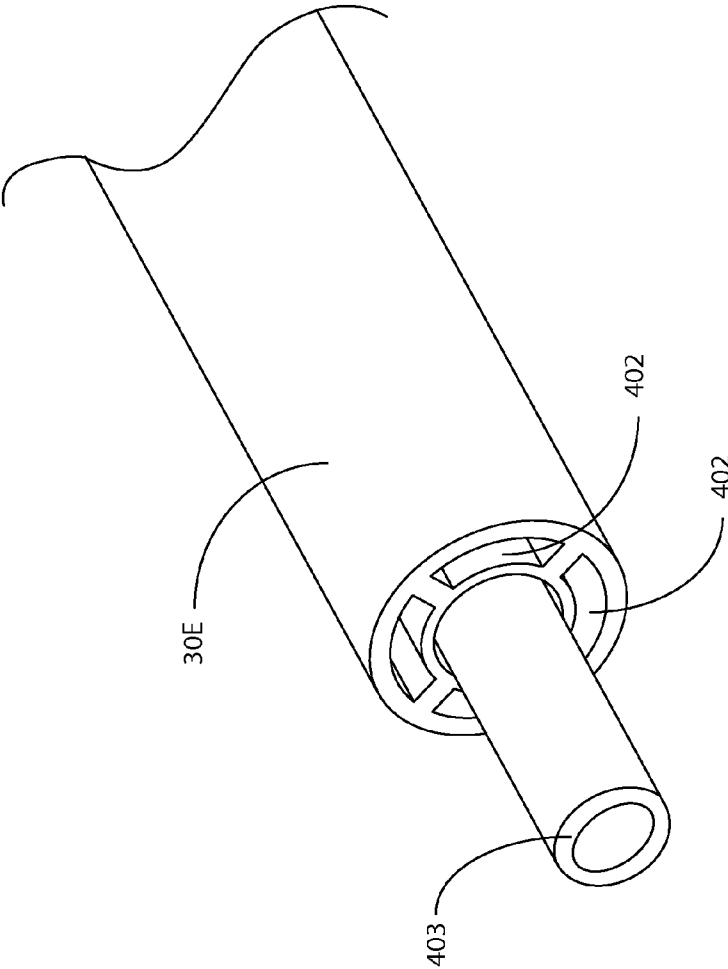


FIG. 29B

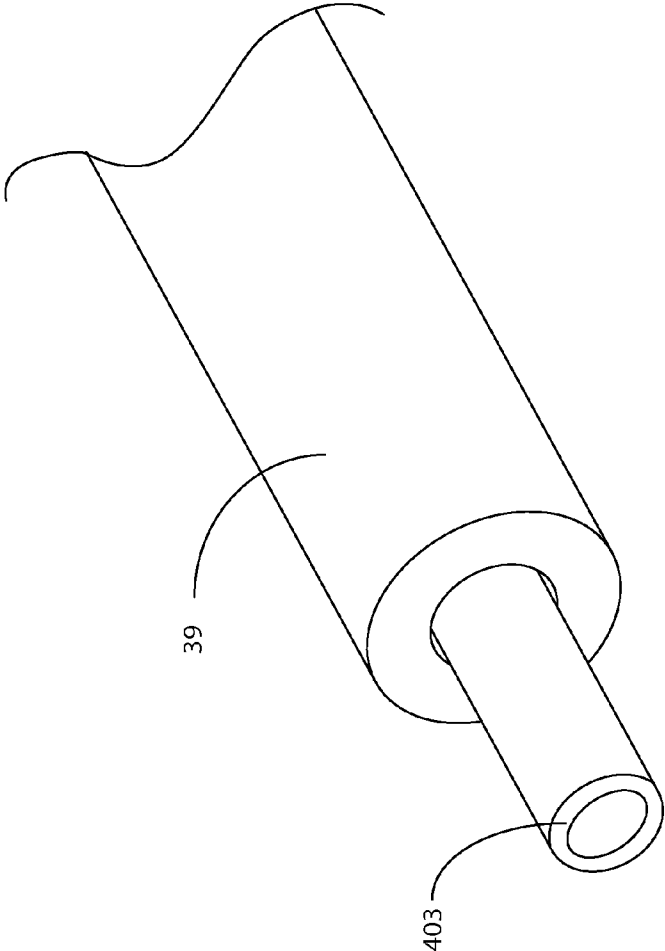


FIG. 29C

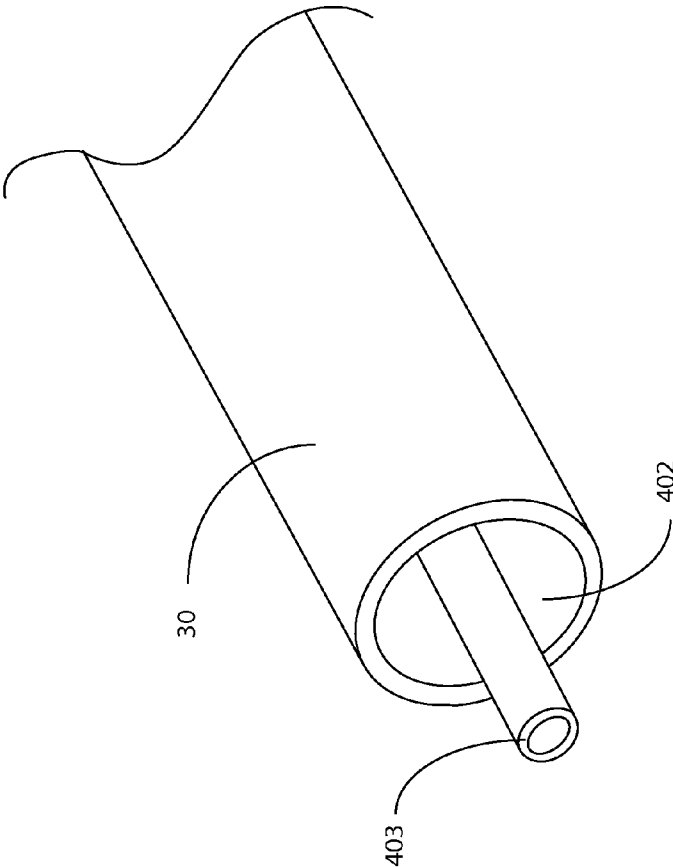


FIG. 29D

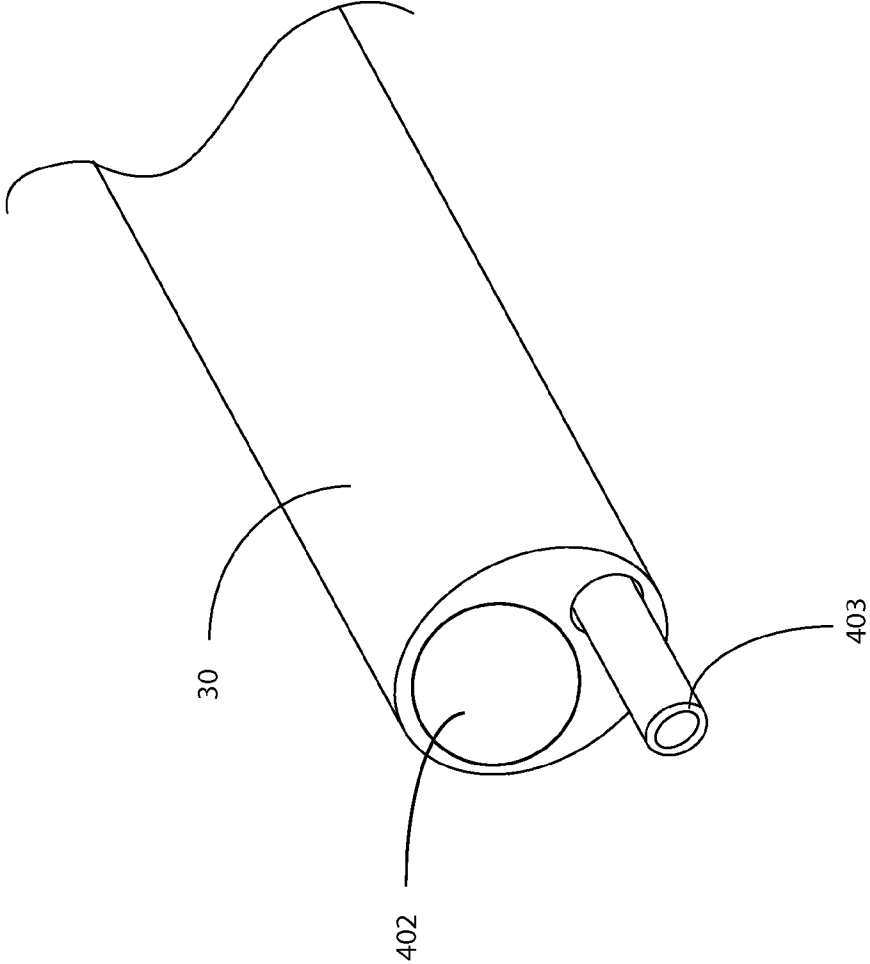


FIG. 29E

DEVICES FOR CLEARING BLOCKAGES IN IN-SITU ARTIFICIAL LUMENS

CROSS-REFERENCE TO RELATED APPLICATIONS

This continuation-in-part application claims the benefit under 35 U.S.C. §120 of U.S. application Ser. No. 12/274, 937, filed on Nov. 20, 2008 entitled FEEDING TUBE CLEANER which in turn claims the benefit under 35 U.S.C. §119(e) of U.S. Provisional Patent Application No. 60/989, 484, filed on Nov. 21, 2007 entitled FEEDING TUBE CLEANER and of U.S. Provisional Patent Application No. 61/099,737, filed on Sep. 24, 2008 entitled DEVICE FOR CLEARING BLOCKAGES IN FEEDING TUBES and all of whose entire disclosures are incorporated by reference herein.

STATEMENT REGARDING FEDERALLY SPONSORED RESEARCH OR DEVELOPMENT

This invention was produced in part using funds from the Federal government under National Science Foundation Award ID nos. IIP-0810029 and IIP-0923861. Accordingly, the Federal government may have certain rights in the invention.

BACKGROUND OF THE INVENTION

1. Field of Invention

The present invention generally pertains to cleaning or clearing devices and methods of using such devices for the in-situ clearing of artificial lumens within a living being including the in-situ clearing of feeding tubes.

2. Description of Related Art

The following is a description of the background of feeding tubes. It should be understood that the device and method of the present invention is not limited to the clearing of feeding tubes but is applicable to a range of artificial lumens such as indwelling catheters and that feeding tubes are being discussed simply by way of example.

A feeding tube is a medical device used to provide nutrition to patients who cannot obtain nutrition by swallowing. The state of being fed by a feeding tube is called enteral feeding or tube feeding. Placement may be temporary for the treatment of acute conditions or lifelong in the case of chronic disabilities. Varieties of feeding tubes are used in medical practice and are usually made of polyurethane or silicone.

A gastric feeding tube, or "G-tube", is a tube inserted through a small incision in the abdomen into the stomach and is used for long-term enteral nutrition. The most common type is the percutaneous endoscopic gastrostomy (PEG) tube. Feeding tubes may also be of the nasogastric type commonly called "NG-tube", which are introduced through the nose, down the esophagus and into the stomach in a procedure called Nasogastric intubation. PEG-tubes on the other hand are placed endoscopically: the patient is sedated, and an endoscope is passed through the mouth and esophagus into the stomach. The position of the endoscope can be visualized on the outside of the patient's abdomen because it contains a powerful light source. A needle is inserted through the abdomen, visualized within the stomach by the endoscope, and a suture passed through the needle is grasped by the endoscope and pulled up through the esophagus. The suture is then tied to the end of the PEG-tube that is to be external, and pulled back down through the esophagus, stomach, and out through the abdominal wall. The tube is kept within the stomach either

by a balloon on its tip (which can be inflated or deflated) or by a retention dome which is wider than the tract of the tube. In the case of NG-tubes, once they are passed through the patient's nostril, a clinician must be careful not to accidentally slip the end of the tube into the patient's lungs. Additionally, upon placing the NG-tube in the patient's gastric system, for example the stomach, it is common for the tubes to slip as the primary securing means is to tape the tube to the patient immediately outside the nostril. Clinicians may pass nutrients to the patient's stomach or remove fluids from the patient via the lumen or NG-tube.

Approximately 410,000 PEG-tubes and 5 million NG-tubes are placed each year in the U.S. A down-side of the life-sustaining feeding tube is that they can become clogged. Based on a 35% clogging rate, US civilian medical facilities, treat over 1.7 million NG clogs and 140k PEG clogs annually.

Numerous conditions that may necessitate enteral nutrition over long periods of time include but are not limited to traumatic injury or elderly illness such as Alzheimer's, Parkinson's, or Cancer. When long-term enteral access is needed, gastrostomy—(G), jejunostomy—(J) or gastrojejun—(GJ) tubes are often surgically inserted. J- and GJ-tubes are employed when gastric complications are present and improved nutrient uptake is necessary. Therefore, the J-tube distal end is positioned in the bowels. Reported clogging rates of GJ and J-tubes have been as high as 35% mainly due to the small bore, considerable length, and convoluted geometries of the tubes once placed. As the discussion below suggests, standard nursing protocols to clear tube occlusions are time consuming at best and are often unsuccessful. GJ- and J-tubes are especially challenging due to the curvature associated with placement.

When a patient's enteral feeding tube becomes clogged, the process of clearing it can be time-consuming and expensive, especially if the tube must be replaced. Additionally, a clog can interrupt the patient's supply of nutrients and cause him discomfort. Many nursing policies recommend flushing feeding tubes with water every four to six hours, and before and after administering medications or checking gastric residuals. Even with these policies, the rate of feeding tube occlusion is approximately 12.5%. Small-bore tubes are even more prone to clogging than are large-bore tubes, and clogging of these tubes has been shown to be a major cause of feeding downtime. A patient with an occluded tube may miss several hours of feeding and receiving nutrients before the tube is unclogged or replaced. This concern, along with patients' discomfort and the expense incurred by having to replace tubes that could not be unclogged, identifies problems to be corrected by the present invention.

Over time, feeding tubes become brittle and need to be replaced. A major cause of this is the accumulation of fungus inside the feeding tube. Standard feeding tube maintenance is to "flush" feeding tubes with water; however, this does not remove debris and fungus from the inner walls. Once a tube clogs, it is prone to reclogging.

Medications are the number one reason for tubes getting clogged. Certain medications, such as Metamucil or liquid pain reliever, build up on the inner walls of the tube and promote clogging. Other medications need to be crushed and mixed with water. If these medications are not adequately flushed or crushed finely, they will clog the tube. Older patients receive an average of 8-11 medications regularly throughout the day. Due to medical restrictions on fluid intake, or if the care-giver is rushed, an adequate flush may not occur. A clogged tube can leave an already compromised

patient without medication or nutrition for hours, or even days, and is extremely frustrating to both the patient and the caregiver.

Patients with long-term feeding tubes are generally cared for at home or in a long term nursing facility. Advancements in technology and home nursing have allowed the utilization of home enteral nutrition to dramatically increase over the last few decades. While this is certainly positive, the down side is that when a feeding tube becomes clogged such that it cannot be unclogged with conventional methods, the patient must be transported to a specialty hospital to have the tube surgically removed and replaced. For persons recovering in rural areas, this could be even more problematic as an extensive car ride—several hours—may be necessary to reach the specialty hospital. This disruption is a time consuming, expensive, and agonizing experience for the patient and family members. Numerous hours without nutrients and medication could have significant adverse effects on recovery of wounded soldiers, elderly and chronically ill patients.

One product which claims the ability to assist in restoring feeding tubes by degrading the clogged matter is the CLOG ZAPPER™ available through CORPAK® MedSystems of Wheeling, Ill. and is disclosed in part in U.S. Pat. No. 5,424,299 (Monte). This product relies on a chemical solution being injected into an enteral feeding tube to clear remnant food from the tube and decontaminate the tube. The chemical solution mixture comprises maltodextrin, cellulase, alpha-amylase, potassium sorbate, papain, ascorbic acid, disodium phosphate, sodium lauryl sulfate, disodium EDTA, and citric acid. While the solution provides some assistance in degrading the clogged matter, some patients may be allergic to at least one of these ingredients and the system for introducing the chemical solution is not always successful.

The current state of science includes three approaches to remove a clog: (1) syringe flush, (2) chemical and enzymatic treatment, and (3) mechanical devices.

Syringe Flush

The most recommended approach is to use a ‘flushing syringe’. The first step is to insert the syringe into the tube and pull back on the plunger to attempt to dislodge the clog. If not successful, warm water is placed into the tube and pressure, alternating with syringe suction, is performed. This may need to be repeated for up to 30 or more minutes. However, this may not always be done with enough efficiency or regularity and a high percentage of tubes remain clogged.

Chemical and Enzymatic Treatment

Chemical approaches to clog removal involve a nurse flushing the tube with a variety of reported substances, such as enzymes, meat tenderizer, soda, and fruit juices. More recently developed chemical approaches include using a dose of pancrelipase (Viokase®) and sodium bicarbonate mixed with water. The Clog Zapper uses a syringe filled with an unclogging powder with a variety of ingredients. Product directions state to allow the solution to set for an hour before flushing the tube. The InTRO-ReDUCER is a catheter that allows the solution to be introduced directly at the clog site, which has been reported to be more effective than introducing the solution at the external end of the feeding tube. Chemical approaches to clog removal are not effective. Patients can also be allergic to the ingredients in the chemical approaches, or adversely affected by the high sodium content.

Mechanical Devices

Mechanical devices to remove clogs are also available. Tiny brushes on wires can be used to break up the clog, but have been reported to pack the material in some clogs even more densely. The Enteral Feeding Tube DeClogger® by Bionix is a plastic, flexible rod with a spiral tip on the end. The

DeClogger can be twisted to break through or pull out obstructions. Even when successful, these approaches can take up to 30 minutes to several hours per patient, do not leave the tube walls clear, and do not progress through tortuous paths well.

What is needed is an apparatus capable of mechanically breaking up the clogged material from the sidewalls and inner portions of indwelling artificial tubes and catheters, and especially enteral feeding tubes. In addition, a regular maintenance schedule is preferred for using the apparatus to clean the walls of the tube. This regular maintenance cleans the tube walls of debris while stopping potential nucleation sites in which new clogs can grow from.

All references cited herein are incorporated herein by reference in their entireties.

BRIEF SUMMARY OF THE INVENTION

These and other features of this invention are described in, or are apparent from, the following detailed description of various exemplary embodiments of this invention.

It is hereby noted that the term “in situ” is defined as performing an act on an element while the element is being utilized for its commonly known function. For example, performing the act of clearing a clog or blockage from a feeding tube in situ refers to cleaning or clearing a clog or blockage in a feeding tube while the feeding tube is connected to the digestive system of a being, human or other.

It should be understood that it is the Applicant’s belief that where the clearing member of the embodiments disclosed herein utilizes a brush or brush function along any portion of the clearing member that makes any entry into the artificial lumen, the clearing member also cleans that interior portion of the artificial lumen.

A device is disclosed for the in situ clearing of blockages in artificial tubes (e.g., feeding tubes, including pediatric feeding tubes, PEG-tubes, NG-tubes, GJ-tubes, NJ-tubes, etc.) completely or partially disposed within a living being. The device comprises: a controller that remains outside of the living being, and wherein the controller comprises an actuator (e.g., voice coil motor; DC motor; piezoelectric actuator such as amplified piezoelectric actuators and Langevin transducers; solenoid motor; pneumatic motor, etc.) for generating repetitive motion (e.g., reciprocating, rotating, etc.); a clearing member having a first end that is releasably coupled to the actuator and having a second working end that is insertable into an opening in the artificial tube; wherein the second working end has a portion that comes into repetitive contact with a blockage in the artificial tube for clearing the blockage therein, wherein the clearing member comprises a flexible material that permits the clearing member to make repetitive contact with the blockage while the clearing member is positioned within a straight portion or within a curved portion of the artificial tube.

A method is also disclosed for the in situ clearing of blockages in artificial tubes (e.g., feeding tubes, including pediatric feeding tubes, PEG-tubes, NG-tubes, GJ-tubes, NJ-tubes, etc.) completely or partially disposed within a living being. The method comprises: coupling a first end of a releasably-securable flexible clearing member to a controller and wherein the controller remains outside of the living being; inserting a second working end of the flexible clearing member into an opening in the artificial tube; energizing the controller such that the flexible clearing member experiences repetitive motion (e.g., reciprocating, rotating, etc.) and positioning the flexible clearing member such that the second working end of the flexible clearing member comes into

5

repetitive contact with the blockage for clearing the blockage therein; and wherein the flexible clearing member clears the blockage when positioned within a straight portion or within a curved portion of the artificial tube.

BRIEF DESCRIPTION OF SEVERAL VIEWS OF THE DRAWINGS

Exemplary embodiments of this invention will be described with reference to the accompanying figures.

FIG. 1 is an isometric view of the control box and clearing stem of the present invention resting on a table;

FIG. 1A is an isometric view of the control box and clearing stem of the present invention disposed on another device support (e.g., a pole cart, bed, etc.), shown in partial, adjacent the patient;

FIG. 2 is a top plan view of another control box with the lid removed;

FIG. 2A is a cross-sectional view of the control box taken along line 2A-2A of FIG. 2;

FIG. 2B is a top plan view of an alternate embodiment of the control box of FIGS. 1 and 1A with the lid removed;

FIG. 3 is a side view of the clearing stem of the present invention;

FIG. 3A is a cross-sectional view of the clearing stem taken along line 3A-3A of FIG. 3;

FIG. 3B is a partial view of the sheath depicting both integer and periodic length markings;

FIG. 3C is a side view of an alternate clearing stem that is the preferred embodiment of the present invention;

FIG. 3D is a cross-sectional view of the alternate clearing stem of FIG. 3C taken along line 3D-3D of FIG. 3C;

FIG. 4 is top plan view shown in cross-section depicting the clearing stem inserted within an artificial lumen in a living being showing the clearing stem clearing a blockage and depicting the stem's radius of curvature;

FIG. 5A is a partial view of the clearing stem whose distal end includes a plastic clearing tip on the distal end of the wire;

FIG. 5B is a partial cross-sectional view of the clearing stem whose distal end includes an alternative hollow cylindrical clearing tip on the distal end of the wire including a tip compression spring (TCS);

FIG. 5C is a partial cross-sectional view of the clearing stem whose distal end includes an alternative clearing tip on the distal end of the wire including a gripping or chopping mechanism;

FIG. 5D is a partial view of the clearing stem whose distal end includes an alternative clearing tip on the distal end of the wire includes a welded ball;

FIG. 6 is a partial view of the clearing stem whose distal end includes a brush mounted on the wire tip;

FIG. 7 is a partial view of the clearing stem whose distal end includes a brush mounted on the distal end of the sheath;

FIG. 8 is a partial view of the clearing stem whose distal end includes a brush mounted on the distal end of the sheath with bristles swept toward the extreme distal end of the stem;

FIG. 9A is a top view of the tube depth-control collar;

FIG. 9B is a side view of the tube depth-control collar;

FIG. 9C is a cross-sectional view of the depth-control collar taken along line 9C-9C of FIG. 9A;

FIG. 9D is a partial isometric view of a fixed tube depth-control collar with the clearing stem inserted into a feeding tube;

FIG. 10 is a plan view of an exemplary voice coil motor (VCM) for use in the present invention;

FIG. 10A is a cross-sectional view of the VCM taken along line 10A-10A of FIG. 10;

6

FIG. 11 is a top plan view of another exemplary motor of the present invention with the lid removed and depicting a DC motor that drives a scotch yoke;

FIGS. 11A-11C depict a sequence of the scotch yoke operation of FIG. 11;

FIG. 12 is a top plan view of another exemplary motor of the present invention with the lid removed and depicting an amplified piezoelectric actuator (APA);

FIG. 12A is a cross-sectional view of the APA control motor taken along line 12A-12A of FIG. 12;

FIG. 12B is a cross-sectional view of Langevin transducer control motor;

FIG. 12C is a functional diagram depicting the first four overtones of clearing stem motion introduced by the Langevin transducer;

FIG. 13 is a top plan view of another exemplary motor of the present invention with the lid removed and depicting a solenoid;

FIG. 13A is a cross-sectional view of the solenoid motor taken along line 13A-13A of FIG. 13;

FIG. 14 is a top plan view of another exemplary motor of the present invention with the lid removed and depicting a pneumatic actuator;

FIG. 14A is a cross-sectional view of the control motor taken along line 14A-14A of FIG. 14;

FIG. 15 is a cross-sectional view of the magnetic pattern used in the VCM showing driving members having opposite pole directions;

FIG. 16A is a partial end view of the drive side of the control box depicting a sealing diaphragm;

FIG. 16B is a partial end view of the drive side of the control box depicting an alternative clearing stem coupling and sealing diaphragm configuration;

FIG. 16C is a partial end view of the drive side of the control box of FIG. 16 showing the clearing stem being engaged with the control box of FIG. 16B;

FIG. 17A is a block diagram of the control box electronics for the reciprocating tube clearer (TC1) configuration;

FIG. 17B is an operational flow diagram of the microprocessor of the control box electronics of FIG. 17A;

FIG. 18A depicts a hand-held version of the present invention showing the handset being gripped by the operator and including a tube depth control-collar on the clearing member;

FIG. 18B depicts an alternative hand-held version of the present invention;

FIG. 18C is a side view of the alternative hand-held version showing the hand grip in cross-section;

FIG. 19 is a cross-sectional view of the hand-held version of FIG. 18A;

FIG. 20 is a cross-sectional view of the DC motor using a planetary gear train configuration;

FIG. 21 is a cross-sectional view of the DC motor using a compound gear train configuration;

FIG. 22 is an enlarged cross-sectional view of the clearing member and its components;

FIG. 23 is an enlarged cross-sectional view of the distal end of the clearing member which uses a helical design;

FIG. 24 is an enlarged cross-sectional view of the push-button actuated tube depth-control collar;

FIG. 25 is an enlarged cross-sectional view of a torque-limiter that is designed to slip once a certain applied torque is exceeded;

FIG. 26 is a cross-sectional view of the hand-held version of the present invention depicting the multi-nodal harmonics while the clearing member is spinning;

7

FIG. 27 is a cross-sectional view of a prior-art hand-held device that generates rotatable motion depicting undesired operation with only a nodal point at the proximal end of the clearing stem;

FIG. 28 is a block diagram of the control box electronics for the rotating tube clearer (TC2) configuration;

FIG. 29 is a partial isometric view of the distal end of the sheath of the tube clearers TC1 and TC2 showing aspiration/irrigation ports;

FIG. 29A is a partial isometric view of the distal end of the sheath of the tube clearers TC1 and TC2 showing aspiration/irrigation ports;

FIG. 29B is a partial isometric view of the distal end of the sheath showing a lumen or wire that is hollow;

FIG. 29C is a partial isometric view of the clearing stem using only a hollow lumen or a wire only, without a sheath, effectively using the indwelling lumen as the sheath;

FIG. 29D is a partial isometric view of the distal end of the sheath of the tube clearers TC1 and TC2 showing a very narrow hollow wire allowing aspiration/irrigation along sides of wire; and

FIG. 29E is a partial isometric view of the distal end of the sheath of the tube clearers TC1 and TC2 showing a small sheath channel for a very narrow hollow wire and a larger channel for aspiration/irrigation.

DETAILED DESCRIPTION OF THE INVENTION

The preferred embodiments of this present invention are illustrated in FIGS. 1-29E with the numerals referring to like and corresponding parts.

The present inventions are portable devices, as well as methods for such devices, for effectively removing, moving or breaking up a clog from the internal portions of an artificial tube or catheter, enteral tube, and preferably a feeding tube, including pediatric feeding tubes. The action of removing clogs and clearing artificial tubes can also be referred to as a "maintenance action".

As will be discussed in detail later, there are basically two types of tube clearers (TC) disclosed herein, both of which are mechanical tube clearers. The first type of tube clearer TC1 includes several embodiments that generate reciprocating motion of a clearing member for removing, moving or otherwise breaking up a clog in the artificial tube. This tube clearer TC1 is preferred for use in nasogastric (NG) feeding tubes, although it should be understood that TC1 is not limited for only clearing NG feeding tubes. FIGS. 1-17B, 29, 29B, 29C, 29D and 29E are directed to TC1.

The second type of tube clearer TC2 involves the generation of rotational motion of a clearing member for removing, moving or otherwise breaking up a clog. This tube clearer TC2 is preferred for use in percutaneous endoscopic gastric (PEG) feeding tubes, although it should be understood that TC2 is not limited for only clearing PEG feeding tubes. FIGS. 5A, 5D, 18A-28, and 29A-29D are directed to TC2.

Both types of tube clearers TC1 and TC2 are unique to feeding tube clearing and overcome major obstacles in critical and long-term care medicine by clearing clogged feeding tubes quickly and efficiently. As will be discussed in detail later, the tube clearer TC1 and TC2 can remove a clog much faster (e.g., in less than 6 minutes) and at a much greater success rate than other currently-available clearing methodologies/devices, while at the same time, resulting in cleaner tube walls. Existing methodologies/devices simply do not work at all, do not clear the clogs properly, or they take a considerable time to do so.

8

In both tube clearers TC1 and TC2, an activation unit or controller remains external to the artificial tube and therefore the patient. The activation unit or controller delivers energy to a clearing stem (also referred to as a "clearing member") which is inserted into the artificial tube and whereby the clearing stem destroys the clog (e.g., clogs of food and/or ground medication, etc.) and cleans the tube walls. As a result, the activation units in these clearers TC1 and TC2 are reusable devices and the clearing stems are disposable. The clearing stems of TC1 and TC2 operate in narrow tube diameters, through several radial curves sufficient to reach, e.g., the bowel. Thus, the tube clearers TC1 and TC2 clear safely and with greater efficiency for NG-, PEG-, GJ- and NJ-tubes. Both tube clearers TC1 and TC2 require no complicated set up, e.g., no tuning is required.

Reciprocating Tube Clearer TC1

As shown in FIG. 1, the tube clearer TC1 comprises an activation unit (also referred to as the "control box" or "controller") 1 which remains external to the artificial tube 39 (see FIG. 4) being cleared, and therefore is also external to the patient (not shown). The activation unit 1 delivers energy to a clearing stem 26 which clears as it moves through the tube inner lumen 41 of the indwelling artificial tube 39, destroying the clog 40 and clearing the walls of the artificial tube 39, viz., the tube inner lumen 41 walls. Where feeding tubes are being cleared by the tube clearer TC1, the tube clearer TC1 breaks up clogs of food and ground medication in a short time (e.g., less than 6 minutes). The reusable control box 1 includes a motor which drives (actuates) the disposable clearing stem 26. The control box 1 is positioned and releasably secured onto a table, tray, or nursing cart 38, such as shown in FIG. 1. Alternatively, the control box 1 can be positioned on a pole cart 38A (see FIG. 1A), or bed rail or any other type of support that is adjacent, or which can be moved adjacent to the patient or living being.

As shown most clearly by way of example in FIG. 3A, the clearing stem 26 comprises a wire 28 running concentrically through a sheath 30. The wire 28 protrudes from the end of the sheath 30 and is actuated while the sheath 30 remains stationary and is secured to a non-moving portion of the control box 1. The motion at the wire tip 29 clears the occlusion or clog 40.

Control Box 1

As shown in FIGS. 2-2B, the control box 1 comprises a motor 14, drive electronics 10, electrical connectors, wiring, and clearing stem connectors. The control box 1 is preferably constructed of polymer, although metallic, rubber, or a combination of all three materials may be used. The preferred polymer is flame-retardant ABS plastic, although other polymers such as polyurethane, polypropylene, and nylon, but not limited to such, may be used for, among other things, their lightweight composition and structural integrity. Metals such as aluminum, titanium, steel, brass in sheet or machined form may also be used, especially where certain motor technologies (e.g., amplified piezoelectric actuators (APAs)) are used; to maintain efficiency of APAs, the non-moving portion of them needs to be effectively clamped or else too much deflection on the side that should be clamped will greatly reduce the APAs' efficiency; a metal control box provides sufficient rigidity to properly clamp. The control box 1 has a releasable securing mechanism such as rubber feet, mechanically actuated suction cup, screws, rubber stops, or magnetic feet, etc. that facilitates its use on a table or nursing cart. As such, the control box 1 remains portable but is stationary during use. The motor 14 drives a motor shaft 15 that generates the reciprocating motion.

It should be understood that FIG. 2B depicts the preferred control box 1 because it comprises a novel clearing stem-control box interface, as will be discussed in detail later with regard to FIGS. 16B-16C. FIG. 2B also depicts, by way of example only, the use of a counter balance mechanism 14A to counteract vibration caused by the reciprocation of an actuating motor 14, as will also be discussed later.

In another embodiment, the electronic circuit and componentry for example power indicator 3, fault indicator 4, enable switch 72 can be incorporated into a membrane switch such as XYMOX Technologies, Inc. Model No. 54894.

Clearing Stem/Member and Connectors

The clearing stem 26 comprises a sheath 30 which is fed into the clogged artificial tube. The preferred sheath material is polytetrafluoroethylene (PTFE) although other tube materials may also be used such as, but not limited to, nylon, polyvinyl chloride (PVC), polyurethane, polyethylene, polypropylene, fluoropolymer, Viton, Hytrel. As mentioned previously, within the sheath 30 is a wire 28, which is attached to the motor 14. The motor 14 supplies reciprocating (also referred to as "oscillating") motion to the wire 28, causing the wire 28 and its wire tip 29 to reciprocate back and forth. As can be seen most clearly in FIGS. 3-3A, the wire 28 protrudes beyond the end of the sheath 30, and into the clog 40 (FIG. 4) which causes the disruption of the clog 40. The length of the wire protrusion 28A beyond the end of the sheath 30 strongly impacts the effectiveness of the clearing. In addition, the roundness of the wire tip 29 strongly impacts the ease of insertion of the clearing stem 26 into the artificial tube 39.

The clearing stem 26 may comprise a length of 60 cm to 250 cm, but preferably 180-220 cm, and most preferably, 203 cm. In addition, the wire 28 may comprise a flexible wire most preferably stainless steel twisted wire, but could also be helical wrapped wire or a flexible stainless steel wire encased in a polymer wrapping, such as shrink wrap. The wire 28 protrudes from the end of the sheath a distance of 0 to 13 cm, but preferably 1 to 5 cm and most preferably 2.54 cm. The clearing stem 26 releasably secures to the control box 1 via a Luer clearing stem connector 6.

It should be noted that that, alternatively, the wire 28 may be hollow to enable other features such as irrigation or aspiration of the artificial lumen, as will be discussed later.

FIGS. 3-3A depict the clearing stem 26 which uses a magnetic-based and Luer lock connection to the control box motor 14, a stem stiffener 31 at a proximal end of the clearing stem 26, the amount that the wire tip 29 extends beyond the sheath 30 (referred to as the "protrusion" or "wire protrusion") 28A, a wire stop 27, and tube depth-control collar 22.

In particular, the proximal end of the clearing stem 26 comprises a clearing stem magnet 33 and a Luer clearing stem fitting 32 (FIGS. 3-3A). The control box 1 includes a Luer clearing stem connector 6 (FIGS. 2-2A) along with a motor magnetic coupler 13 which itself includes an internal magnet 12 in the coupler bore. To releasably secure the clearing stem 26 to the control box 1, the clearing stem magnet 33 is passed through the Luer clearing stem connector 6, through a diaphragm 9 and into the motor magnetic coupler 13 where the clearing stem magnet 33 and magnet 12 come into contact to form the magnetic coupling. The Luer clearing stem fitting 32 and Luer clearing stem connector 6 are then engaged to form the Luer lock configuration. Advantages to this magnetic connector include: the omission of threads (which can suffer from stripping), the avoidance of any special tools to facilitate connection, reduced occurrence of bio-contamination, and the avoidance of having to disassemble any portion of the control box 1 in order to switch clearing stems 26. The design of the mechanical components and the strength of the two

magnets 33/12 are critical to avoid detaching the clearing stem 26 when the motor 14 is reciprocating. By way of example only, the magnets 12/33 may comprise rare earth magnets (e.g., neodymium) for holding the clearing stem wire 28 to the motor shaft 15. The appropriately-sized magnets may provide from 0.5 to 3.0 lbs of holding force. The sheath 30 is held fast to the control box 1 by the Luer lock connector/receptacle combination. It should be understood that clamping of the sheath 30 needs to have a certain force to secure the sheath 30, but not crush the sheath 30. The stiffness of the sheath 30 must be adequate to preserve the inner diameter cross section during operation. This is necessary to ensure the wire 28 is not pinched by the operator and its motion impeded. The wire 28 must also be flexible enough to navigate a small radius of curvature, such as 2.54 cm radius, while maintaining operation, as can be seen in FIG. 4. In particular, FIG. 4 depicts a clog 40 blocking the tube inner lumen 41 of an artificial tube 39 and wherein the clearing stem 26 navigates a tight radius of curvature, R, and clears the clog 40 which is located past the radius of curvature R. The magnets 33/12 may be cylindrical in shape and the magnet 12 within the motor magnetic coupler 13 is recessed within the motor magnetic coupler 13 that fits over the motor shaft 15. The magnet recess 12A keeps the magnet from sliding along its surface plane and becoming detached while it is reciprocating. A sensor (magnetic or contact, not shown) may also be implemented to illuminate an indicator 75A (e.g., an LED, see FIGS. 2A and 17A) on the control box 1 to confirm that the magnetic connection is securely made. This feature also alerts the user if the connection becomes broken during use.

In an alternate embodiment, the magnet 33 (or 12) may only be located on one of the mating pieces, and a disc or cylinder of magnetic material, be located on the other.

It should be understood that this magnetic Luer lock coupling is by way of example only. It is within the broadest scope of the invention to include other types of releasably securable connector mechanisms, such as, but not limited to, threaded couplings.

As mentioned previously, the control box 1 includes a diaphragm 9 which seals the control box 1 from contamination from the outside. As can be seen most clearly in FIGS. 2-2A, the diaphragm 9 permits magnetic attachment of the clearing stem 26 so that the magnets 33/12 can make contact while at the same time sealing the box 1 such that no debris, biological or other, enters the control box 1. FIG. 16A is an enlarged partial view showing the sealing diaphragm 9 that does not interfere with motor shaft 15 motion. The diaphragm 9 prevents, among other things, the ingress of liquids into the control box 1. The diaphragm 9 may also be located externally or on the boundary of the control box 1 so that it can be cleaned more easily.

As also mentioned previously, the preferred control box 1 is that shown in FIG. 2B wherein a preferred novel clearing stem-control box interface is used. In particular, FIGS. 16B-16C depict the drive side of the control box 1 which includes a sheath attachment bracket 83, an alternate diaphragm 9A, a diaphragm sealing ring 84 (see also FIG. 2B), the motor (e.g., voice coil motor, VCM) shaft 15 along with an alternate motor magnetic coupler 13A (e.g., a magnetic coupler for a VCM). As can be seen from FIG. 16B, the alternate diaphragm 9A contains no holes or apertures through which the clearing stem 26 passes. The diaphragm sealing ring 84 secures the compliant alternate diaphragm 9A in place. To facilitate coupling the clearing stem 26 to this control box, as can be seen most clearly in FIG. 16C, the proximal end of the clearing stem 26 comprises an alternate clearing stem fitting 32A and an alternate clearing stem magnet 33A positioned

11

within an alternate clearing stem magnetic fitting 33B. In order to couple the clearing stem 26 to the control box motor 14, the alternate clearing stem magnet fitting 33B is brought into close proximity with the alternate diaphragm 9A such that the two magnets 12 and 33A are magnetically coupled and abutting through the alternate diaphragm 9A. Thus, there is no breach of the seal of the control box 1 because the alternate diaphragm 9A remains closed. Simultaneously, the alternate clearing stem fitting 32A is secured in the sheath attachment bracket 83. As a result, reciprocation of the motor shaft 15 can occur without passing through any aperture or opening in the alternate diaphragm 9A. FIGS. 1 and 1A depict a drive-end view of the clearing stem 26 coupled to the control box 1.

As can be appreciated from FIG. 3A, the wire stop 27 limits the amount of travel of the wire 28 to the right (i.e., towards the motor 14) during operation. In an alternate embodiment, as shown in FIGS. 3C and 3D, the wire stop 27 has been removed and instead an alternate wire stop 27A is used closer to the proximal end of the clearing stem 26. This alternate wire stop 27A comprises a stretchable/pliant (e.g., silicon) tube whose ends are bonded to the alternate clearing stem fitting 32A on one side and to the alternate clearing stem magnet fitting 33B on its other side. This alternate wire stop 27A supports the wire 28 that passes through it. During operation, the alternate wire stop 27A compresses and expands accordingly without interfering with wire 28 oscillation/travel. This alternate wire stop 27A is preferred because it is located externally of the artificial tube 39 and thereby avoids having a stop at the working end of the wire 28 that could interfere with operation. Thus, the alternate wire stop 27A serves to keep the wire 28 from sliding out of the sheath 30.

As shown in FIG. 3A, the wire tip 29 of the wire is rounded to allow the wire 28 to break up a clog 40 (FIG. 4), and to resist penetrating an organ (e.g., stomach or other tissue/organ, etc.) should the wire tip 29 ever make its way close to an organ. The wire protrusion 28A may also be given added flexibility by design compared to that of the rest of the wire 28, to further reduce the risk of the clearing stem wire tip 29 having enough force to penetrate an organ (e.g., the stomach) and/or to increase displacement at the wire tip 29 and facilitate clearing of the clog 40. As mentioned previously, the length of the wire protrusion 28A beyond the end of the sheath 28 and the roundness of the wire tip 29 strongly impact the ease of insertion into an artificial tube. Ideally, the wire tip 29 radius is 0.5 to 2.0 times the overall wire 28 diameter. The stiffness of the sheath 30 comprises a balance between being stiff enough to prevent the operator from clamping down on the wire 28 and stopping wire 28 motion versus being flexible enough to enter an artificial (e.g., feeding) tube 39 and to navigate curves in the tube inner lumen 41 of the artificial tube 39.

Another safety feature of the present invention TC1 is that the force generated at the end of the wire tip 29 is less than 5% of the force generated at the motor 14 and therefore, this force reduction provides a safety feature of avoiding puncturing an organ accidentally but yet providing sufficient force to break up the clog 40 and helping to clear the walls of the tube.

As mentioned previously, a stem stiffener 31 (FIGS. 3-3A) is provided at the proximal end of the clearing stem 26 which prevents the operator from over-bending the clearing stem 26 and thereby stopping the reciprocation. The stem stiffener 31 may be constructed of the same material (of a larger diameter than the wire 28 or sheath 30), may be integrated into the sheath 30 via custom extrusion, or may be constructed of a different material, such as any polymer or metal.

12

To prevent the "over-insertion" of the clearing stem 26, a tube depth-control collar 22 (FIGS. 3-3A and 9A-9C) is provided. The tube depth-control collar 22 comprises a tube depth-control collar body 24 which includes an internal spring 25. A tube depth-control collar push button 23 is provided to lock or unlock the tube depth-control collar 22. In particular, as shown most clearly in FIG. 9A, the depth control collar push button 23 has a central passageway of push button 23A and the tube depth-control collar body 24 has a central passageway of collar body 24A. A spring 25 acts to misalign these two passageways 23A/24A. Thus, to re-position the tube depth-control collar 22 along the length of the sheath 30 (not shown), the depth control collar push button 23 is depressed which momentarily relieves any clamping force on the sheath 30 and the tube depth-control collar 22 can then be moved. When the operator wishes to lock the tube depth-control collar 22 in position, he/she releases the tube depth-control collar push button 23 which results in the sheath 30 being clamped between an upper portion of collar body 24B of the tube depth-control collar body 24 and a lower portion 23B of the tube depth-control collar push button 23. The force applied by the depth-control collar to the sheath 30 needs to be compressive enough to hold the tube depth-control collar body 24 in place against the sheath 30, but not to clamp the sheath 30 onto wire 28. Sheath length markings 30A (FIG. 3B) and integer markings 30B (FIG. 3B) are provided to facilitate positioning the tube depth-control collar 22 along the length of the sheath 30 depending on the length of the artificial tube 39 being cleared. The markings 30A/integers 30B are in ascending or descending order from the distal end 30C of the sheath 30 to the proximal end 30D. Along with the stiffness of the sheath 30, the spring constant of the spring 25 comprises a balance between the force necessary to maintain the tube depth-control collar body 24 in place on the sheath 30 while avoiding the tube depth-control collar body 24 from clamping down on the wire 28 and stopping wire 28 motion.

It should be understood that it is within the broadest scope of the present invention to include fixed tube depth-control collars 22A, such as that shown in FIGS. 3C, 3D and 9D. In particular, a plurality of clearing stems 26 may be provided, each having a fixed tube depth-control collar 22A fixed at a predetermined length (e.g., 35 inches, 44 inches, etc.) along the sheath 30. FIG. 9D shows the fixed tube depth-control collar 22A abutting the proximal end of the feeding tube FT thereby preventing the sheath 30 from entering any further within the feeding tube FT. Using this embodiment, the operator selects one clearing stem 26, from a plurality of clearing stems 26, having a particular fixed tube depth-control collar 22A and clearing stem 26 length that is appropriate for the particular feeding tube FT that contains a clog that is to be cleared.

To facilitate clearing, a brush may be included on the wire tip 29 or on the distal end of the sheath 30. For example, FIG. 6 depicts a wire tip brush 35 on the end of the wire 28 whereas FIGS. 7 and 8 depict respective brushes with sheath tip brush 36 and forward swept sheath tip brush 37 on the end of the sheath 30. Therefore, as the wire protrusion 28A reciprocates, the wire tip brush 35 cleans the tube walls or when the sheath 30 is inserted into the artificial tube 39, the insertion motion causes the brush 36 or 37 to clean the tube walls, as well as facilitate the movement of the dislodged blockage and/or its pieces. In particular, the small brush (e.g., polyester, foam, or twisted in wire) on the distal end of sheath (36 or 37) or wire (35) provides more thorough clearing of tube walls. With particular regard to brush 36 or 37, mounted on the distal end of the sheath 30, the brush 36 or 37 is non-moving in this embodiment, which helps to clear excess particles from tube

walls after the wire protrusion **28A** has cleared the clog **40** and as the sheath **30** is retracted and moved out of the artificial tube **39**. The advantage of the brush **36** or **37** on the sheath **30** is that the brush **36** or **37** does not impede the wire **28** motion at all. It should be noted that the forward swept sheath tip brush **37** on the distal end of the sheath **30** shown in FIG. **8** includes bristles that are swept in the distal direction. This makes clearing effective as the forward swept sheath tip brush **37** is inserted into the tube, but also allows for a smoother retraction because the sweep-direction of the bristles reduces the resistance of the forward swept sheath tip brush **37** when the operator is removing the clearing stem **26** from the artificial tube **39**. This reduced resistance minimizes the chance of dislodging the artificial tube **39** from the patient when the clearing stem **26** is removed.

Other configurations of the clearing stem **26** include a range of wire tip **29** designs. For example, a sphere (e.g., metal or plastic) anywhere along the length of the wire protrusion **28A** may be included, such as the ball tip **34E** in FIG. **5D**. If the sphere is included at the wire tip **29**, this helps prevent the inadvertent insertion into an organ (e.g., stomach) wall, and also prevents the inadvertent retraction of the wire protrusion **28A** into the sheath **30** during use, setup or clearing illustrated in FIG. **5D**. Another alternative end may comprise a plastic end wherein a plastic tip is fused or ultrasonically welded to the wire tip **29** and which may comprise the shape of a point, helix, or radius, etc., illustrated in FIG. **5a**. In addition, these alternative tips may further comprise ridges or a pattern designed to sweep broken debris away from the clog **40** site. FIG. **5A** depicts the distal end of the wire **28** with a plastic wire tip **34**. An alternative tip design may include a spring guide wire design possibly exemplified by Lake Region Medical Paragon Pre-coat guidewires. Another alternative tip could be flexible such as a Tecoflex® tip which causes the tip to slide across contacted tissue rather than puncturing tissue, thus providing an additional safety feature.

FIG. **5B** depicts another alternative end which may comprise a small spring mechanism which provides increased displacement and protection against an over-insertion puncture. In particular, a plastic or metal alternate tubing tip **34A** is positioned over the distal end of the wire **28**. The rear end of the alternate tubing tip **34A** is secured to one end of a tip compression spring TCS that is slid onto the wire **28**. A fixed member **34B** is secured to the wire **28** and to the other end of the tip compression spring TCS. Thus, the alternate tubing tip **34A** acts as a further protection against accidental contact with soft tissue, since the alternate tubing tip **34A** can only be retracted when it encounters a solid object, e.g., a clog, and whereby the wire tip **29** is then exposed to the solid object. Once the clog is cleared, the alternate tubing tip **34A** springs back in position ahead of the wire tip **29** to shield it from contact with bodily tissue or organs. Moreover, the wire tip **29** may also comprise a small gripping mechanism wherein the wire tip **29** contains a small cable-actuated gripping mechanism to dislodge clogs **40** or retrieve samples of clog material. In particular, FIG. **5C** depicts gripping/chopping mechanism **34C** that are hinged or pivoted at pivot point **34D**. By actuating a control member (not shown, e.g., a cable, rod, electromechanical motor, piezoelectric motor etc.), the gripping/chopping mechanism **34C** can be closed around a clog specimen or used to tear away the clog material to dislodge clogs or retrieve a sample of the clog material.

An alternative design to the wire **28** is the provision of a flexible portion of wire **28** located between the end of the sheath **30** and the wire tip **29**. Thus, the wire protrusion **28A** may comprise a material that is more flexible than the remaining part of the wire **28** that couples to the motor shaft **15**.

Control Box Motor for TC1

As mentioned previously, the motor **14** drives the wire **28**, creating linear displacement. The back and forth displacement of the wire **28** allows it to break up and clear clogs **40** in artificial tubes (e.g., enteral feeding tubes and especially NG feeding tubes), while simultaneously cleaning debris from the tube walls. The wire tip **29** of the wire **28** has a linear displacement, preferably, in the range of 0.25 to 25 mm, more preferably 2-10 mm from the distal end of the sheath **30**. The frequency of operation of the motor shaft **15** preferably varies from 10 to 100 Hz but more preferably in the 15-40 Hz range. The motor **14** has a range of displacement preferably from 1-40 mm and more preferably in the range of 10-30 mm. The motor blocking force (i.e., the maximum force output) has a preferable range of 2-25N and more preferably 6-14N.

The reciprocating motion of the clearing stem **26** of the present invention TC1 can be achieved using a variety of motor technologies, such as, but not limited to, voice coil motors (VCMs) as illustrated for the motor **14** (FIGS. **2-2B**, **10-10A** and **15**), DC motors **49** (FIG. **11**, **11A-11C**), piezoelectric transducers, including amplified piezoelectric actuator motors **59** (APA, such as those disclosed in U.S. Pat. No. 6,465,936 (Knowles, et. al), whose entire disclosure is incorporated by reference herein) (FIGS. **12-12A**), piezoelectric actuators, active polymer compound actuators, solenoid motors **55** (FIGS. **13-13A**), pneumatic motors **42** (FIGS. **14-14A**), magnetostrictive transducers, electrostrictive transducers, etc.

As shown in FIGS. **2-2A**, **10-10A**, and **15** the motor **14** may comprise a voice coil motor (VCM) having a VCM body **16** mounted within end bearings **18**, a displaceable motor shaft **15**, dampers or spring **19**, and magnets **20** mounted to the motor shaft **15**, with pole pieces **21A**, **21B** and **21C** (FIGS. **2A**, **10A** and **15**) located at the ends and within the center of the magnets **20**. Coil windings **17** are wound around the VCM body **16** and thus do not interfere with VCM motor shaft **15** displacement. Motor mounts **7** and motor mount dampers **8** secure the motor **14** within the control box **1** while avoiding direct coupling against the bottom surface of the control box **1**. A motor printed circuit board (PCB) **11** distributes the current commands from the electronics **10** to the coil windings **17** through wires **53**. When an electric current is applied through the coil windings **17**, a magnetic field, due to Ampere's Law, is produced inside the coil windings. The non-uniform magnetic field at the ends exerts a force on the permanent magnets **20**. Alternating the current alternates the direction of the magnetic field gradients and results in a reciprocating motion of the motor shaft **15** with respect to the VCM body **16**. The magnitude of the force is determined by the magnetic flux density, which is proportional to the number of turns per length of the coil, current magnitude, cross-sectional area of the coil, as well as the strength of the permanent magnets **20**. The springs **19** absorb the energy associated with abrupt changes in the direction of the inertial force of the magnets **20** and VCM body **16** when actuated, resulting in a lowering of vibration and increasing the tube clearer TC1 usability and efficiency.

By way of example only, the spring constant of the springs **19** can range from 0.5-5 lb/in, and more preferably 1.5-2.5 lb/in.

A soft stop **SS** may be installed at the free end of the VCM motor shaft **15** because the shaft tends to drift off center during use.

A further variation of the use of a plurality of magnets is to arrange the plurality of magnets into two "driving members" disposed between the pole pieces **21A-21C**, mentioned previously. Pole pieces **21A-21C** are typically ferromagnetic and

15

are preferably stainless steel. As shown most clearly in FIG. 15, the south poles of the first magnetic driving member 20N and the south poles of the second magnetic driving member 20S are fixedly secured to the opposing faces of the pole piece 21B in order to provide a zone of maximum magnetic flux density which extends radially outwardly from the central portion of the pole piece 21B, similar to the configuration disclosed in U.S. Pat. No. 4,363,980 (Peterson) whose entire disclosure is incorporated by reference herein. Alternatively, each magnetic driving member 20N and 20S may be replaced with a single elongated permanent magnet, rather than using a plurality of magnet elements as shown in FIG. 15. In either case, the driving members 20N and 20S have opposite pole directions.

It is within the broadest scope of the present invention that the relative positions of the coil windings 17 and the magnets 20 are reversed (not shown), i.e., the coil windings 17 are wound directly around the motor shaft 15 and the magnets 20 are positioned around the VCM body 16 and thus do not interfere with the motor shaft's 15 reciprocation.

Alternatively, a dual coil motor or actuator (also not shown) is also within the broadest scope of the present invention. In particular, instead of using magnets 20, two coil windings are used wherein one coil is wound directly around the motor shaft 15 and a second or outer coil is wound around the first or inner coil but without interfering with shaft displacement. Each coil is supplied with respective alternating current sources which generate respective electromagnetic fields that also generate a reciprocating motion of the motor shaft 15. The inner coil may conduct direct current DC while the outer coil conducts alternating current AC. Alternatively, the inner coil may conduct alternating current AC while the outer coil conducts direct current DC, or both the inner coil and the outer coil may conduct alternating current AC.

Moreover, to reduce vibration caused by the oscillating motion of the motor shaft 15, a secondary VCM or counter balance mechanism 14A of similar size (also referred to as a "countermass" or "counterbalance") may be included and driven at an opposite phase (e.g., 180° phase lag) for cancelling vibration caused by the motor 14. See FIG. 2B. Thus, when the tube clearer TC1 is operated such that the first VCM is activated to cause the motor shaft 15 to move, a first momentum vector is produced. The second VCM is operated such that it creates a second momentum vector equal in magnitude but opposite in direction to the first momentum vector, such that the net sum of the first and second momentum vectors is minimized and preferably equal to zero. In particular, to maximize vibration reduction, the moving parts (shaft, magnets, pole pieces, attachments, etc.) of the counter balance mechanism 14A should have a moving mass and velocity (frequency and displacement) equal to that of the moving parts of the actuating motor 14. This is based on the principle of Conservation of Momentum. The sine waves that actuate both VCMs must have a 180 degree phase lag between them. This causes their forces to be opposite and (ideally) equal, cancelling each other out. As such, operation of the tube clearer TC1 does not cause "chatter" and therefore there is no irritation to the operator or patient.

DC Motor 49

The motor may also comprise DC or DC brushless motor 49 for creating reciprocating displacement via a scotch yoke SY or similar mechanism. FIG. 11 depicts the control box 1 using a DC motor 49 and scotch yoke SY as the actuating mechanism. No signal generating electronics are needed for this application since the DC motor 49 is simply turned on to cause a rotating crank CR to drive the scotch yoke slider 50 and the scotch yoke shaft 52 in reciprocating motion. The

16

adapter 51 transmits the scotch yoke SY motion to the scotch yoke shaft 52. FIGS. 11A-11C show three still frames as an example of scotch yoke SY motion. FIG. 11A and FIG. 11B show Scotch yoke forward displacement direction 50A and FIG. 11C shows Scotch yoke rearward displacement direction 50B are moving in a reciprocating motion.

APA Motor 59

An amplified piezoelectric actuator (APA) 60 creates reciprocating displacement in the lower range, preferably (0.1 to 2.0 mm), anchored to the control box 1. One or more APA motors 59 can be used in series, as this increases displacement. FIGS. 12-12A depict the control box 1 with an APA as the actuating mechanism. In particular, the APA actuator 60 is mounted to the control box via an actuator mount 61 which is indirectly coupled to the control box 1 bottom via motor mount damper 8. An actuator shaft 62 conveys the reciprocating motion, from APA actuator 60 expansion and contraction, to the clearing stem (not shown) via the magnetic coupling discussed earlier for the other embodiments.

Langevin Transducer 77

A Langevin transducer 77 can be used for the motor 14. As shown in FIG. 12B, the Langevin transducer comprises a plurality of piezoelectric elements 78 are arranged to cause a horn 81 to vibrate to form the reciprocating motion. The horn 81 is secured to an actuator mount 61 using a pre-stress bolt 79. The Langevin transducer 77 includes a tail mass 80 for bolt-clamping the Langevin transducer 77 to the actuator mount 61. The forward end of the horn 81 is tapered such that a distal end of the horn passes through the control box alternate diaphragm 9A. A clearing stem attachment 82 is provided to receive/mate with the clearing stem 26 as discussed previously. A power source (not shown) that provides the proper activation energy is coupled through the power plug 5 and via electronic control wires 53.

It should be noted that activation of the Langevin transducer 77 creates reciprocating motion with the introduction of several overtones (viz., first-fourth overtones), shown in FIG. 12C. As part of the design of the present invention, the lateral displacement caused by these overtones is kept to a minimum. In particular, the piezoelectric elements 78 (e.g., a plurality of piezoelectric ceramic discs) are held in compression between the tail mass 80 and horn 81; and the pre-stress bolt 79 passing from a proximal end of the tail mass 80 and threading into the horn 81. Vibratory motion is caused by the activation of the piezoelectric elements 78 upon being exposed to an alternating electric field such as from an AC electrical current applied to electrical contacts (not shown) formed on opposing sides of each of the piezoelectric elements 78. The vibratory motion is translated as a standing harmonic wave spanning longitudinally across the horn 81 and to the clearing stem (not shown). Therefore, when operated at ultrasonic frequencies, the Langevin transducer 77 translates the ultrasonic energy as a reciprocating vibration to the clearing stem 26, and produces a standing wave within the flexible member. The horn 81 and tail mass 80 are made of a metal such as titanium, stainless steel or, preferably, aluminum. The pre-stress bolt 79 is generally of stainless steel, but not limited thereto.

Solenoid Motor 55

The solenoid motor 55 shown in FIGS. 13-13A mounted in the control box 1 operates in a very similar manner as does the motor 14, discussed previously. A return spring 58 is required with the solenoid 56 since it has one-way actuation. In particular, the electronics 10 are configured to pulse the solenoid 56 such that during the pulse, the solenoid shaft 57 is driven to the left in FIGS. 13-13A and when the pulse is terminated,

the return spring 58 restores the solenoid shaft 57 to the right. This action is repeated at the frequencies discussed previously.

Pneumatic Motor 42

FIGS. 14-14A depict a pneumatic motor 42 for creating the reciprocating motion. In particular, the pneumatic motor shaft 44 is driven by the pneumatic motor 42 which receives pneumatic pulses from a pneumatic pulse generator (not shown) via an air supply inlet 54 on the control box 1 and through internal tubing 47. The pneumatic motor 42 is positioned within a pneumatic motor housing 43 which includes a pneumatic motor diaphragm 46 for distributing the pneumatic pulse evenly to the pneumatic motor shaft 44, thereby maintaining its alignment, while at the same time providing a tightly-sealed motor configuration. The pneumatic pulse causes the pneumatic motor shaft 44 to be driven to the left while compressing a return spring 58. Once the pneumatic pulse is terminated, the return spring 58 restores the pneumatic motor shaft 44 to the right. This action is repeated at the frequencies discussed previously.

Electronics

FIG. 17A provides a block diagram of the electronic system 63 contained within the electronics 10. A microprocessor (e.g., MSP430F2618TPMR) controls the power electronics 73 to the motor 14. Although not shown, a power supply (e.g., an Autodyne UL medically-approved power supply AMP6301-08) converts the 120VAC from the wall outlet to 24VDC. A microprocessor power unit MPU 69 (e.g., a voltage regulator circuit, such as the LM317/LM337) reduces the incoming (e.g., +24VDC) power 67 to a lower power (e.g., +3.3VDC indicated by 70) for use by the microprocessor 71. The microprocessor 71 controls the motor 14 via power electronics 73, as well as all of the associated indicators, such as LED indicators 3, 4, 75 and 75A. The power electronics 73 convert the microprocessor 71 commands into a power signal to motor 76 (24Vp-p AC) using internal inverters to activate the motor 14. An enable switch 72 is provided to permit the clearing stem to be continuously reciprocated for a predetermined period of time (e.g., 4-20 minutes), which avoids running the device TC1 for too long but provides sufficient time to effect clearing the clog. A control box power switch 2 is coupled to the microprocessor power unit (MPU) 69 via a fuse 66. A power indicator (e.g., LED) 3 is provided on the control box 1. When the control box 1 is externally powered, e.g., from 120VDC, 60 Hz wall power, a power-cord (not shown) is supplied with the control box 1, and which includes an AC/DC converter. It should be understood that this does not limit the operation of the present invention to wall power in any manner and that the control box 1 can be operated off any type of power source, including battery power.

The electronic system 63 may also include a displacement sensor DS (e.g., an LVDT (e.g., Macro Sensors CD 375-500) or force sensor/load cell (e.g., Futek LPM 200); or eddy current sensor (e.g., Micro-Epsilon eddy NCDT 3010), etc.) for accomplishing closed loop motor control as well as detecting changes in the clearing process. For example, the sensor DS forms a closed loop with microprocessor 71 for maintaining the motor shaft 15 in a centered position, which maintains the motor 14 where the force is the greatest and provides optimum control. Alternatively, the sensor DS may comprise a displacement/force feedback sensor or even an optical displacement sensor (e.g., Variohm Eurosens). The DS sensor output may also be used for self-centering of the wire 28 during operation. As part of the closed loop control, it may be advantageous to also change any DC offset to alter the force profile at the wire tip 29 and to provide more power to one side.

In addition, an impedance sensor/current sensor IS may be included for detecting the change in voltage/current of the motor 14 and communicating with the microprocessor 71 for determining the status of the clearing process, such as initial contact with blockage, passage therethrough, etc. This status can be conveyed through a display or clearing status indicator 75 (e.g., LEDs, 7-segment displays, audible indicators, etc.) or a series of differently-colored LEDs 75 (e.g., from green to yellow to red). Alternatively, where the displacement sensor DS comprises a displacement/force feedback sensor, this sensor's output can be used to detect when the clog 40 is contacted and when it is penetrated.

As mentioned earlier, in order to indicate that the clearing stem magnet 33 and the control box magnet 12 are coupled properly, a magnetic/conductive sensor to determine if a solid clearing stem connection has been made which can then be provided to an indicator 75A. By way of example only, a magnetic sensor could be implemented to determine safe connectivity between magnets in operation, such as a Honeywell Magnetometer, HMR2300. These magnetometers measure both magnetic field intensity and direction using their Anisotropic Magneto-Resistive sensors. The ability to acquire this information can be utilized by the microprocessor 71 to ensure the magnet polarities are correct, and that the magnets field intensity is at a safe level (e.g., they have not been de-magnetized). Similarly, an anti-tamper circuit may also be included in the electronic system 63 which interrupts operation if the control box 1 is attempted to be opened. A corresponding tamper sensor may also be provided that causes the indicator 75A on the control box 1 to indicate if someone has opened, or attempted opening the lid of the control box 1. Furthermore, control box screws can be configured to disable operation of the control box 1, if they are attempted to be removed during activation.

The microprocessor 71 can be programmed to drive the electronic system 63 at the needed voltage and frequency, converting 120V 60 Hz wall power to needed parameters to drive the motor 14 at, for example 15-40 Hz (e.g., 25 Hz). In particular, several fault conditions are programmed into the microprocessor 71 for which it interrupts device TC1 operation:

$V_{input} < 20VDC;$

$V_{input} > 25VDC;$

Overtemperature condition pertaining to the amplifier IC;

Short circuit condition pertaining to the amplifier IC;

Should any of these fault conditions occur, the microprocessor 71 activates a fault indicator 4. Also, as discussed earlier, the enable switch 72 permits the operator to initiate the reciprocating motion without the need to hold any trigger. The enable switch 72 permits the control box 1 to maintain the reciprocating motion for a predetermined period of time (e.g., 4-20 minutes) before the reciprocating motion is terminated.

FIG. 17B provides a flow diagram of the microprocessor 71 operation: at step power up 85, the microprocessor 71 is powered up following activation of the power switch 2 by the operator. The microprocessor 71 then conducts a one second step initialization 86. Once the initialization 86 is completed the microprocessor 71 activates the power indicator 3 (e.g., typically a green light (GL) or indication). At this point, device TC1 remains in a disabled state until the enable switch 72 is activated by the operator; "enable button pressed" step 89 of the flow diagram represents activation of the enable switch 72 resulting in the enabled state 88 of the device where the clearing stem 26 is being reciprocated as described previously. The microprocessor 71 then maintains operation of this reciprocation for the predetermined period (e.g., 4-20 minutes) shown as time interval 93 in the flow diagram. At the

end of the predetermined period, the microprocessor 71 terminates the reciprocating movement of the clearing stem 26 and returns to step disabled 87. In addition, upon activation of the enable switch 72 by the operator, the microprocessor 71 monitors the device TC1 for the faults described above, indicated by the paths—fault detected 90 of the flow diagram. If a fault 91 is detected by the microprocessor 71, the microprocessor 71 terminates clearing stem reciprocation and activates the fault indicator 4 (e.g., typically a yellow light (YL) or indication). The microprocessor 71 then shuts down (step power cycle 92) the device TC1.

Operation of the present invention tube clearer TC1 is as follows: if wall power is being used, the connector end of the power cord (not shown) is inserted into power plug 5 (FIGS. 2-2A) on the control box 1 and the other end of the power cord is coupled to a power supply which is coupled to a standard 120V RMS/60 Hz three-prong outlet. The control box 1 is turned on using the power switch 2 which turns on the power indicator 3 which verifies that the control box 1 is operating properly.

A new clearing stem 26 is removed from its packaging (but not discarded since the contaminated clearing stem 26 will be placed in the packaging and then discarded). If a plurality of clearing stems 26 are provided with tube depth-control collars fixed at different positions, the operator needs to select the clearing stem which has the appropriate fixed collar position; if, the tube depth-control collar is adjustable, the operator needs to position the collar appropriately along the clearing stem.

The following discussion of the operation is based upon the control box shown in FIGS. 2-2A, it being understood that this is by way of example only. The wire end of the wire 28 comprising the clearing stem magnet 33 is gently pulled out from within the sheath 30 and then the clearing stem magnet 33 is inserted into the bore of the Luer clearing stem connector 6 until the operator feels the pull of the clearing stem magnet 33 to the other magnet 12 and/or hears the magnets connect. The sheath 30 is then pushed until the Luer clearing stem fitting 32 is flush with the Luer clearing stem connector 6 on the control box 1. The Luer clearing stem fitting 32 is then twisted onto the Luer clearing stem connector 6. Next, the distal end wire tip 29 of the clearing member 26 is inserted a few inches into the artificial tube. The enable switch 72 is pressed to activate the reciprocating motion. While holding the artificial tube 39 in one hand, the clearing stem 26 is held in the other hand while the clearing stem 26 is advanced into the artificial tube. When the clog is initially encountered, the clearing status indicator 75 changes to alert to the initial contact, and the operator begins to apply a slight force to the clearing stem 26. Facilitating clog clearance can be achieved by the operator moving the clearing stem 26 back and forth slightly to clear the clog. These steps are repeated until the clog has cleared, in which case, the clearing status indicator 75 showing that the clog has been cleared activates. If the clog is cleared before the predetermined period (e.g., 4-20 minutes) is reached, the operator can depress the enable switch 72 again to stop the reciprocating movement and then depress the power switch 2 to shut off power to the device TC1. The clearing stem 26 can then be removed from the artificial tube (e.g., feeding tube FT) and then the working end of the clearing stem 26 can be inserted into the packaging. The artificial tube should be flushed with water to verify that the clog has been cleared; if not, the working end of the clearing stem 26 should be removed from the packaging and the clearing procedure repeated. If the clog is verified as being cleared, the clearing stem 26 is disengaged from the control box 1 in accordance with the version of the control box 1 being used.

For example, if the preferred control box 1 (e.g., FIG. 16C) is being used, the alternate clearing stem fitting 32A is disengaged from the sheath attachment bracket 83 and the alternate clearing stem magnet 33A is pulled away from the alternate diaphragm 9A; alternatively, where the Luer fitting version of the control box 1 (e.g., FIG. 16A) is used, the operator twists the Luer clearing stem fitting 32 and removes the clearing stem magnet 33 end of the clearing stem 26 from the control box 1. In either situation, the clearing stem 26 is placed back in the packaging and this is discarded in a suitable biohazard container.

FIG. 29 provides a partial isometric end view of a working end 401 of the wire 28 of the clearing stem 26 which utilizes a sheath with channels 30E that includes ports 402 which can be used for irrigation and/or aspiration. These ports 402 form the end of conduits in the sheath with channels 30E whose other ends are coupled to an aspiration source (not shown, e.g., a vacuum source, etc.) and/or an irrigation source (also not shown, e.g., a saline solution source, or other liquid source). During clog break-up, broken pieces of the clog can be aspirated out of the artificial tube using the sheath with channels 30E and where irrigating the clog vicinity is required, the sheath with channels 30E can be used to deliver such liquids. When aspirating and irrigating simultaneously, aspiration flow should equal irrigation flow rate. The appropriate flow rates are preferably 1-15 ml/min.

Another alternate clearing stem configuration is replacing the wire 28 with a hollow lumen or wire 403 to allow aspiration or irrigation down the hollow lumen or wire 403 to achieve the same purposes discussed with regard to FIG. 29. This alternative configuration is shown in FIG. 29B. Thus, the sheath ports 402 and the hollow lumen or wire 403 may cooperate in different configurations to achieve irrigation/aspiration alternatively or simultaneously. By way of example, the sheath ports 402 can be irrigating while the hollow lumen or wire 403 is suctioning, or vice versa. Alternatively, all of the ports 402 and the hollow lumen or wire 403 can be operating as irrigators or aspiration.

Another alternate clearing stem configuration is to use the indwelling artificial tube 39 effectively as the sheath, as illustrated in FIG. 29C. In this case, a wire 28 or hollow lumen or wire 403 is inserted directly into an artificial tube 39 without the sheath 30. The motor 14 drives the wire 28 or hollow lumen or wire 403 with motion as described previously, to disrupt the clog 40. Although not shown, the tube depth-control collar 22 may also be secured at the desired length to prevent over-insertion of the wire 28 or hollow lumen or wire 403, with the collar 22 impacting the end of open proximal end of the artificial tube 39 during operation. Alternatively, the wire 28 or hollow lumen or wire 403 may include the fixed tube depth-control collar 22A to also limit over-insertion. Using this configuration, the hollow lumen or wire 403 can achieve irrigation or suction alternatively. An advantage of this configuration is that elimination of the sheath can allow access to narrower lumens. The phrase “completely exposed” when used with the device TC1 means a device TC1 that does not use a sheath.

Another alternate clearing stem configuration is a very narrow hollow lumen or wire 403 compared to the sheath 30 such that the areal differential between the hollow lumen or wire 403 and sheath 30 allows for aspiration/irrigation as illustrated in FIG. 29D.

Another alternate clearing stem configuration is the sheath 30 has two ports. One is quite small and is possibly used for a very narrow hollow lumen or wire 403 and the port 402 is used for aspiration/irrigation as illustrated in FIG. 29E.

Rotating Tube Clearer TC2

As with TC1, tube clearer TC2 is a mechanical tube clearer but instead of generating reciprocating motion, tube clearer TC2 generates rotating motion to achieve artificial tube clearing, preferably for PEG feeding tubes. FIG. 18A depicts the tube clearer TC2 which comprises a reusable handset 115 (which remains outside the artificial tube and the patient) having a motor 108 (e.g., a DC motor) that drives (rotates) a disposable or limited-reuse clearing member 114. The handset 115 is held by the operator's hand 136 during the clearing procedure.

It should be noted that, alternatively, clearing member 114 may also be hollow for irrigation or aspiration, or other features.

The tube clearer TC2 (FIG. 19) comprises a clearing member 114 that includes a magnetic connector 103 at one end which attaches to a torque limiter 105 of the handset 115. Attached at the distal end of the clearing member 114 is a narrow flexible rod, preferably a polymer piece of tubing with a clearing brush 101 located on its distal end. The clearing member 114 can be solid or hollow. In the solid embodiment, the distal end of the clearing member 114 is attached to the clearing brush 101 and the proximal end of the clearing member 114 is attached to a magnetic connector 103. In the hollow embodiment, the wire holding the clearing brush 101 may extend the central length of the clearing member 114 to the magnetic connector 103. The clearing member 114 is flexible in order to conform to various radius of curvatures R. It is rotated by the motor 108 within the handset 115. The rotary motion of the clearing brush 101 clears the clog, occlusion, or debris from the tube (not shown).

Clearing Member and Connectors

The clearing member 114 comprises a polymer tube with a clearing brush 101 inset at its distal end. The preferred polymer materials are nylon and polyurethane, although other materials may be used, such as polytetrafluoroethylene (PTFE), Polyvinyl chloride (PVC), polyethylene, polypropylene, and fluoropolymer. The length of the clearing member 114 is equal to the length of the feeding tube +/-one inch, depending on application. FIG. 22 shows the layout of the clearing member 114. At the proximal end of the clearing member 114 is a polymer magnetic connector 103 which includes a clearing member magnet adapter 104 in its inner bore and which sits flush to the proximal end of the clearing member 114. To attach the clearing member 114 to the handset, as shown in FIG. 19, the magnetic connector 103 is inserted into a receiving bore 105A within the torque limiter 105 of the handset 115. Disposed within the bore end is a magnetic element 105B and wherein when the magnetic connector 103 is inserted into the receiving bore 105A, the clearing member magnet adapter 104 and magnetic element 105B contact. To facilitate a tight connection, the magnetic connector 103 comprises a hexagonal-shape, or other non-round shape, that fits into a correspondingly-shaped receiving bore 105A. DC motor 108 output is conveyed to the clearing member stem 102 through a gear train 107 and gear train output shaft 106.

The clearing brush 101 at the distal end has several unique features. It could be a twisted-in-wire type clearing brush 101 with a negative taper NT, as shown in FIG. 23. By way of example only, the clearing brush 101 may comprise a twisted-in wire type; alternatively, the brush 101 may comprise a helical-wound wire or other type brush design. "Negative taper" implies that the clearing brush 101 bristles are wider in diameter at the distal end than at the proximal end of the clearing brush 101. There are several reasons for this configuration in the clearing member's 114 design. Most conven-

tional brushes have a taper smaller at the distal end and larger at the proximal end. However, for this application it would require over-insertion to clear the full bore of the end of the artificial tube (e.g., feeding tube) 119. The negative taper NT also allows the helix-type wound clearing brush 101 to be extended rearward, as shown by the path of freed clog particles arrow 120 in FIG. 23. When rotating (indicated by the rotation of brush arrow 121), this clearing brush 101 design forces wicking of the loosened clog debris away from the clog 122 also in the direction of the path of freed clog particles arrow 120. This is important for fast, effective clearing. If the clog 122 was not removed from the clog site, it could be compacted further, making the clog 122 even more difficult to remove. The negative taper NT also allows for contact with the tube walls (in order to clean them), but only in a limited area. Having contact only in a limited area reduces the amount of drag on the artificial tube 119 and the torque transmitted to it and thus this minimizes any chance of dislodging the artificial tube 119 from within the patient when the clearing member 114 is removed from the artificial tube 119. The shape of the (distal) tip of the clearing brush 101 is also important for this application. Unlike many standard twisted-in-wire brushes, which are cut at the ends after twisting, the TC2 clearing brush 101 could possibly be wound with a rounded tip—the wire bends 180 degrees. This bend prevents any sharp end from coming into contact with the stomach, intestine, or other organs/tissues if over-inserted past the end of the artificial tube 119. Thus, the clearing brush 101 transfers minimal torque due to its unique geometry, but its helical design is also able to remove loosened debris from the clog 122.

In another embodiment, the brush tip 101A (FIG. 19) radius of the clearing brush 101 can be modified, e.g., rounded to allow the clearing brush 101 to break up a clog, but to not penetrate an organ (e.g., stomach or other tissue/organ, etc.) should the brush tip 101A ever make its way close to an organ. The clearing brush 101 may also be retracted from the distal end of the clearing member to decrease the chance of the clearing brush 101 catching in stomach or other tissue. In another embodiment, the brush tip 101A can be modified by the addition of a flexible tip such as a Tecoflex® tip. In another embodiment, brush tip 101A can be modified by the addition of ball tip 34E as illustrated in FIG. 5D.

Handset 115

Preferably, the handset 115 is shaped like a pistol, with contours to fit the user's fingers comfortably while he/she is using it, as shown by the operator's hand 136 (FIG. 18A). An index finger trigger 109 controls operation. The trigger 109 is a momentary power switch that only provides power when being pressed. The handset 115 is composed of three parts, one battery cover and two halves which are fastened together by screws or built-in snap fit connectors to form a handset housing 113. It also contains an isolated battery compartment 112 to facilitate battery 111 changes without exposing any components to contaminants that could cause device failure or reduce reliability. A control circuit 110 (FIG. 19) conveys power to the DC motor 108.

In this embodiment the handset contains an isolated compartment in which a common battery size is used. For example, the handset 115 can be designed to accommodate any battery size such as 9V, AA, AAA, or a specialty size and a plurality of batteries where required. Alternatively, the handset 115 may comprise a rechargeable battery such that there is no need to remove any batteries. A charger (not shown) may accompany the handset 115 such that the rechargeable battery can be inductively charged and this configuration has advantages over the battery operated setup,

including: no panels are removable on the handset **115** which eliminates the possibility of contamination; and also reduces cost and disposal of batteries. The inductive charger may comprise a base unit, rechargeable battery, and circuitry. The base unit may comprise an enclosure with a slot or depression or cradle into which the handset **115** is positioned. The base unit plugs into a standard 120V outlet. A coil in the base unit transmits a magnetic field to a coil in the handset **115**, and a charging circuit would transform the signal to the correct voltage and route it to the rechargeable battery located in the handset **115**.

Motor

The motor **108** of the tube clearer TC2 is preferably a DC motor or a brushless DC motor and gear combination. The gear mechanism may be a precision gear head, such as one utilizing a planetary gear train **116** or a compound gear train **118** utilizing two or more standalone gears. Motor and gear output speed ranges from 600 RPM to 1800 RPM, more preferably 740 to 1140 RPM. The torque limiter **105** is also preferred in this embodiment. The maximum output torque can preferably range from 20 mNm to 40 mNm with a more preferable torque of 24 to 34 mNm. A voltage of less than or equal to 9 volts DC is preferred to drive the motor **108**, such that standard commercially-available batteries can be used. FIG. **20** shows a DC motor **108** with a planetary gear train **116** whereas FIG. **21** shows a DC motor **108** with a compound gear train **118** configuration that is coupled to the motor output shaft **117**. Thus, torque, speed and geometry of the clearing stem define the optimal operation of the device TC2. Alternatively, the motor **118** itself may have a torque output of preferably 20 mNm to 40 mNm, with a more preferable torque of 24 to 34 mNm, in which case the torque limiter **105** would not be necessary.

In another embodiment, a DC or brushless DC motor **108** and gear combination is used in combination with a torque limiter **105**. The torque limiter **105** is attached in-line with the motor output shaft **117** and allows slippage once the maximum output torque is reached. In another embodiment, a DC or brushless DC motor **108** and gear combination is used in combination with a hammering device, similar to that found in hammer drills (U.S. Pat. No. 5,653,294 (Thurler, et al.) and whose entire disclosure is incorporated by reference herein). This device creates an oscillatory motion along with the rotary motion to clear the clog. In another embodiment, the DC or brushless DC motor in all examples above is replaced with a piezoelectric motor with similar specifications.

Tube Depth-Control Collar

As with TC1, tube clearer TC2 comprises a tube depth-control collar **133**, as shown in FIG. **24**. This depth-control collar permits one-handed operation using no special tools. The tube depth-control collar **133** mounts along the rod portion of the clearing member **114**. The tube depth-control collar **133** is formed to be well-balanced and lightweight so as to not cause unwanted harmonics in the clearing member **114** during rotation. The tube depth-control collar **133** comprises a lightweight, circular tube depth-control collar housing **129** which includes a displaceable tube depth-control collar push button **130** that acts against a preloaded spring **132** bias and which locks against the clearing member **114** which passes through the opening for clearing member **131**. FIG. **18A** depicts the tube depth-control collar **133** on the clearing member **114**.

Motor Torque Limiting

In a preferred embodiment of the handset **115**, the torque applied to the clearing member **114** is limited by controlling the voltage and current applied to the DC motor and ultimately to the gears. These voltage and current limits are

established by testing and determining the minimum angle of twist that are unacceptable when the clearing brush **101** is in a locked condition within tubes under test. An alternative method involves the use of a DC motor with a torque limiter **105** as depicted in FIGS. **19** and **25**. The torque limiter **105** is a two-piece patterned disc, preloaded by a preload spring **125**. The spring force controls torque at which disc slippage occurs. In particular, the torque limiter **105** comprises an input coupler **123**, a torque limiter output shaft **135**, a preload collar **134** and a torque limiter profile **124**. The input coupler **123** couples to the gear train **107** and the torque limiter output shaft **135** couples the clearing member **114**. As can be appreciated, when a certain applied torque is exceeded, the torque limiter **105** is designed to slip at the interface or torque limiter profile **124** to disengage and thereby prevent the clearing member **114** from exceeding the torque limit.

Clearing Member Control

The tube clearer TC2 must control harmonics so that the clearing member **114** does not become uncontrollable and cause injury/damage. During device activation, the tube clearer TC2 rotates the clearing member **114** with a displacement diameter that is preferably from 0 mm to 40 mm and a more preferred diameter of 25.4 mm or less. FIG. **26** shows multi-nodal harmonics (i.e., node points **126**) occurring in the clearing member **114** while spinning and also depicts the maximum desired displacement **127A**. This is preferred as its shape limits the displacement by geometry. The distance between the first two nodal points **126** is indicated by distance between nodal points **128**, and as can be seen in FIG. **26**, this distance decreases for subsequent nodal points **126**. The maximum desired displacement **127A** of the clearing stem is preferred to be kept to 25.4 mm or less. In contrast, FIG. **27** depicts a commercially-available rotary tool **115A** (e.g., a hand-held drill) rotating the clearing member **114**, showing the undesirable profile of rotating stem **129A** (and its undesirable corresponding maximum radial displacement **127B**) of the clearing stem motion because there is only one nodal point at the proximal end of the clearing member **114**. This type of deformation is not preferred because it is more likely to be unstable.

FIG. **28** depicts a block diagram of the electronics of the device TC2. In particular, a DC motor **108** provides the rotational motion to the clearing stem **114**. The motor **108** receives its input voltage **140** from a voltage regulator **137** which in turn receives power **139** from a power source or battery **111** (e.g., 9V battery, a rechargeable battery, etc.) when the trigger **109** is activated by the operator. A power indicator **138** (see FIG. **18A** also), driven by the voltage regulator, is also provided.

FIG. **29A** provides a partial isometric end view of the device TC2 showing the clearing brush **101** coupled to the clearing member stem **102** which utilizes a sheath with channels **30E** that includes ports **402** which can be used for irrigation and/or aspiration. These ports **402** form the end of conduits in the sheath with channels **30E** whose other ends are coupled to an aspiration source (not shown, e.g., a vacuum source, etc.) and/or an irrigation source (also not shown, e.g., a saline solution source, or other liquid source). During clog break-up, broken pieces of the clog can be aspirated out of the artificial tube using the sheath with channels **30E** and where irrigating the clog vicinity is required, the sheath with channels **30E** can be used to deliver such liquids. When aspirating and irrigating simultaneously, aspiration flow should equal irrigation flow rate. The appropriate flow rates are preferably between 1-15 mL/min. The clearing brush **101** can also be placed back along the clearing member stem **102** away from the distal end of the clearing member **114** to decrease the

potential for the clearing brush 101 grabbing or interacting with the stomach or other organ or tissue. Alternatively, the various configurations shown in FIGS. 29 and 29B-29E can also be used with the device TC2. The phrase “completely exposed” when used with the device TC2 means a device TC2 that does not use a sheath.

FIGS. 18B-18C depict an alternative voice coil motor tube clear device TC2. Instead of using a “pistol-style” housing, the device TC2 of FIGS. 18B-18C comprise an elongated hand grip 301. In addition, unlike the rotational motion of the TC2 device shown in FIG. 18A, the alternative voice coil motor tube clear device 300 generates reciprocating motion (as discussed previously with regard to the TC1 devices). In particular, within the hand grip 301 is positioned a voice coil motor 305 that, when energized, causes the clearing stem 303 to reciprocate. The tip of the clearing stem 303 includes a clearing brush 304. As shown most clearly in FIG. 18B, a clearing stem adapter 302 is provided on an end of the hand grip 301 for securing the clearing stem 303 to the voice coil motor 305 in the hand grip 301. A power indicator 138 is also provided to indicate when power is being provided to the clearing stem 303 for reciprocating motion. A power switch/trigger 109A is provided so that the user can manually control the activation of the device, similar to the pistol-style embodiment.

It should be noted that, alternatively, clearing stem 303 may also be hollow for irrigation or aspiration, or other features and may have similar configurations as shown in FIGS. 29-29E.

It should be further understood that the preferred embodiments of the present invention are for the in-situ clearing of artificial lumens in a living being, but that these embodiments can be used for clearing lumens located outside of the living being, as well as for clearing other types of lumens not associated with living beings.

Now that exemplary embodiments of the present invention have been shown and described in detail, various modifications and improvements thereon will become readily apparent to those skilled in the art. Accordingly, the spirit and scope of the present invention is to be construed broadly and limited only by the appended claims, and not by the foregoing specification.

APPENDIX

Reference Characters and Their Associations		
APA	Amplified Piezoelectric Actuator	TC1
CR	Crank	TC1
DS	Displacement Sensor	TC1
IS	Impedance Sensor	TC1
NT	Negative Taper Angle	TC2
SY	Scotch Yoke	TC1
SS	Soft Stop	TC1/TC2
TC1	Tube Clearing Device 1	TC1
TC2	Tube Clearing Device 2	TC2
TCS	Tip Compression Spring	TC1
GL	Green Light	TC1
YL	Yellow Light	TC1
FT	Artificial/Feeding Tube	TC1/TC2
R	Radius of Curvature	TC1
1	Control Box	TC1
2	Power Switch	TC1
3	Power Indicator	TC1
4	Fault Indicator	TC1
5	Power Plug	TC1
6	Clearing Stem Connector	TC1
7	Motor Mount	TC1
8	Motor Mount Damper	TC1
9	Diaphragm	TC1
9A	Alternate Diaphragm	TC1

APPENDIX-continued

Reference Characters and Their Associations			
10	Electronics		TC1
11	Motor PCB		TC1
12	Magnet		TC1
12A	Magnet Recess		TC1
13	Motor Magnetic Coupler		TC1
13A	Alternate Motor Magnetic Coupler		TC1
14	Motor		TC1
14A	Counter Balance Mechanism		TC1
15	Motor Shaft		TC1
16	VCM Body		TC1
17	Winding		TC1
18	End Bearing		TC1
19	Spring		TC1
20	Magnets		TC1
20N-20S	Magnetic Driving members		TC1
21A-21C	Pole Pieces		TC1
22	Tube Depth-Control collar		TC1
22A	Fixed Tube Depth-Control collar		TC1
23	Depth Control Collar Push Button		TC1
23A	Central passageway of push button		TC1
23B	Lower portion of press button		TC1
24	Tube Depth-Control Collar Body		TC1
24A	Central passageway of collar body		TC1
24B	Upper portion of collar body		TC1
25	Spring		TC1
26	Clearing Stem		TC1
27	Wire Stop		TC1
27A	Alternate Wire Stop		TC1
28	Wire		TC1
28A	Wire Protusion		TC1
29	Wire Tip		TC1
30	Sheath		TC1
30A	Sheath length markings		TC1
30B	Integer markings		TC1
30C	Distal End		TC1
30D	Proximal End		TC1
30E	Sheath with Channels		TC1
31	Stem Stiffener		TC1
32	Clearing Stem Fitting		TC1
32A	Alternate Clearing Stem Fitting		TC1
33	Clearing Stem Magnet		TC1
33A	Alternate Clearing Stem Magnet		TC1
33B	Alternate Clearing Stem Magnet Fitting		TC1
34	Plastic Wire Tip		TC1
34A	Alternate Tubing Tip		TC1/TC2
34B	Fixed Member		TC1/TC2
34C	Gripping/Chopping Mechanism		TC1
34D	Pivot Point		TC1
34E	Ball Tip		TC1/TC2
35	Wire Tip Brush		TC1
36	Sheath Tip Brush		TC1
37	Forward Swept Sheath Tip Brush		TC1
38	Nursing Cart		TC1
38A	Pole		TC1
39	Artificial Tube		TC1
40	Clog		TC1
41	Tube Inner Lumen		TC1
42	Pneumatic Motor		TC1
43	Pneumatic Motor Housing		TC1
44	Pneumatic Motor Shaft		TC1
46	Pneumatic Motor Diaphragm		TC1
47	Internal Tubing		TC1
48	Scotch Yoke Motor		TC1
49	DC Motor		TC1
50	Scotch Yoke Slider		TC1
50A	Scotch Yoke Forward Displacement direction		TC1
50B	Scotch Yoke Rearward Displacement direction		TC1
51	Adapter		TC1
52	Scotch Yoke Shaft		TC1
53	Wires		TC1
54	Air Supply Inlet		TC1
55	Solenoid Motor		TC1
56	Solenoid		TC1
57	Solenoid Shaft		TC1
58	Return Spring		TC1

APPENDIX-continued

Reference Characters and Their Associations		
59	APA Motor	TC1
60	Actuator	TC1
61	Actuator Mount	TC1
62	Actuator Shaft	TC1
63	Electronic System	TC1
66	Fuse	TC1
67	Power	TC1
69	Micro Processor Power Unit (MPU)	TC1
70	+3.3 VDC	TC1
71	Microprocessor	TC1
72	Enable Switch	TC1
73	Power Electronics	TC1
75	Clearing Status Indicator	TC1
75A	Indicator	TC1
76	power signal to motor	TC1
77	Langevin Transducer motor	TC1
78	Piezoelectric elements	TC1
79	Pre-stress bolt	TC1
80	Tail Mass	TC1
81	Horn	TC1
82	Clearing Stem Attachment	TC1
83	Sheath Attachment Bracket	TC1
84	Diaphragm Sealing Ring	TC1
85	Power Up	TC1
86	Initialization	TC1
87	Disabled	TC1
88	Enabled	TC1
89	Enable Button Pressed	TC1
90	Fault Detected	TC1
91	Fault	TC1
92	Power Cycle	TC1
93	Time Interval	TC1
101	Clearing Brush	TC2
101A	Brush tip	TC2
102	Clearing Member Stem	TC2
103	Magnetic Connector	TC2
104	Magnetic Adapter	TC2
105	Torque Limiter	TC2
105A	Receiving Bore	TC2
105B	Magnetic Element	TC2
106	Gear Train Output Shaft	TC2
107	Gear Train	TC2
108	Motor	TC2
109	Trigger	TC2
109A	Power Switch/trigger	TC2
110	Control Circuit	TC2
111	Battery	TC2
112	Battery Compartment	TC2
113	Handset Housing	TC2
114	Clearing Member	TC2
115	Handset	TC2
115A	Commercial Available Rotary Tool	TC2
116	Planetary Gear Train	TC2
117	Motor Output Shaft	TC2
118	Compound Gear Train	TC2
119	Artificial Tube	TC2
120	Path of Freed Clog Particles	TC2
121	Rotation of Brush Arrow	TC2
122	Clog	TC2
123	Input Coupler	TC2
124	Torque Limiter Profile	TC2
125	Preload Springs	TC2
126	Nodal Points	TC2
127A	Maximum Desired Displacement	TC2
127B	Undesirable Displacement	TC2
128	Distance between nodal points	TC2
129	Tube depth-control collar housing	TC2
129A	Undesired Profile of Rotating Stem	TC2
130	Tube Depth-Control Collar Push Button	TC2
131	Opening for Clearing Member	TC2
132	Preloaded Spring	TC2
133	Tube depth-control collar	TC2
134	Preload Collar	TC2
135	Torque Limiter Output Shaft	TC2
136	Operator's Hand	TC2
137	Voltage Regulator	TC2
138	Power Indicator	TC2

APPENDIX-continued

Reference Characters and Their Associations		
5	139	Power TC2
	140	Input Voltage TC2
	300	Voice Coil Motor (VCM) Tube Clear TC2
	301	Hand Grip TC2
	302	Clearing Stem Adapter TC2
	303	Clearing Stem TC2
	304	Clearing Brush TC2
10	305	Voice Coil Motor TC2
	401	Working End TC1/TC2
	402	Port TC1/TC2
	403	Hollow Lumen or Wire TC1/TC2
15	While the invention has been described in detail and with reference to specific examples thereof, it will be apparent to one skilled in the art that various changes and modifications can be made therein without departing from the spirit and scope thereof.	
20	What is claimed is:	
	1. A device for the in situ clearing of blockages in artificial tubes completely or partially disposed within a living being, said device comprising:	
25	a controller that remains outside of the living being, said controller comprising an actuator for generating repetitive motion; and	
	a clearing member having:	
30	a first end that is releasably coupled to the actuator,	
	a second working end that is insertable into an opening in the artificial tube, said second working end having a portion that comes into repetitive contact with a blockage in the artificial tube for clearing the blockage therein, and	
35	a flexible material that permits said clearing member to make repetitive contact with the blockage while said clearing member is positioned within a straight portion or within a curved portion of the artificial tube,	
	a stationary sheath,	
40	a movable wire disposed within said stationary sheath, said wire having a first wire end that is coupled to said actuator and a second wire end that comprises a tip that contacts said blockage, and	
	a pliant wire stop through which said wire is disposed, said pliant wire stop coupled between said first wire end and a stationary fitting through which said wire passes,	
45	wherein said actuator provides reciprocating motion to said clearing member that occurs along a longitudinal axis of the artificial tube.	
50	2. The device of claim 1 further comprising a tube depth-control collar along a length of said clearing member.	
	3. The device of claim 2 wherein said tube depth-control collar is fixed on said clearing member.	
55	4. The device of claim 2 wherein said tube depth-control collar is adjustable along said clearing member.	
	5. The device of claim 1 wherein said controller further comprises a displaceable diaphragm that seals a portion of a controller housing, said actuator comprising a shaft having a first shaft end that is in contact with a first side of said diaphragm, and	
60	wherein said first wire end is in contact with a second side of said diaphragm, said first shaft end and said first end of said movable wire being magnetically coupled together through said diaphragm.	
65	6. The device of claim 1 wherein said tip comprises a tip compression spring that is coupled to a first end of a tip tubing	

29

having a second end that is open to permit an extreme tip end of said wire to extend out of, or be retracted within, said tip tubing.

7. The device of claim 1 wherein said tip comprises a gripper that forms an extreme end of said tip for gripping the blockage or portions thereof.

8. The device of claim 1 wherein said tip comprises a brush at an extreme tip end.

9. The device of claim 1 wherein said tip comprises an extreme tip end in the form of a sphere.

10. The device of claim 1 wherein said stationary sheath comprises a first sheath end having a brush positioned thereat.

11. The device of claim 10 wherein said brush comprises bristles that are swept forward.

12. The device of claim 1 wherein said movable wire is hollow for permitting irrigation or aspiration at said second wire end.

13. The device of claim 1 wherein said stationary sheath comprises a first sheath end having ports that permit irrigation or aspiration at said second working end through said ports when said ports are coupled to irrigation/aspiration sources.

14. The device of claim 1 wherein said actuator operates in the range of approximately 15-40 Hz.

15. The device of claim 1 wherein the tip of the movable wire has a linear displacement of approximately 2-10 mm.

16. The device of claim 1 further comprising a clearing status indicator that detects changes in power input to said actuator for providing the operator with an indication as to progress with a clearing process.

17. The device of claim 1 further comprising a clearing status indicator, said clearing status indicator detecting that a first shaft end of said actuator and said first end of said movable wire are magnetically coupled together.

18. The device of claim 1 further comprising a tamper indicator, said tamper indicator providing an alert to the operator of tampering with the controller or if said clearing member is not properly connected to an actuator shaft.

19. The device of claim 1 further comprising a negative taper brush disposed on an extreme tip of said clearing member, said negative taper brush comprising bristles that are wider in diameter at a distal end of said brush than at a proximal end of said brush.

20. The device of claim 1 wherein the stationary sheath further comprises a first working end having ports that permit irrigation or aspiration at said second working end through said ports when said ports are coupled to irrigation/aspiration sources.

21. The device of claim 1 wherein said clearing member is a hollow element that permits irrigation or aspiration at said second working end.

22. A device for the breaking-apart of occlusions within a living being, comprising:

a controller comprising an actuator for generating repetitive motion;

a clearing member having a first end that is releasably coupled to the actuator and having a second working end that comes into repetitive contact with an occlusion, and wherein the clearing member comprises a flexible wire configured to be caused to reciprocate by the actuator's repetitive motion and a sheath configured to not be caused to reciprocate by the actuator's repetitive motion; wherein the clearing member further comprises, at its first end,

a stationary fitting;
a tube; and

30

a second fitting configured to be caused to reciprocate by the actuator's repetitive motion; and wherein the wire passes through the tube and is coupled to the second fitting.

23. A device for the in situ clearing of blockages in artificial tubes completely or partially disposed within a living being, said device comprising:

a controller that remains outside of the living being, said controller comprising an actuator for generating repetitive motion; and

a clearing member having:

a first end that is releasably coupled to the actuator, a second working end that is insertable into an opening in the artificial tube, said second working end having a portion that comes into repetitive contact with a blockage in the artificial tube for clearing the blockage therein, and

a flexible material that permits said clearing member to make repetitive contact with the blockage while said clearing member is positioned within a straight portion or within a curved portion of the artificial tube, a stationary sheath, and

a movable wire disposed within said stationary sheath, said wire having a first wire end that is coupled to said actuator and a second wire end that comprises a tip that contacts said blockage,

wherein said controller comprises a displaceable diaphragm that seals a portion of a controller housing, wherein said actuator provides reciprocating motion to said clearing member that occurs along a longitudinal axis of the artificial tube and comprises a shaft having a first shaft end that is in contact with one side of said diaphragm, and

wherein said first wire end is in contact with the other side of said diaphragm, said first shaft end and said first end of said movable wire being magnetically coupled together through said diaphragm.

24. The device of claim 23 further comprising a tube depth-control collar along a length of said clearing member.

25. The device of claim 24 wherein said tube depth-control collar is fixed on said clearing member.

26. The device of claim 24 wherein said tube depth-control collar is adjustable along said clearing member.

27. The device of claim 23 further comprising a rigid wire stop coupled to said wire adjacent said second wire end, said rigid wire stop having an outer dimension that is greater than an outer dimension of said sheath.

28. The device of claim 23 wherein said tip comprises a tip compression spring that is coupled to a first end of a tip tubing having a second end that is open to permit an extreme tip end of said wire to extend out of, or be retracted within, said tip tubing.

29. The device of claim 23 wherein said tip comprises a gripper that forms an extreme end of said tip for gripping the blockage or portions thereof.

30. The device of claim 23 wherein said tip comprises a brush at an extreme tip end.

31. The device of claim 23 wherein said tip comprises an extreme tip end in the form of a sphere.

32. The device of claim 23 wherein said stationary sheath comprises a first sheath end having a brush positioned thereat.

33. The device of claim 32 wherein said brush comprises bristles that are swept forward.

34. The device of claim 23 wherein said movable wire is hollow for permitting irrigation or aspiration at said second wire end.

31

35. The device of claim 23 wherein said stationary sheath comprises a first sheath end having ports that permit irrigation or aspiration at said second working end through said ports when said ports are coupled to irrigation/aspiration sources.

36. The device of claim 23 wherein said actuator operates in the range of approximately 15-40 Hz.

37. The device of claim 23 wherein the tip of the movable wire has a linear displacement of approximately 2-10 mm.

38. The device of claim 23 further comprising a clearing status indicator that detects changes in power input to said actuator for providing the operator with an indication as to progress with a clearing process.

39. The device of claim 23 further comprising a clearing status indicator, said clearing status indicator detecting that a first shaft end of said actuator and said first end of said movable wire are magnetically coupled together.

40. The device of claim 23 further comprising a tamper indicator, said tamper indicator providing an alert to the operator of tampering with the controller or if said clearing member is not properly connected to an actuator shaft.

41. The device of claim 23 further comprising a negative taper brush disposed on an extreme tip of said clearing member, said negative taper brush comprising bristles that are wider in diameter at a distal end of said brush than at a proximal end of said brush.

42. The device of claim 23 wherein the stationary sheath further comprises a first working end having ports that permit irrigation or aspiration at said second working end through said ports when said ports are coupled to irrigation/aspiration sources.

43. The device of claim 23 wherein said clearing member is a hollow element that permits irrigation or aspiration at said second working end.

44. A device for the in situ clearing of blockages in artificial tubes completely or partially disposed within a living being, said device comprising:

a controller that remains outside of the living being, said controller comprising an actuator for generating repetitive motion; and

a clearing member having:

a first end that is releasably coupled to the actuator,

a second working end that is insertable into an opening in the artificial tube, said second working end having a portion that comes into repetitive contact with a blockage in the artificial tube for clearing the blockage therein, and

a flexible material that permits said clearing member to make repetitive contact with the blockage while said clearing member is positioned within a straight portion or within a curved portion of the artificial tube, a stationary sheath, and

a movable wire disposed within said stationary sheath, said wire having a first wire end that is coupled to said actuator and a second wire end that comprises a tip that contacts said blockage, said tip comprising a tip compression spring that is coupled to a first end of a tip tubing having a second end that is open to permit an extreme tip end of said wire to extend out of, or be retracted within, said tip tubing,

wherein said actuator provides reciprocating motion to said clearing member that occurs along a longitudinal axis of the artificial tube.

45. The device of claim 44 further comprising a tube depth-control collar along a length of said clearing member.

46. The device of claim 45 wherein said tube depth-control collar is fixed on said clearing member.

32

47. The device of claim 45 wherein said tube depth-control collar is adjustable along said clearing member.

48. The device of claim 44 further comprising a rigid wire stop coupled to said wire adjacent said second wire end, said rigid wire stop having an outer dimension that is greater than an outer dimension of said sheath.

49. The device of claim 44 wherein said tip comprises a gripper that forms an extreme end of said tip for gripping the blockage or portions thereof.

50. The device of claim 44 wherein said tip comprises a brush at an extreme tip end.

51. The device of claim 44 wherein said tip comprises an extreme tip end in the form of a sphere.

52. The device of claim 44 wherein said stationary sheath comprises a first sheath end having a brush positioned thereat.

53. The device of claim 52 wherein said brush comprises bristles that are swept forward.

54. The device of claim 44 wherein said movable wire is hollow for permitting irrigation or aspiration at said second wire end.

55. The device of claim 44 wherein said stationary sheath comprises a first sheath end having ports that permit irrigation or aspiration at said second working end through said ports when said ports are coupled to irrigation/aspiration sources.

56. The device of claim 44 wherein said actuator operates in the range of approximately 15-40 Hz.

57. The device of claim 44 wherein the tip of the movable wire has a linear displacement of approximately 2-10 mm.

58. The device of claim 44 further comprising a clearing status indicator that detects changes in power input to said actuator for providing the operator with an indication as to progress with a clearing process.

59. The device of claim 44 further comprising a clearing status indicator, said clearing status indicator detecting that a first shaft end of said actuator and said first end of said movable wire are magnetically coupled together.

60. The device of claim 44 further comprising a tamper indicator, said tamper indicator providing an alert to the operator of tampering with the controller or if said clearing member is not properly connected to an actuator shaft.

61. The device of claim 44 further comprising a negative taper brush disposed on an extreme tip of said clearing member, said negative taper brush comprising bristles that are wider in diameter at a distal end of said brush than at a proximal end of said brush.

62. The device of claim 44 wherein the stationary sheath further comprises a first working end having ports that permit irrigation or aspiration at said second working end through said ports when said ports are coupled to at least one of irrigation and aspiration sources.

63. The device of claim 44 wherein said clearing member is a hollow element that permits irrigation or aspiration at said second working end.

64. A device for the in situ clearing of blockages in artificial tubes completely or partially disposed within a living being, said device comprising:

a controller that remains outside of the living being, said controller comprising an actuator for generating repetitive motion; and

a clearing member having:

a first end that is releasably coupled to the actuator,

a second working end that is insertable into an opening in the artificial tube, said second working end having a portion that comes into repetitive contact with a blockage in the artificial tube for clearing the blockage therein, and

33

a flexible material that permits said clearing member to make repetitive contact with the blockage while said clearing member is positioned within a straight portion or within a curved portion of the artificial tube, a stationary sheath comprising a first sheath end having a brush positioned thereat, a movable wire disposed within said stationary sheath, said wire having a first wire end that is coupled to said actuator and a second wire end that comprises a tip that contacts said blockage, and wherein said actuator provides reciprocating motion to said clearing member that occurs along a longitudinal axis of the artificial tube.

65. The device of claim 64 wherein said brush comprises bristles that are swept forward.

34

66. The device of claim 64 further comprising a tube depth-control collar along a length of said clearing member.

67. The device of claim 66 wherein said tube depth-control collar is fixed on said clearing member.

68. The device of claim 66 wherein said tube depth-control collar is adjustable along said clearing member.

69. The device of claim 64 wherein said stationary sheath comprises a first sheath end having ports that permit irrigation or aspiration at said second working end through said ports when said ports are coupled to irrigation/aspiration sources.

70. The device of claim 64 wherein said clearing member is a hollow element that permits irrigation or aspiration at said second working end.

* * * * *

UNITED STATES PATENT AND TRADEMARK OFFICE
CERTIFICATE OF CORRECTION

PATENT NO. : 8,262,645 B2
APPLICATION NO. : 12/964252
DATED : September 11, 2012
INVENTOR(S) : Roger B Bagwell et al.

Page 1 of 1

It is certified that error appears in the above-identified patent and that said Letters Patent is hereby corrected as shown below:

In the Specifications

In col. 1 line 25, the phrase “may have” should be replaced with the word --has--.

Signed and Sealed this
Twentieth Day of May, 2014



Michelle K. Lee
Deputy Director of the United States Patent and Trademark Office